

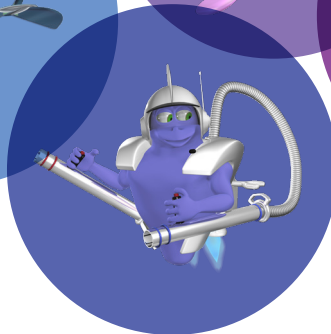
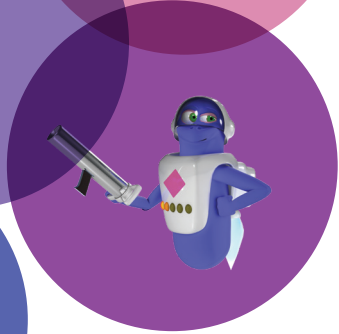
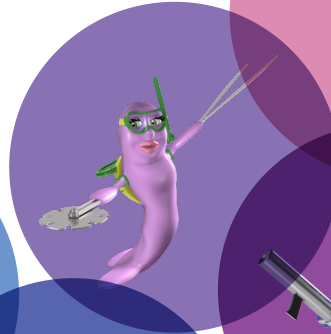
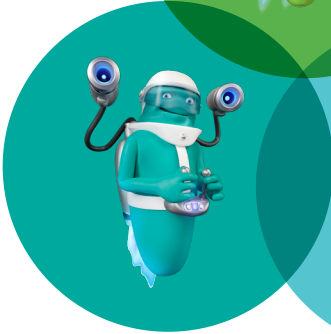
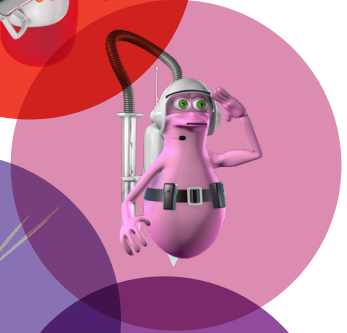
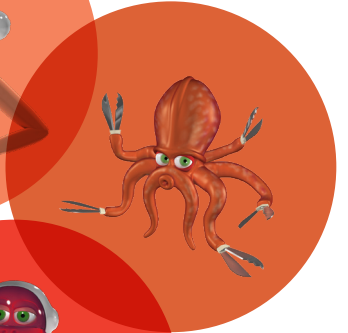


# SmartEnzymes™ for Better Biologics

Antibody Digestion  
Antibody Deglycosylation  
Antibody Conjugation  
Proteomics  
Glycan Profiling  
Affinity Purification  
Research Antibodies

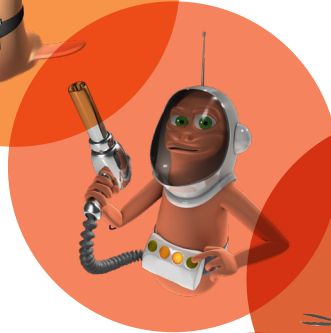

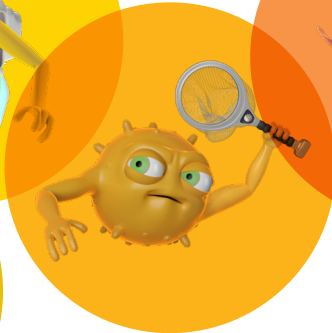

Nature offers a vast source of enzymes, perfected through evolution to perform defined reactions.

At Genovis, we believe that enzymes with unique properties can be used as biological tools to support the research and development of complex biopharmaceuticals to help bring safe and effective medicines to patients in need.

Our task is to identify new enzymes and give them names.

We call them *SmartEnzymes*.

## Antibody Digestion

IgG-specific proteases generate subunits suitable for analytical and preparative workflows for the characterization, development and production of biopharmaceuticals.

<b>FabRICATOR®</b> (IdeS) Below hinge digestion of IgG	6
<b>FabALACTICA®</b> (IgdE) Above hinge digestion of human IgG1	14
<b>FabDELLO™</b> Above hinge digestion of human IgG1	18
<b>FabRICATOR® Z</b> (IdeZ) Below hinge digestion of mouse IgG	20
<b>FabULOUS™</b> (SpeB) Above hinge digestion of IgG	22
<b>GingisKHAN®</b> (Kgp) Above hinge digestion of human IgG1	24

## Antibody Deglycosylation

IgG-specific endoglycosidases trim Fc glycoforms on native antibodies, including high-mannose, hybrid, complex and bisecting type glycans.

<b>GlycINATOR®</b> (EndoS2) Hydrolysis of IgG Fc glycans	26
<b>IgGZERO®</b> (EndoS2) Hydrolysis of complex Fc glycans	28

## Antibody Conjugation

Technologies for site-specific conjugation of antibodies and glycoengineering of IgG with defined glycoforms.

<b>GlyCLICK®</b> Site-specific conjugation of native IgG	30
<b>TransGLYCIT™</b> Transglycosylation of IgG	34

## Proteomics

Protease with unique specificity generates peptides suitable for proteomics and MS analysis.

<b>GingisREX®</b> (RgpB) Arginine-specific protein digestion	38
---	----

## Glycan Profiling

The analysis of N- and O-linked glycosylated proteins are facilitated using a suite of glycan-specific proteases, endo- and exoglycosidases.

<b>PNGase F</b> Hydrolysis of N-glycans	40
<b>OpeRATOR®</b> O-glycan-specific protein digestion	44
<b>OglyZOR®</b> Hydrolysis of core-1 O-glycans	46
<b>SialEXO®</b> Hydrolysis of sialic acids	48
<b>FucosEXO™</b> Hydrolysis of $\alpha$ 1-2,3,4 fucose	50
<b>GalactEXO™</b> Hydrolysis of $\beta$ 1-3,4-linked galactose	52
<b>GalNAcEXO™</b> Hydrolysis of $\alpha$ -linked GalNAcs	54
<b>OmniGLYZOR™</b> Hydrolysis of $\alpha$ -linked GalNAcs	56

## Affinity Purification

Tool for affinity purification enriches proteins and peptides carrying O-glycans from complex samples.

<b>GlycOCATCH®</b> Enrichment of O-glycoproteins	58
---	----

## Research Antibodies

Antibodies to detect SmartEnzymes in various assays.

<b>Anti-FabRICATOR®</b> Detection of the FabRICATOR enzyme	60
---	----



**FabRICATOR®**  
Below hinge  
digestion of IgG



**FabALACTICA®**  
Above hinge  
digestion of human IgG1



**FabDELLO™**  
Above hinge  
digestion of human IgG1

ANTIBODY DIGESTION



**PNGase F**  
Hydrolysis of  
N-glycans



**OpeRATOR®**  
O-glycan-specific  
protein digestion



**OglyZOR®**  
Hydrolysis of  
core-1 O-glycans



**SialEXO®**  
Hydrolysis of  
sialic acids

GLYCAN PROFILING



**GlycINATOR®**  
Hydrolysis of  
IgG Fc glycans



**IgGZERO®**  
Hydrolysis of  
complex Fc glycans



**GlyCLICK®**  
Site-specific  
conjugation of native IgG



**TransGLYCIT™**  
Transglycosylation  
of IgG

ANTIBODY DEGLYCOSYLATION

ANTIBODY CONJUGATION





**FabRICATOR® Z**  
Below hinge  
digestion of mouse IgG



**FabULOUS™**  
Above hinge  
digestion of IgG



**GingisKHAN®**  
Above hinge  
digestion of human IgG1



**FucoseEXO™**  
Hydrolysis of  
 $\alpha$ 1-2,3,4 fucose



**GalactEXO™**  
Hydrolysis of  $\beta$ 1-3,4-  
linked galactose



**GalNAcEXO™**  
Hydrolysis of  
 $\alpha$ -linked GalNAcs



**GlycOCATCH®**  
Enrichment of  
O-glycoproteins

— AFFINITY PURIFICATION —



**GingisREX®**  
Arginine-specific  
protein digestion



**Anti-FabRICATOR®**  
Detection of the  
FabRICATOR enzyme


— PROTEOMICS —

— RESEARCH ANTIBODIES —



# FabRICATOR®

## Below Hinge Digestion of IgG

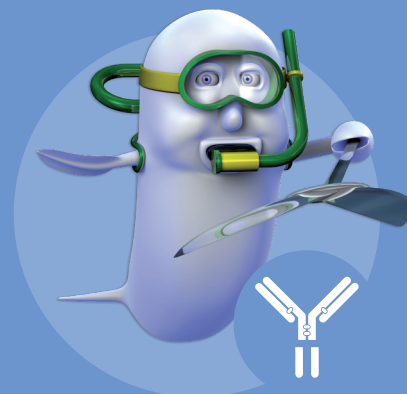

**FabRICATOR (IdeS) is an IgG-specific cysteine protease that digests antibodies at a single amino acid site below the hinge, generating homogenous F(ab')<sub>2</sub> and Fc fragments. The enzyme is widely used in characterization, quality control, stability testing, production monitoring and clone selection of antibody-based therapeutics such as mAbs, ADCs, biosimilars and Fc-fusion proteins.**

- 


Human IgG1-4, IgG of some classes from monkey, rat, rabbit, sheep
- 

CPAPELLG / GPSVF
- 

30 min reaction
- 

No reducing agents or co-factors needed

### Antibody Subunit Workflow

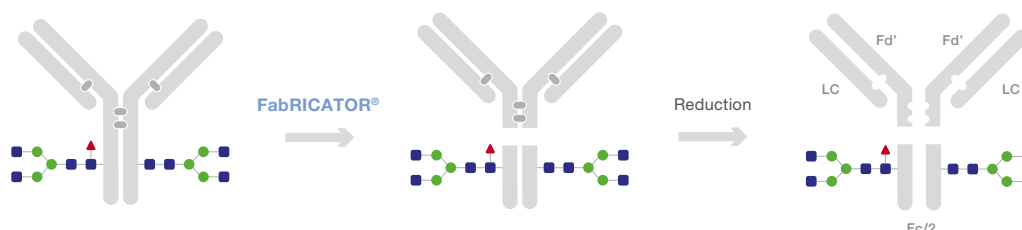

The FabRICATOR enzyme digests antibodies in 30 minutes at 37°C to generate intact F(ab')<sub>2</sub> and Fc subunits without the need for optimization and no risk of

over-digestion, otherwise common using proteases such as papain, ficin or pepsin. For sample preparation in middle-level LC-MS workflows, the subunits can be

reduced to generate 23-25 kDa fragments allowing increased mass resolution for accurate determination of critical quality attributes.

### Product Formats



**FabRICATOR®**  
Lyophilized enzyme for digestion of IgG



**FragIT™**  
Immobilized enzyme for digestion of IgG in spin columns



**FragIT™ kit**  
Immobilized enzyme and CaptureSelect™ Fc for purification of IgG subunits



**FabRICATOR®-HPLC**  
Immobilized enzyme on HPLC-compatible resin for on-column digestion



**FabRICATOR® MagIC**  
Immobilized enzyme on magnetic beads for automated workflows



**FabRICATOR® Validation kit**  
Three different batches for validation of FabRICATOR-based analytical methods

Thermo Scientific™ CaptureSelect™ resin from Thermo Fisher Scientific Inc. and its subsidiaries.  
Thermo Scientific and CaptureSelect are trademarks of Thermo Fisher Scientific Inc. and its subsidiaries.




## Determining Antibody Quality Attributes using FabRICATOR®

The homogenous fragments generated with FabRICATOR can be mono-isotopically resolved for in-depth mapping of attributes that may be unresolved at the intact level. This middle-level workflow can be used for amino acid sequence verification and oxidation analysis of quality attributes, including deamidation, pyroglutamination and glycan profiling. The relative abundance of eleven Fc

glycans were quantified on a therapeutic antibody after FabRICATOR digestion for increased mass resolution, and no additional digestion or labeling were needed for similarity assessment.


Digestion of antibodies into homogenous subunits allows for middle-level analysis for rapid monitoring of changes that may occur during storage and formulation, such as methionine oxidation.

A therapeutic antibody subjected to forced oxidation was site-specifically digested with FabRICATOR, deglycosylated by GlycINATOR and analyzed at the subunit level. Domain specific data was obtained in the analysis of oxidation events as shown in the deconvoluted mass spectra for the Fc/2 fragment. The analysis was performed in a high throughput format.

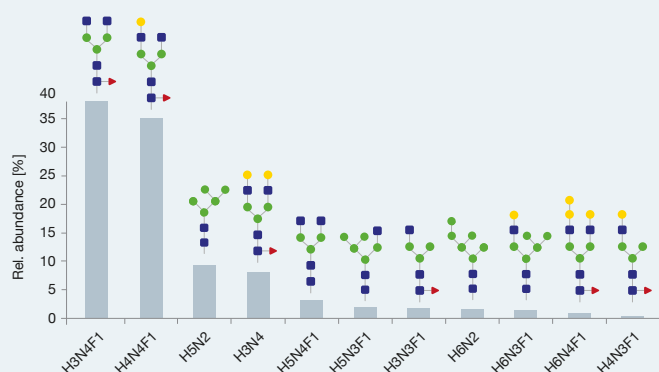

**Determining antibody quality attributes.** Relative abundance of glycans on cetuximab Fc/2 fragments, generated with FabRICATOR and analyzed with LC-MS at the subunit level by Ayoub and colleagues (Ayoub, 2013).

### Adalimumab (IgG1)

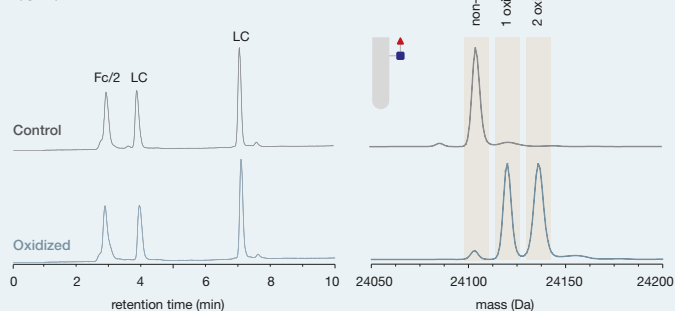


**Quality control of antibody oxidation.** Oxidation analysis of adalimumab showing UV (280 nm) chromatograms before (control) and after (oxidized) incubation with 0.25% H<sub>2</sub>O<sub>2</sub>. A deconvoluted mass spectrum of the deglycosylated Fc/2 fragments demonstrates the degree of oxidation.

## Preparative Generation of F(ab')<sub>2</sub> and Fc/2 Subunits with FragIT™ Kit



For applications that require fast and convenient digestion of antibodies or Fc-fusion proteins without residual enzyme in the preparation, FragIT spin columns that contain FabRICATOR immobilized on agarose beads can be used. With FragIT kit, the homogenous F(ab')<sub>2</sub> and Fc fragments can also be separated using CaptureSelect™ affinity purification columns.

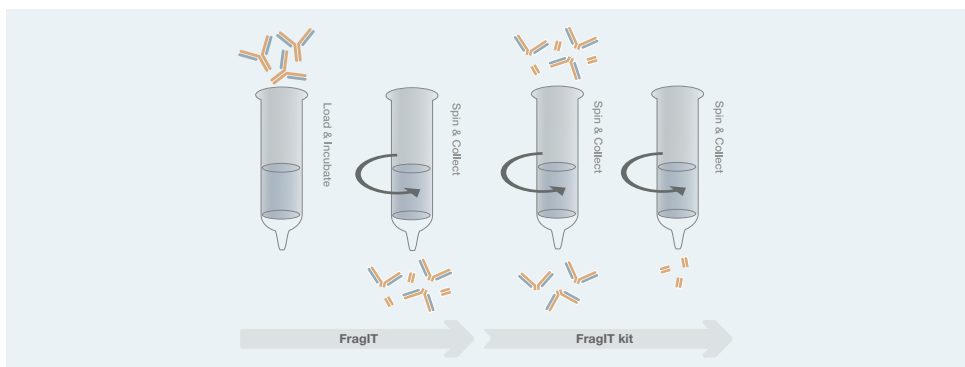

### References

Ayoub, D. et al., 2013. *Correct primary structure assessment and extensive glyco-profiling of cetuximab by a combination of intact, middle-up, middle-down and bottom-up ESI and MALDI mass spectrometry techniques.* mAbs, 5(5), pp.699–710.




## On-column Digestion for Middle-Level Workflows with FabRICATOR®-HPLC

FabRICATOR-HPLC contains the FabRICATOR enzyme immobilized on an HPLC-compatible resin for consistent and reproducible on-column antibody digestion. The FabRICATOR-HPLC col-

umn can be used in a standard LC-MS routine analysis, but more advanced configurations are possible, for example with a 2D-LC setup. Ultimately, a bioreactor can be connected directly to the

MS in an automated online middle-level workflow. This would significantly reduce operator time, sample handling errors and increase throughput.

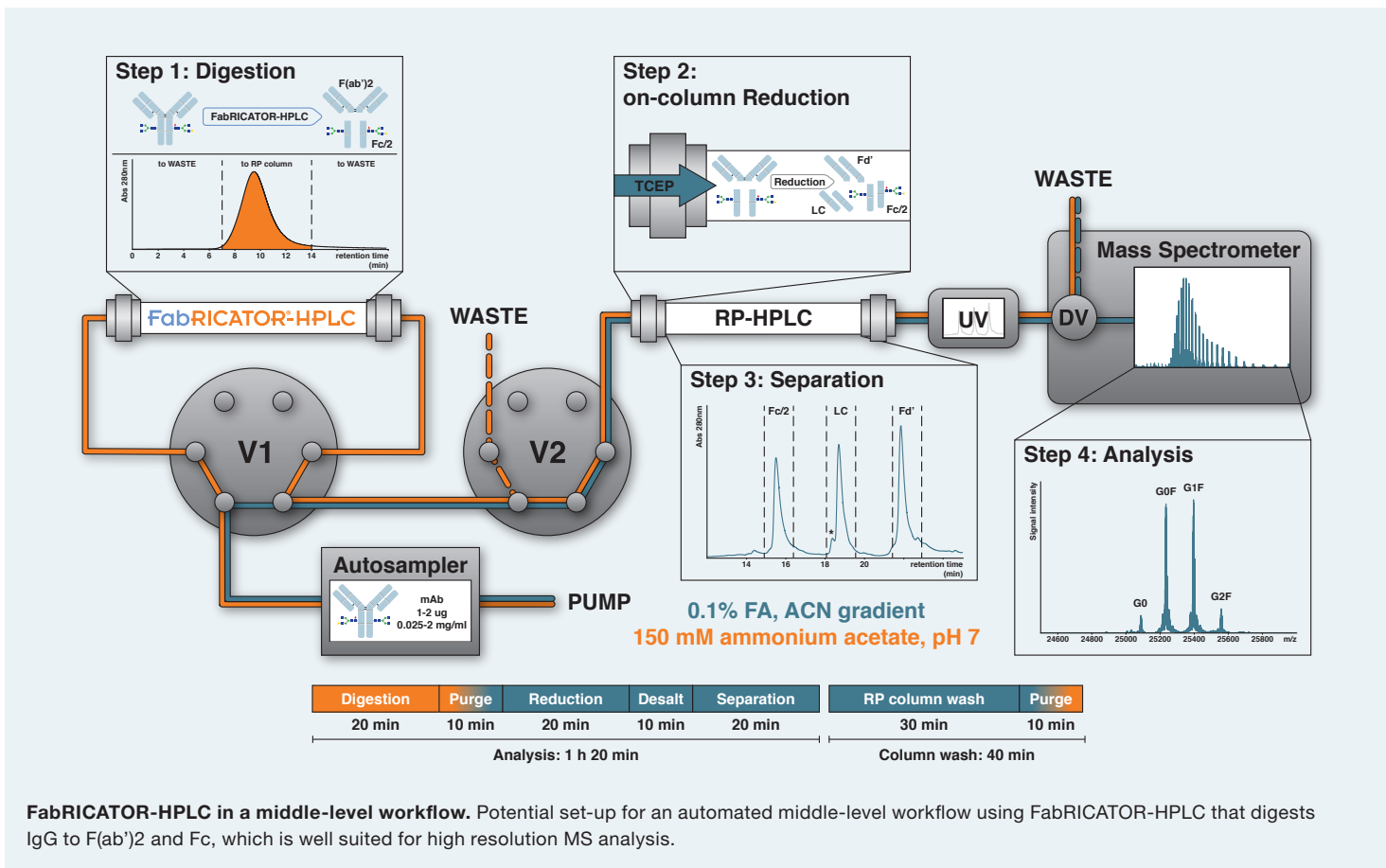




## Robust On-column Performance

The on-column digestion efficiency was maintained by FabRICATOR-HPLC during continuous operation at 37°C for up to 14 days according to tests, where 2 µg trastuzumab was injected every 4 h and analyzed by LC-MS. The generated antibody fragments were reduced on-column and analyzed using an automated workflow with more than 95 % antibody digestion yield during the 14-day period.

### Performance testing of FabRICATOR-HPLC.

Comparison between digestion with standard in-solution FabRICATOR protocol (orange) and digestion using the automated FabRICATOR-HPLC workflow (blue), and quantification of the column's performance over a period of 14 days (bottom). The asterisk marks LC fragments that are not completely reduced with one intramolecular disulfide bridge intact.

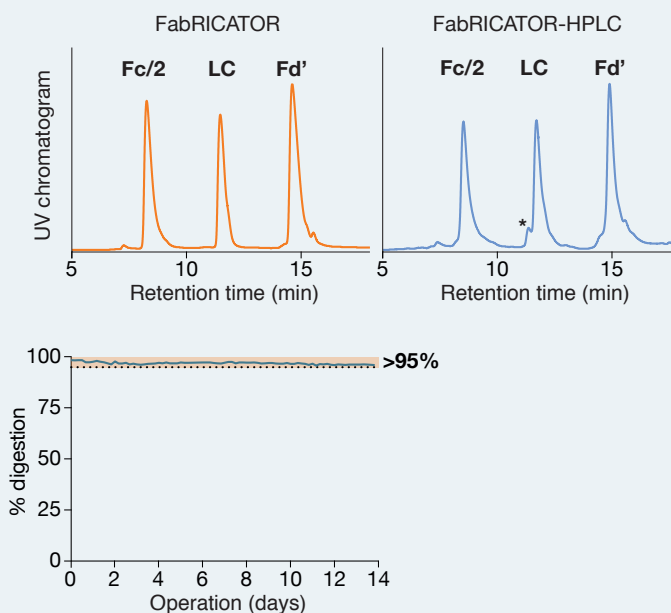


## Automated Antibody Fc Glycan Analysis

The performance and operational stability of the FabRICATOR-HPLC column are demonstrated here using analysis of trastuzumab Fc glycosylation. Automated middle-level analysis using FabRICATOR-HPLC yielded mass spectra virtually indistinguishable from the one obtained by a standard in-solution digestion workflow. The analytical results of the Fc glycosylation profiles were stable and reproducible during 14 days of continuous operation with standard deviations of less than 0.5 % for all glycoforms.

### Glycan analysis using FabRICATOR-HPLC.

Deconvoluted mass spectra of trastuzumab Fc/2 from in-solution FabRICATOR digestion (top, orange), FabRICATOR-HPLC (top, blue), and glycosylation profiles of trastuzumab (bottom, orange, n=10) or (bottom, blue, n=28).

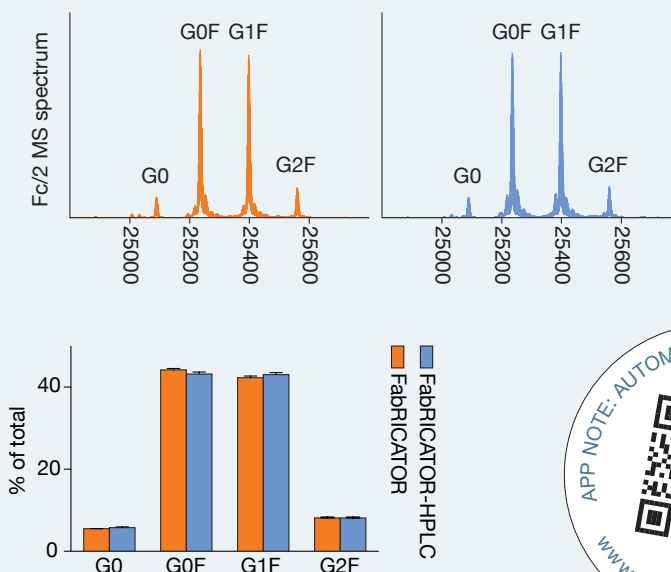

## FabRICATOR®-HPLC

Column hardware: PEEK/biocompatible

Column dimensions: 2.1 mmD × 50 mL

Support resin: POROS®\*

Typical flow rate: 0.025-0.05 mL/min

Maximum Pressure: 100 bar

Operating pH: 6.5 - 8.0

Operating temperature: 37°C

Storage conditions: +4-8°C  
(Do not freeze!)

Number of days of continuous operation: >10\*\*

Injections per column: >200\*\*

Start material: Human IgG1-4,  
Fc-fusion proteins

\* See legal and disclaimers. \*\* Depending on specific application




## Parallel Subunit Generation of Antibodies for Automated Middle-level Workflows with FabRICATOR® MagIC



For automated processing of antibodies, FabRICATOR MagIC contains the FabRICATOR enzyme immobilized on

magnetic agarose beads. The magnetic format enables parallel sample preparation and analysis of multiple antibodies

in middle-level workflows. It can be used in both manual or automated setups for routine analysis of IgG-based biologics.

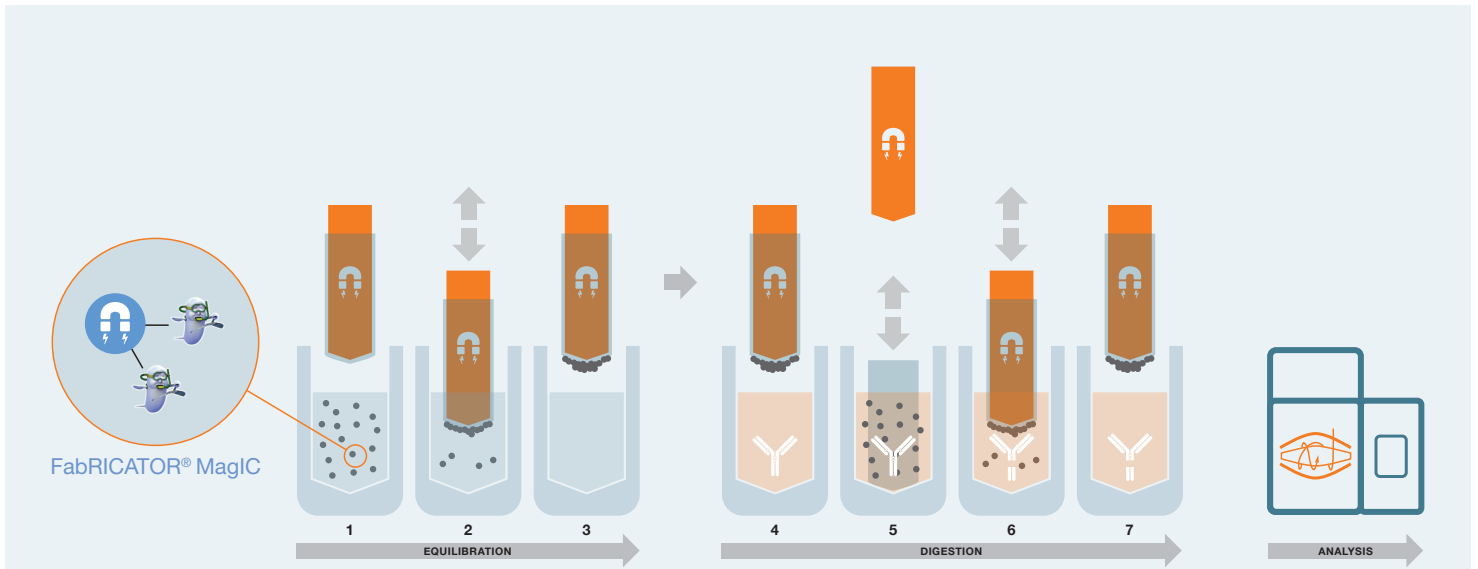


**FabRICATOR MagIC workflow on Thermo Scientific KingFisher™ Purification System.** IgG, FabRICATOR MagIC beads and PBS are dispensed into a 96-deepwell plate according to the KingFisher workflow. The beads are (1) transferred to wells containing PBS for (2) equilibration and (3) collection. The beads are (4) added to the samples and (5) incubated with mixing for 10-20 min at 37°C, prior to (6) collection and the pure F(ab')<sub>2</sub> and Fc are left in the sample wells (7). The trademark KingFisher is the property of Thermo Fisher Scientific.

### Automated Middle-Level Workflow

FabRICATOR MagIC generates intact F(ab')<sub>2</sub> and Fc subunits and the disulfide bonds can easily be reduced using 5 mM TCEP in a one-pot reaction.

The generated Fc/2, LC and Fd' fragments allow detailed analysis and monitoring using middle-level workflow, and results in a better resolution compared to intact level approaches.

Reversed-phase LC-MS analysis of the resulting fragments from an antibody processed with FabRICATOR MagIC demonstrates the complete digestion.

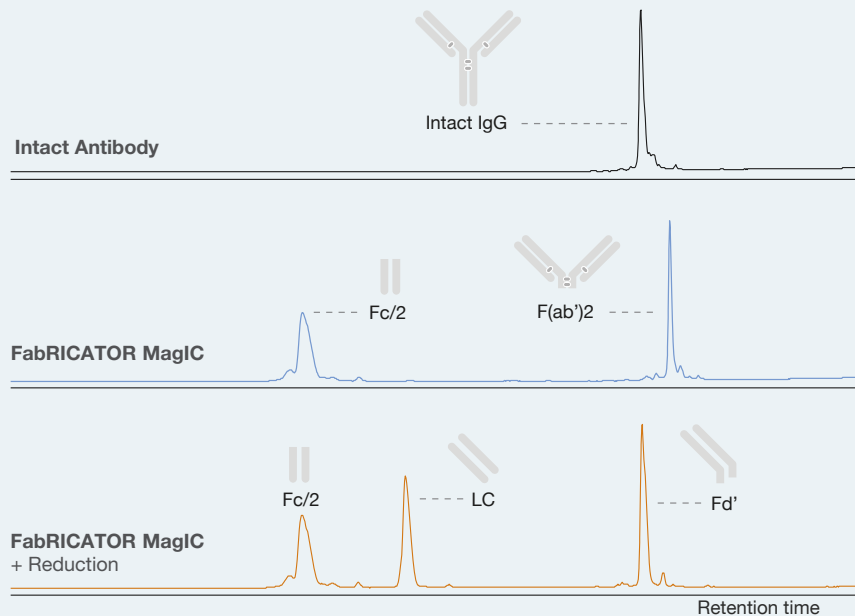

**Generation of antibody subunits with FabRICATOR MagIC.** Analysis of a monoclonal antibody by RPLC-MS at the intact level (top), after FabRICATOR MagIC digestion (middle) and after treatment with FabRICATOR MagIC and 5 mM TCEP (bottom). Data was obtained in collaboration with Symphogen A/S (Ballerup, Denmark).


### Fast Automated Monitoring of CQAs

Using FabRICATOR MagIC for automated digestion of multiple mAbs, different CQAs such as oxidation can be analyzed. A mix of six different mAbs subjected to


forced degradation was analyzed by LC-MS at the intact or the subunit level. The results show the ability to quickly detect and monitor such modifications in a complex mixture of mAbs using

FabRICATOR MagIC. A detailed Fc-glycan profiling was also performed on a mAb by middle-level analysis after FabRICATOR MagIC digestion and reduction with 5 mM TCEP.


#### Intact Antibodies

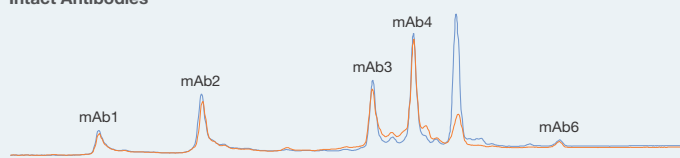


#### FabRICATOR MagIC + Reduction

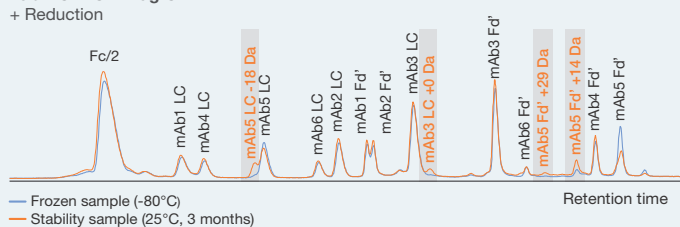

**Oxidation studies using FabRICATOR MagIC.** Six mAbs were analyzed for tryptophan oxidation during forced degradation studies by RPLC-MS at the intact level (top panel) and after FabRICATOR MagIC with 5 mM TCEP (bottom panel). New peaks detectable by the subunit analysis are highlighted. Data was obtained in collaboration with Symphogen A/S (Ballerup, Denmark).

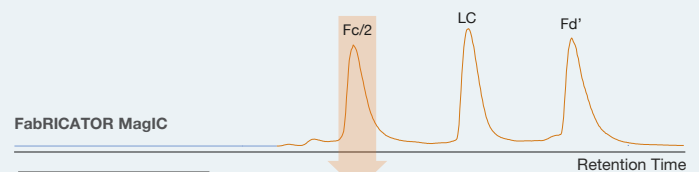

Glycoform	Abundance
Man5	2.5%
G0	4.7%
G0F	47.6%
G1F	38.4%
G2F	6.8%

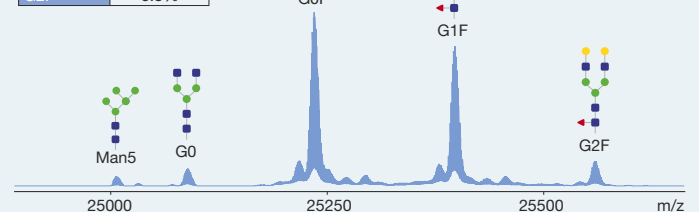

**Glycan analysis using FabRICATOR MagIC.** Analysis of Fc glycosylation by middle-level LC-MS using FabRICATOR MagIC (top) and quantification of glycoforms from a deconvoluted mass spectrum of the Fc/2 fragment (bottom).



**FabRICATOR®**

Lyophilized enzyme for digestion of IgG.



PRODUCT	DESCRIPTION	ID
FabRICATOR, 2000 units	Digests 2 mg IgG	A0-FR1-020
FabRICATOR, 5000 units	Digests 5 mg IgG	A0-FR1-050
FabRICATOR, 5 x 5000 units	Digests 5 x 5 mg IgG	A0-FR1-250
FabRICATOR, 96 x 100 units	Digests 96 x 100 µg IgG	A0-FR1-096
FabRICATOR, 8 x 100 units	Digests 8 x 100 µg IgG	A0-FR1-008
FabRICATOR LE (low endotoxin), 2000 units	Digests 2 mg IgG	A0-FR8-020
FabRICATOR LE (low endotoxin), 5000 units	Digests 5 mg IgG	A0-FR8-050

**FabRICATOR® Validation kit**

Three different batches of the FabRICATOR enzyme for validation of FabRICATOR-based analytical methods.



PRODUCT	DESCRIPTION	ID
FabRICATOR, 3 x 2000 units	Digests 3 x 2 mg IgG	A0-FR4-060

**FragIT™**

Immobilized FabRICATOR enzyme for digestion of IgG in spin columns.



PRODUCT	DESCRIPTION	ID
FragIT, Microspin 2 x 0.5 mg	Digests 2 x 0.5 mg IgG	A0-FR6-010
FragIT, Microspin 5 x 0.5 mg	Digests 5 x 0.5 mg IgG	A0-FR6-025
FragIT, Microspin 10 x 0.5 mg	Digests 10 x 0.5 mg IgG	A0-FR6-050
FragIT, Midispin 1-10 mg	Digests 1-10 mg IgG	A0-FR6-100
FragIT, Maxispin 10-100 mg	Digests 10-100 mg IgG	A0-FR6-1000





**FragIT™ kit**

Immobilized FabRICATOR enzyme and CaptureSelect™ Fc affinity spin columns for digestion and purification of IgG subunits.



PRODUCT	DESCRIPTION	ID
FragIT kit, Microspin 0.5 mg	Digests and purifies 0.5 mg IgG	A2-FR2-005
FragIT kit, Microspin 5 × 0.5 mg	Digests and purifies 5 × 0.5 mg IgG	A2-FR2-025
FragIT kit, Microspin 10 × 0.5 mg	Digests and purifies 10 × 0.5 mg IgG	A2-FR2-050
FragIT kit, Midispin 10 mg	Digests and purifies 10 mg IgG	A2-FR2-100
FragIT kit, Maxispin 100 mg	Digests and purifies 100 mg IgG	A2-FR2-1000

**FabRICATOR®-HPLC**

Immobilized FabRICATOR enzyme on HPLC-compatible resin for on-column digestion.



PRODUCT	DESCRIPTION	ID
FabRICATOR-HPLC	Number of injections per column: >200	A0-FRC-050

**FabRICATOR® MagIC**

Immobilized FabRICATOR enzyme on magnetic beads for automated workflows.




PRODUCT	DESCRIPTION	ID
FabRICATOR MagIC, 2 ml	Digests up to 24 × 200 µg IgG	A0-FRM-024
FabRICATOR MagIC, 4 × 2 ml	Digests up to 96 × 200 µg IgG	A0-FRM-096



# FabALACTICA®


## Above Hinge Digestion of Human IgG1

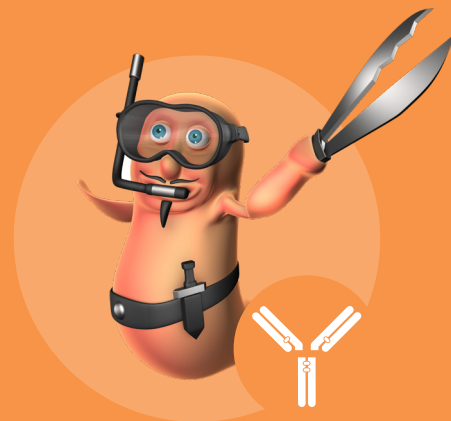

**FabALACTICA (IgdE)** is a cysteine protease that digests human IgG1 at a specific site above the hinge, generating intact and homogenous Fab and Fc fragments. The enzyme facilitates the characterization of bi- and multispecific antibodies and analyses of intact and paired Fc glycosylation. It also simplifies the study of monovalent binding, higher order structures, disulfide scrambling and antibodies with mutated hinge regions.

-  Human IgG1
-  KSCDKT / HTCPPC
-  O/N reaction
-  No reducing agents or co-factors needed

### Antibody Subunit Workflow

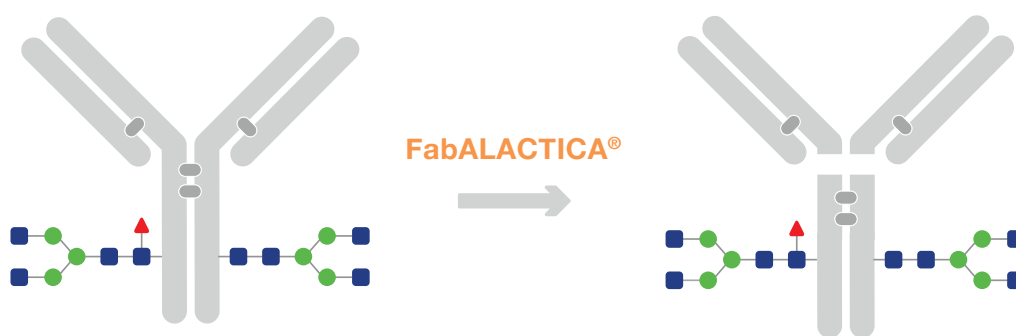

The use of proteases with high specificity for IgG has allowed for subunit profiling of antibody-based therapeutics and studies of key quality attributes using LC-MS.

The FabALACTICA enzyme specifically digests human IgG1 at a single site above the hinge. The generated fragments can be used for a broad range of applications,

including antibody characterization and NMR studies for higher order structures.

### Product Formats



**FabALACTICA®**  
Lyophilized enzyme for digestion of human IgG1



**Immobilized FabALACTICA®**  
Immobilized enzyme for digestion of human IgG1 in spin columns



**FabALACTICA® Fab kit**  
Immobilized enzyme and CaptureSelect™ Fc for purification of human IgG1 subunits

Thermo Scientific™ CaptureSelect™ resin from Thermo Fisher Scientific Inc. and its subsidiaries.  
Thermo Scientific and CaptureSelect are trademarks of Thermo Fisher Scientific Inc. and its subsidiaries.



## Comparison of Enzymatic Strategies for Obtaining Fab Fragments



The SmartEnzymes FabALACTICA and GingisKHAN digests IgG at a single digestion site and generate homogenous Fab and Fc fragments from human IgG1. Traditional enzymes such as papain and Lys-C can also be used to generate Fab fragments, but the unspecific enzymatic activities need careful optimization to minimize over-digestion.

Enzyme	FabALACTICA	GingisKHAN	Papain	Lys-C
Digestion site				
Specificity	IgG-specific/ One digestion site	One digestion site	Unspecific/ Several digestion sites	Unspecific/ Several digestion sites
Selectivity	Human IgG1	Human IgG1	Several proteins	Several proteins
Reducing conditions	No	Yes, 2 mM cysteine	Yes	No
Reaction time	O/N (16-18 h)	1 h	1-24 h	1-24 h


## Generation of Intact Fab and Fc Fragments



In order to demonstrate the activity of FabALACTICA, the three therapeutic monoclonal antibodies cetuximab,

trastuzumab and adalimumab were digested and the resulting fragments were analyzed using RP-HPLC.

The defined peaks display the homogenous fragment generation and the robust enzymatic performance.

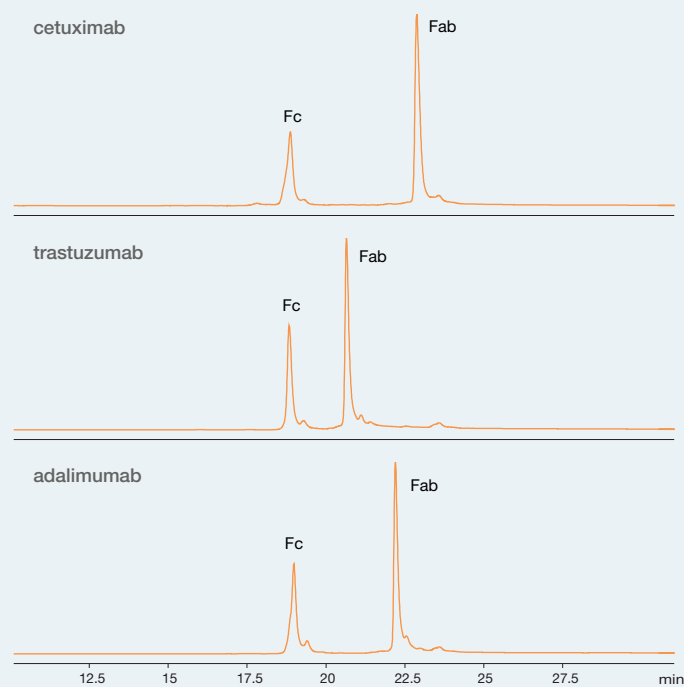

Separation of intact Fab and Fc fragments by RP-HPLC after FabALACTICA digestion. The digestion was performed O/N at 37°C.


## Paired Glycan Analysis using FabALACTICA® and LC-MS



Fc glycosylation is a critical quality attribute of monoclonal antibodies that can affect their stability, solubility, pharmacodynamic-, and pharmacokinetic properties. This post-translational modification

requires close monitoring during both development and manufacturing of a therapeutic antibody. Here, trastuzumab was digested using FabALACTICA, and the resulting intact Fc fragment (~53 kDa)

enabled simultaneous characterization of the two conserved Fc glycosylation sites and their pairing with high mass accuracy.

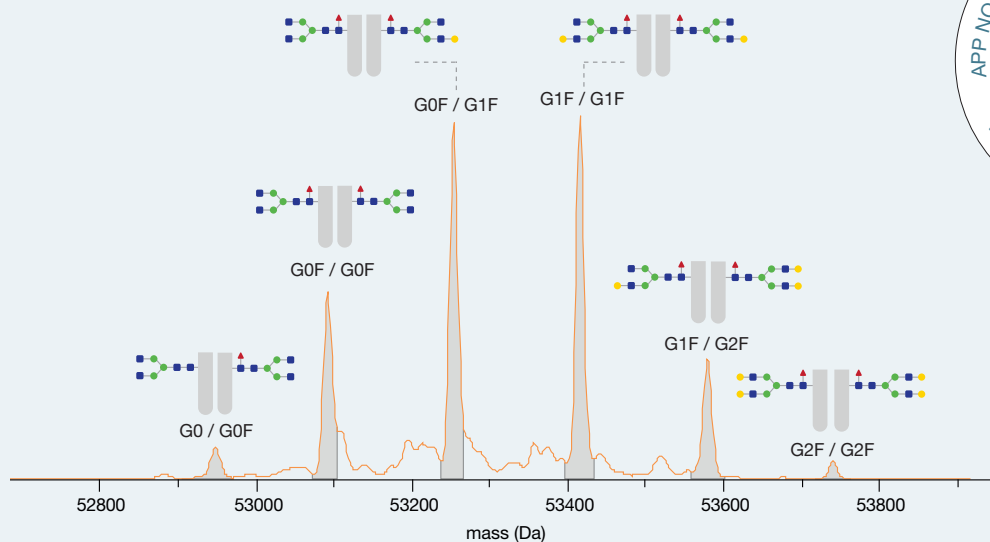

**Paired glycan analysis of trastuzumab Fc fragments.** Trastuzumab was digested using FabALACTICA O/N at 37°C, and intact Fc and Fab fragments were studied using LC-MS.




## 2D-NMR Spectroscopy Analysis using Immobilized FabALACTICA®



Critical quality attributes such as aggregation and unwanted inherent heterogeneity can be assessed by measuring the higher order structure (HOS) of the protein. Two-dimensional nuclear magnetic resonance spectroscopy

(2D-NMR) allows for the precise atomic-level comparison of HOS for IgG-based biopharmaceuticals. Here, infliximab was first treated with H<sub>2</sub>O<sub>2</sub> to induce oxidation before it was digested using Immobilized FabALACTICA.

The NMR spectra of the oxidized sample displayed an altered local environment for the methionines due to oxidation in both Fc and Fab fragments, while the higher order structure remained largely unaffected.

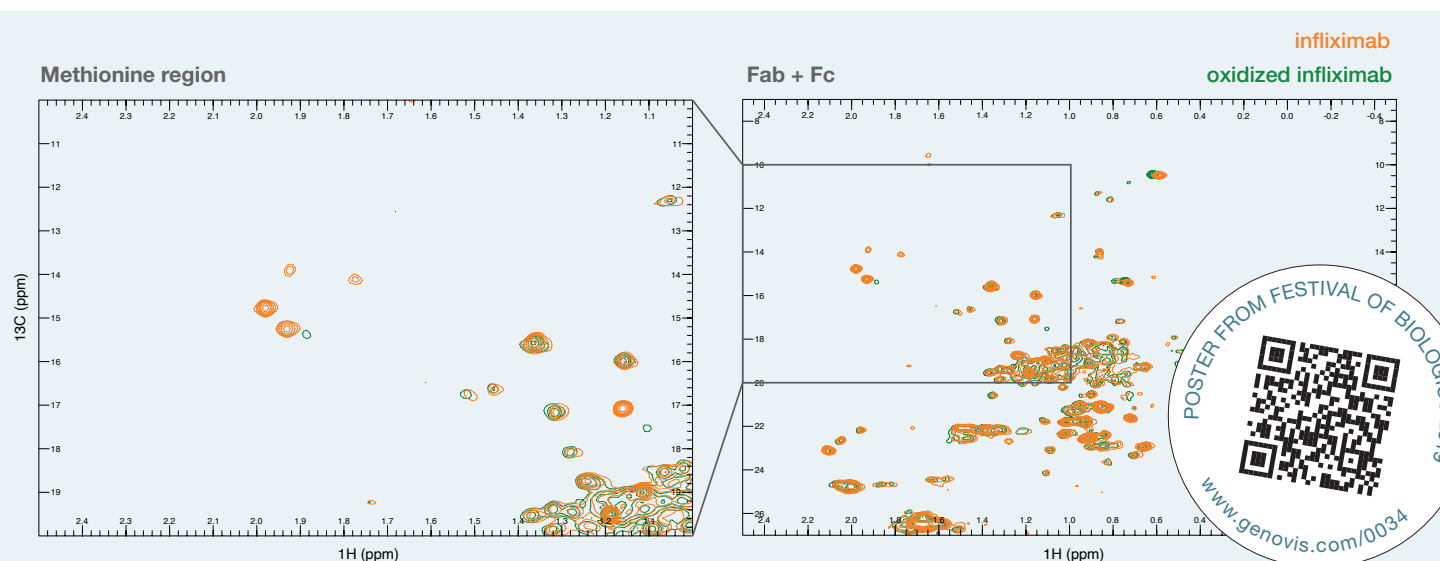

**2D-NMR analysis following FabALACTICA digestion to evaluate the impact of oxidation on infliximab.** 1H-13C methyl HSQC of the combined Fab and Fc fragments of infliximab. Originator (orange) and oxidized originator sample (green). The spectra were acquired at 310 K using an 800 MHz Bruker Avance at the Swedish NMR Center. Data was obtained in collaboration with SARomics Biostructures (Lund, Sweden).






## Preparation of Pure Fabs using FabALACTICA® Fab Kit



The FabALACTICA Fab kit consists of spin columns of Immobilized FabALACTICA for antibody digestion and spin columns of CaptureSelect™ Fc resin for affinity binding of the Fc fragments. Here, the

therapeutic antibodies adalimumab, trastuzumab and cetuximab were incubated with Immobilized FabALACTICA before the Fc fragments were captured in the CaptureSelect™ Fc spin columns.

The Fabs were easily collected by centrifugation, and the resulting Fab preparation was homogenous and pure.

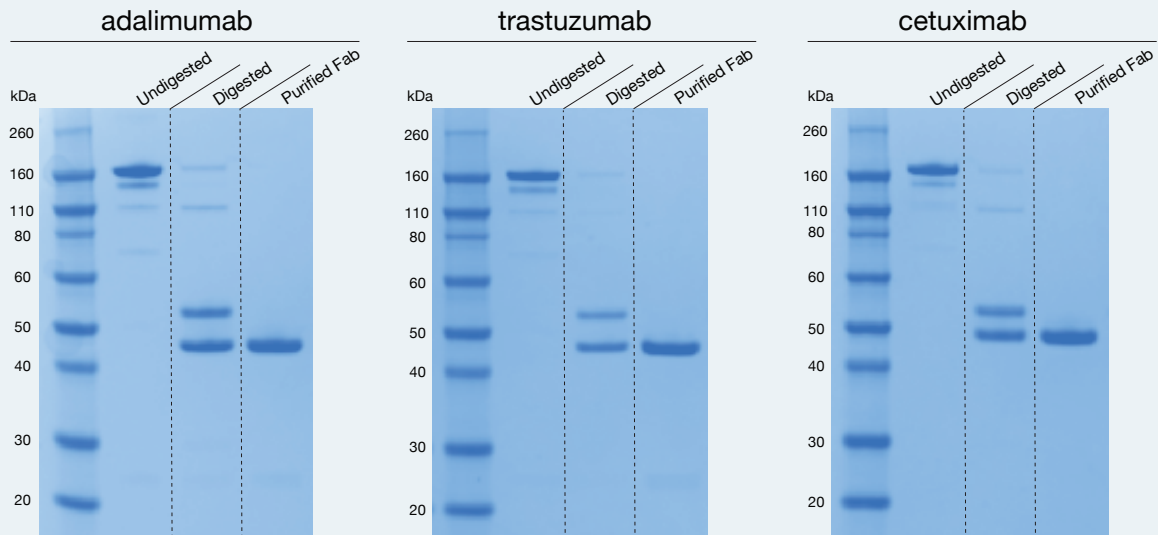

**Preparation of pure Fabs from three monoclonal antibodies using the FabALACTICA Fab kit.** Adalimumab, trastuzumab and cetuximab digested by Immobilized FabALACTICA. Pure Fab fragments were obtained in a high yield from all three mAbs using the FabALACTICA Fab kit.

### FabALACTICA®

Lyophilized enzyme for digestion of human IgG1.



PRODUCT	DESCRIPTION	ID
FabALACTICA, 2000 units	Digests 2 mg human IgG1	A0-AG1-020

### Immobilized FabALACTICA®

Immobilized enzyme for digestion of human IgG1 in spin columns.



PRODUCT	DESCRIPTION	ID
Immobilized FabALACTICA, Microspin 2 × 0.5 mg	Digests 2 × 0.5 mg human IgG1	A0-AG6-010
Immobilized FabALACTICA, Microspin 10 × 0.5 mg	Digests 10 × 0.5 mg human IgG1	A0-AG6-050
Immobilized FabALACTICA, Midispin 5-10 mg	Digests 5-10 mg human IgG1	A0-AG6-100
Immobilized FabALACTICA, Maxispin 10-100 mg	Digests 10-100 mg human IgG1	A0-AG6-1000

### FabALACTICA® Fab kit

Immobilized enzyme and CaptureSelect™ Fc affinity spin columns for generation and purification of human IgG1 subunits.




PRODUCT	DESCRIPTION	ID
FabALACTICA Fab kit, Microspin 0.5 mg	Digests and purifies 0.5 mg human IgG1	A2-AFK-005
FabALACTICA Fab kit, Microspin 5 × 0.5 mg	Digests and purifies 5 × 0.5 mg human IgG1	A2-AFK-025
FabALACTICA Fab kit, Midispin 10 mg	Digests and purifies 10 mg human IgG1	A2-AFK-100
FabALACTICA Fab kit, Maxispin 100 mg	Digests and purifies 100 mg human IgG1	A2-AFK-1000

# FabDELLO™


## Above Hinge Digestion of Human IgG1



**FabDELLO is a protease that specifically digests human IgG1 at a single site above the hinge, generating intact Fab and Fc fragments within two hours with no need for reducing conditions. The enzyme is active on antibodies with mutated hinge regions, such as the LALA mutation, and enables middle-level LC-MS analysis for characterization of critical quality attributes on antibodies with mutated hinges.**

-  Human IgG1
-  KSCDK / THTCPPCP
-  2 hour reaction
-  Requires presence of calcium ions. 200 mM CaCl<sub>2</sub> is included.

### Antibody Subunit Workflow

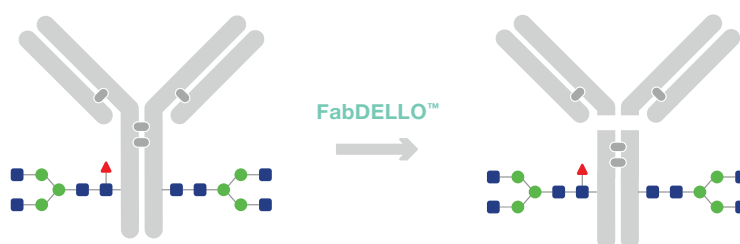

FabDELLO digests human IgG1 at a single site above the hinge (...KSCDK / THTCPPCP...), generating intact Fab and

Fc fragments. A second digestion site may appear on the Fc if the N-glycans are removed. FabDELLO is active under

native conditions and requires the presence of calcium ions. Optimal activity is obtained at 37°C and pH 7-8.5.

### Product Formats



#### FabDELLO™

Lyophilized enzyme for digestion of human IgG1

### FabDELLO

PRODUCT	DESCRIPTION	ID
FabDELLO, 8 × 100 units	Digests 8 × 100 µg human IgG1	B1-BD1-008
FabDELLO, 96 × 100 units	Digests 96 × 100 µg human IgG1	B1-BD1-096






## Middle-level LC-MS Analysis of Mutated Antibodies

Mutated hinge regions are common modifications of therapeutic antibody candidates. The hinge mutations may affect the antibody digestion efficiency of enzymes with high substrate specificity. FabDELLO acts specifically at an exposed lysine residue above the hinge and enables middle-level LC-MS analysis of antibodies with mutated hinge regions. Here, the commercially available IgG1 risankizumab, with a LALA mutation in the lower hinge region, was digested with FabDELLO and reduced before middle-level LC-MS analysis. The defined peaks observed after digestion with FabDELLO show homogenous fragment generation and robust performance of the enzyme.

**Middle-level LC-MS analysis of risankizumab after digestion with FabDELLO.** The digested material was subjected to 20mM DTT for 1 h at 37°C to reduce the cysteine bridges. The sample was analyzed by reversed-phase LC-MS on a Waters™ BioAccord™ system equipped with a Waters™ BioResolve™ RP mAb column (2.1 × 50mm). UV 280 chromatogram showing Fc/2, LC and Fd, and deconvoluted Fc/2, LC and Fd are shown.

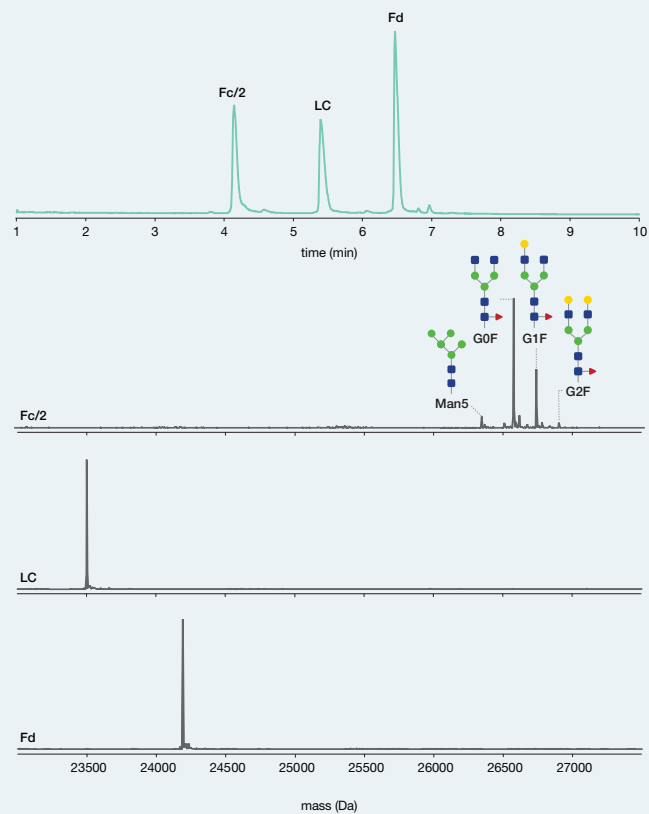


## Complete Digestion of Mutated Antibodies and Fc-fusion Proteins



FabDELLO digests human IgG1 at a single site above the hinge and generates intact Fab and Fc fragments. To demonstrate the specific activity of FabDELLO, a selection of IgG1 and Fc-fusion proteins

were digested and analyzed by non-reduced SDS-PAGE. A complete digestion was achieved for all substrates, including IgG1 with the otherwise difficult LALA and LFLE mutations. The precise and

effective performance of the FabDELLO enzyme on this wide range of substrates validates its use for the development of biopharmaceuticals.

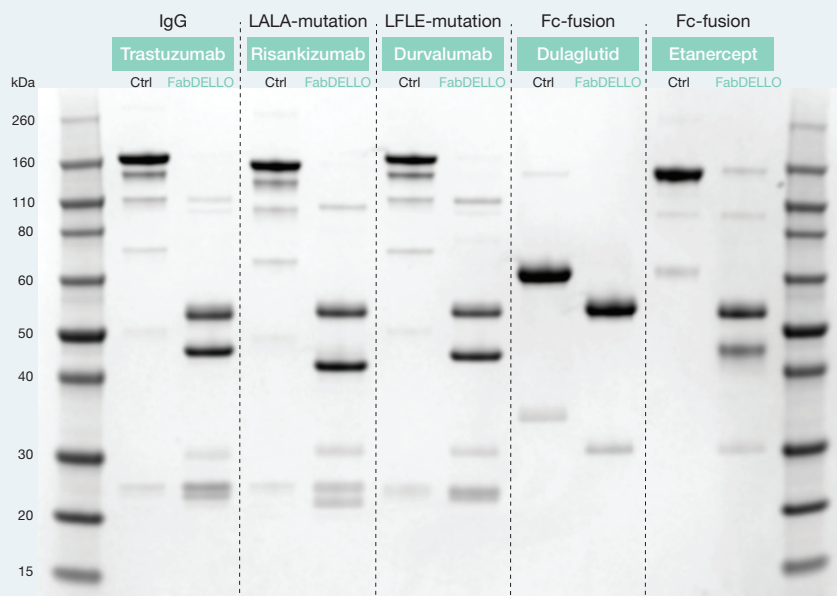

SDS-PAGE analysis under non-reducing conditions of commercial antibodies and Fc-fusion proteins digested with FabDELLO.

# FabRICATOR® Z


## Below Hinge Digestion of Mouse IgG



**FabRICATOR Z (IdeZ)** is a cysteine protease that digests mouse IgG2a and IgG3 at a specific site below the hinge, generating a homogenous pool of F(ab')<sub>2</sub> and Fc fragments after a two-hour incubation.

 Mouse IgG2a and IgG3, human IgG1-4, IgG of some classes from monkey, rabbit and sheep
  CPAPNLLG / GPSVF
  2h reaction
  No need for reducing agents or co-factors

### Antibody Digestion Workflow

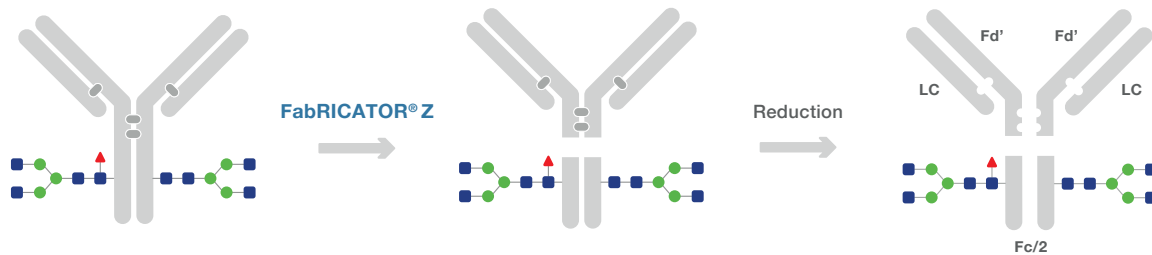

FabRICATOR Z digests mouse IgG2a and IgG3 at a specific site below the hinge region, generating intact F(ab')<sub>2</sub>

and Fc fragments. Some mouse IgG2a that FabRICATOR fails to digest, are readily digested by FabRICATOR Z, but

longer incubation times may be required. There is no risk of over-digestion because of the high specificity of the enzyme.

### Product Formats



#### FragIT™ Z

Immobilized enzyme for digestion of mouse IgG2a and IgG3 in spin columns



#### FabRICATOR® Z

Lyophilized enzyme for digestion of mouse IgG2a and IgG3



#### FragIT™ Z kit

Immobilized enzyme and CaptureSelect™ Fc for purification of IgG subunits

Thermo Scientific™ CaptureSelect™ resin from Thermo Fisher Scientific Inc. and its subsidiaries.  
Thermo Scientific and CaptureSelect are trademarks of Thermo Fisher Scientific Inc. and its subsidiaries.






## Digestion of Mouse IgG2a using FabRICATOR® Z

FabRICATOR Z (IdeZ) and FabRICATOR (IdeS) were here used to digest mouse IgG2a. After 2 hours of incubation, FabRICATOR Z (IdeZ) was found to readily digest mouse IgG2a, whereas FabRICATOR only digested a small amount of the antibody.

### Digestion of mouse IgG2a using FabRICATOR Z and FabRICATOR.

The F(ab')<sub>2</sub> fragment is detected at ~110 kDa and the Fc fragments at ~30 kDa. The enzyme is detected at ~37 kDa.

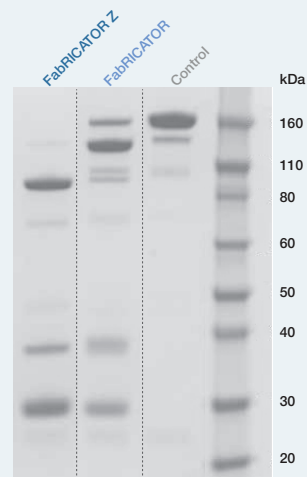


## Preparative Generation of Mouse IgG2a and IgG3 Subunits with FragIT™ Z Kit



For applications that require fast and easy digestion of mouse IgG2a without residual enzyme in the preparation, the FragIT Z spin columns containing FabRICATOR Z immobilized on agarose beads can be used. With FragIT Z kit, the homogenous F(ab')<sub>2</sub> and Fc fragments can be separated with affinity purification columns for preparation of pure antibody fragments.

**Schematic overview of FragIT Z and FragIT Z kit.** Generation and purification of F(ab')<sub>2</sub> and Fc fragments from mouse IgG2a and IgG3 using FragIT Z and FragIT Z kit.

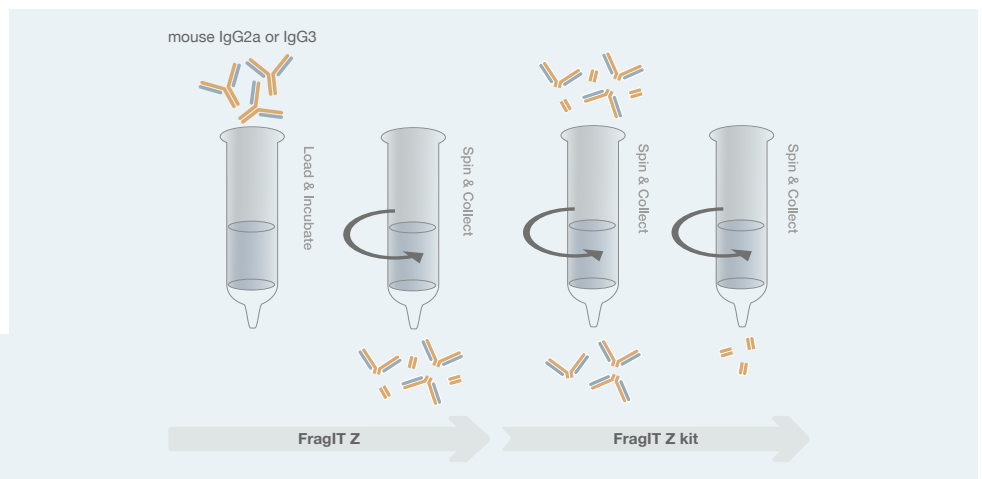

### FabRICATOR® Z



PRODUCT	DESCRIPTION	ID
FabRICATOR Z, 2000 units	Digests 2 mg IgG	A0-FRZ-020

### FragIT™ Z



PRODUCT	DESCRIPTION	ID
FragIT Z, Microspin 2 × 0.5 mg	Digests 2 × 0.5 mg IgG	A0-FZ6-010
FragIT Z, Microspin 5 × 0.5 mg	Digests 5 × 0.5 mg IgG	A0-FZ6-025
FragIT Z, Microspin 10 × 0.5 mg	Digests 10 × 0.5 mg IgG	A0-FZ6-050

### FragIT™ Z kit



PRODUCT	DESCRIPTION	ID
FragIT Z kit, Microspin 0.5 mg	Digests and purifies 0.5 mg IgG	A2-FZ2-005
FragIT Z kit, Microspin 5 × 0.5 mg	Digests and purifies 5 × 0.5 mg IgG	A2-FZ2-025

# GingisKHAN®


## Above Hinge Digestion of Human IgG1



**GingisKHAN (Kgp) is a cysteine protease that digests human IgG1 at a specific site above the hinge, generating intact and homogenous Fab and Fc fragments. The enzyme is used to characterize antibody-based biotherapeutics using LC-MS, to study Fc glycan analysis, bispecific antibodies, affinity and avidity effects and to identify post-translational modifications.**



### Antibody Subunit Workflow

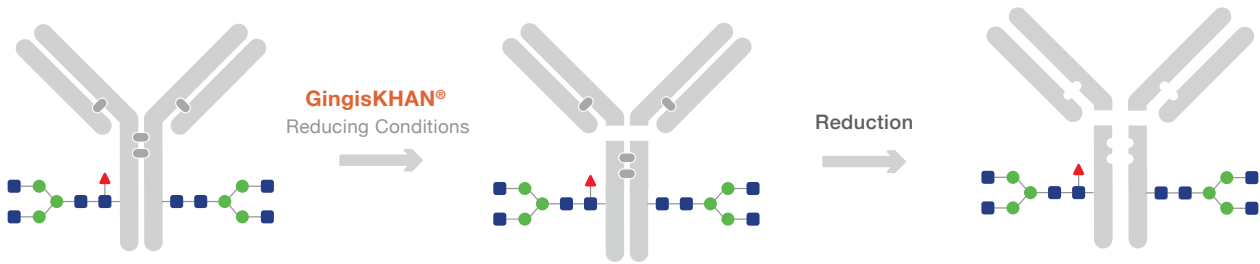

GingisKHAN is a cysteine protease that digests human IgG1 at a specific site above the hinge, generating intact Fab

and Fc fragments. A second digestion site on the Fc may appear if the N-glycans are removed. GingisKHAN

requires mild reducing conditions to be active, and optimal activity is obtained at 37°C and pH 8.

### Product Formats



**GingisKHAN®**  
Lyophilized enzyme for digestion of human IgG1




**GingisKHAN® Fab kit**  
Lyophilized enzyme and CaptureSelect™ CH1 for purification of human IgG1 subunits

Thermo Scientific™ CaptureSelect™ resin from Thermo Fisher Scientific Inc. and its subsidiaries.  
Thermo Scientific and CaptureSelect are trademarks of Thermo Fisher Scientific Inc. and its subsidiaries.



## Digestion of Monoclonal Antibodies and Fc-Fusion Proteins using **GingisKHAN®**

The specific digestion of Fc-fusion proteins or antibodies with mutated hinges into smaller subunits can be challenging. The **GingisKHAN** enzyme digests human IgG1 at a single exposed lysine site above the hinge under reducing conditions. The required reducing condition is mild enough to preserve the intra- and interchain disulphide bridges, resulting in intact Fab and Fc fragments. Here, a selection of monoclonal human IgG1 antibodies and a Fc-fusion protein were incubated with **GingisKHAN** for 1 hour at 37°C and analyzed using HPLC. The resulting fragments displayed the efficient digestion of **GingisKHAN** for the generation of intact Fab fragments from trastuzumab and adalimumab and the separation of the Fc from the fused TNFR domain of etanercept.

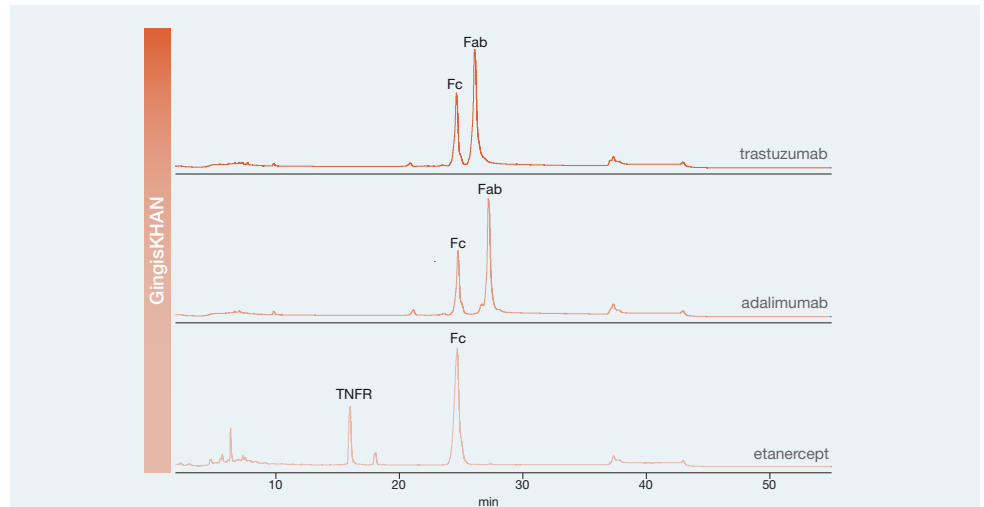


**Specific digestion of human IgG1 using **GingisKHAN**.** LC separation of Fab and Fc fragments from a selection of therapeutic antibodies (trastuzumab and adalimumab) and a Fc-fusion protein (etanercept) after **GingisKHAN** digestion.

## Pure Human IgG1 Fab Fragments using **GingisKHAN®** Fab Kit



Pure and intact human IgG1 Fab fragments can be generated and purified using the **GingisKHAN** Fab kit. Here, trastuzumab Fab fragments were generated after **GingisKHAN** digestion. The fragments were affinity purified using the CaptureSelect™ human CH1 spin column. The enzyme together with the trastuzumab Fc was yielded in the flow-through, and the Fab was easily eluted by lowering the pH.

**SDS-PAGE analysis under non-reducing conditions of the human IgG1 trastuzumab digested by **GingisKHAN**.** After digestion, the Fab fragments were captured in the affinity spin column and could easily be eluted.

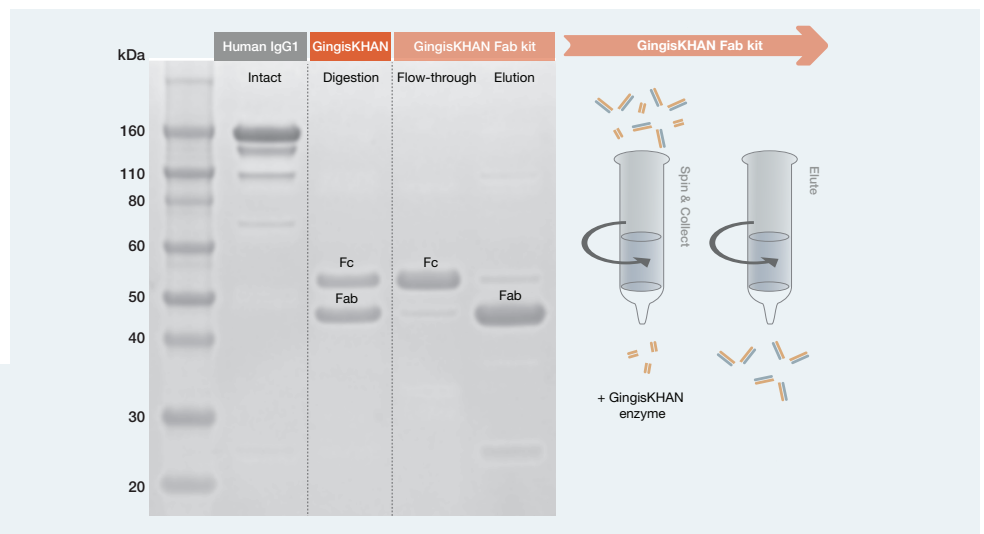

### **GingisKHAN®**



PRODUCT	DESCRIPTION	ID
GingisKHAN, 2000 units	Digests 2 mg human IgG1	B0-GKH-020


### **GingisKHAN®** Fab kit



PRODUCT	DESCRIPTION	ID
GingisKHAN Fab kit	Digests and purifies 2 mg human IgG1	B0-GFK-020

# FabULOUS™

## Above Hinge Digestion of IgG

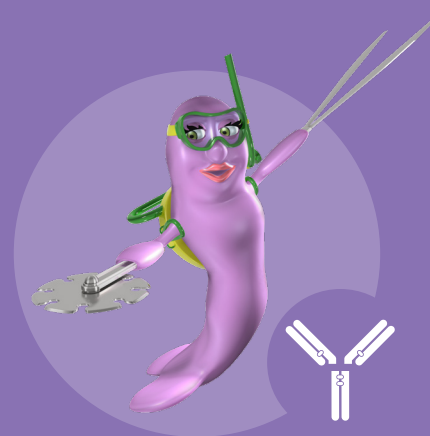

**FabULOUS (SpeB) is a cysteine protease that digests in the hinge region of IgG from several different species and subclasses, generating a homogenous pool of Fab fragments. The prepared Fab fragments can for example be used in affinity studies, studies of Fab glycosylation and structural studies.**

- 


Human IgG and IgG from mouse, rat, goat, sheep and rabbit.
- 

KTHT / CPPCPAP (human IgG1)
- 

60 min reaction
- 

Requires reducing conditions

### Antibody Digestion Workflow

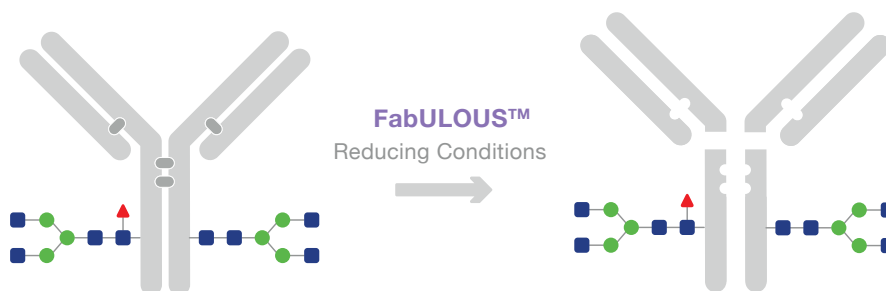

FabULOUS digests IgG and generates Fab and Fc fragments. The primary digestion site on human IgG1 is between the amino acids T225 and C226. The FabULOUS enzyme digests IgG from several species

and subclasses, including mouse (such as mouse IgG1), rat, goat, sheep and rabbit. The enzyme requires reducing conditions for activity on IgG, and if strongly reducing conditions are used,

it is likely that the interchain disulfide bonds will be reduced. The reaction is performed at neutral pH under reducing conditions and the digestion is complete after incubation at 37°C for 1 hour.

### Product Formats



**FabULOUS™**  
Lyophilized enzyme for digestion of IgG




**FabULOUS™ Fab kit**  
Lyophilized enzyme and CaptureSelect™ LC-kappa (mur) for purification of mouse IgG subunits

Thermo Scientific™ CaptureSelect™ resin from Thermo Fisher Scientific Inc. and its subsidiaries.  
Thermo Scientific and CaptureSelect are trademarks of Thermo Fisher Scientific Inc. and its subsidiaries.



## Generating Fab and Fc Fragments using FabULOUS™

The preparation of intact antibody fragments with enzymes can be challenging for many different species and subclasses of IgG due to unspecific enzymatic activity. FabULOUS is a valuable research tool that can generate intact Fab fragments from several species of IgG, including mouse and human. For some species of IgG, FabULOUS is active under mild reducing conditions that allows for the generation of intact Fab fragments without reduction of interchain disulfides. Intact Fab fragments can also be achieved if a buffer exchange is performed after digestion at stronger reducing conditions. This enables reformation of the disulfide bonds in the Fab.

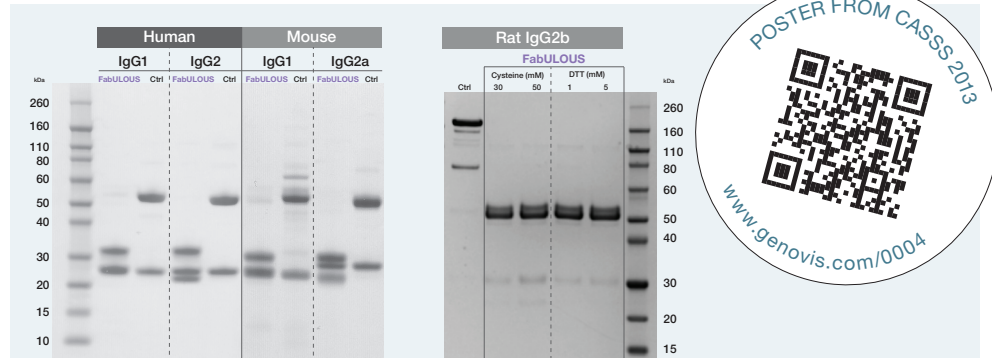


**Generating Fab and Fc Fragments using FabULOUS.** Left: SDS-PAGE analysis of human IgG1 and IgG2 and mouse IgG1 and IgG2a digested in the hinge region using FabULOUS. The generated fragments were analyzed with SDS-PAGE under reducing conditions. Right: SDS-PAGE analysis of rat IgG2b digested in the hinge region using FabULOUS at different reducing conditions. The generated fragments were buffer exchanged and analyzed with SDS-PAGE under non-reducing conditions.

## Generating and Purifying Fab Fragments from Mouse IgG using FabULOUS™ Fab Kit

The FabULOUS Fab kit is designed to generate and purify Fab fragments from mouse IgG. The kit consists of the lyophilized FabULOUS enzyme for antibody digestion and CaptureSelect™ LC-Kappa (mur) affinity spin columns for easy purification of the prepared Fab fragments from mouse IgG. The antibody is digested within 60 minutes using the lyophilized FabULOUS enzyme, and subsequently the prepared Fab fragments bind to the affinity spin columns and are easily eluted. The prepared Fab fragments from mouse IgG1 using FabULOUS Fab kit are shown here.

**Generation and purification of Fab fragments from mouse IgG1.** The top panel shows the intact monoclonal mouse IgG1 antibody, the middle panel shows the analysis of the fragments after FabULOUS digestion, and the bottom panel shows the eluted Fab fragments from the CaptureSelect™ LC-Kappa (mur) column.

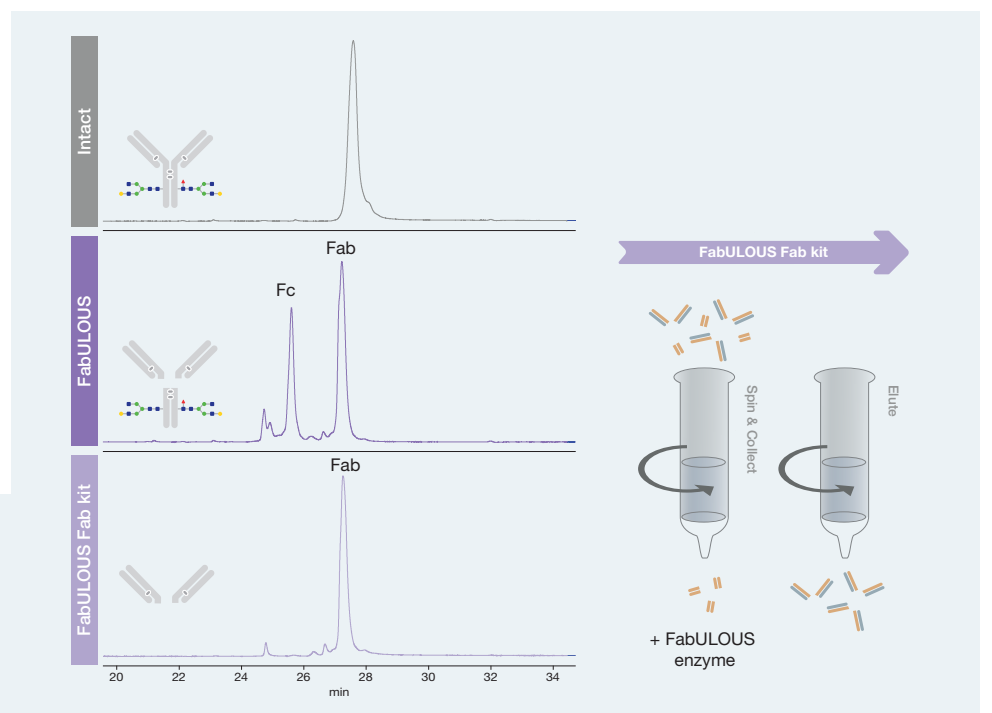

### FabULOUS™

PRODUCT	DESCRIPTION	ID
FabULOUS, 2000 units	Digests 2 mg IgG	A0-PU1-020

### FabULOUS™ Fab kit

PRODUCT	DESCRIPTION	ID
FabULOUS Fab kit, mouse	Digests and purifies 2 mg mouse IgG	A1-PFK-020

# GlycINATOR®

## Hydrolysis of IgG Fc glycans



**GlycINATOR (EndoS2)** is an IgG-specific endoglycosidase that hydrolyzes all glycoforms present at the Fc glycosylation site. The enzyme acts on native IgG and leaves the core GlcNAc with or without fucose. The enzyme is used to reduce the complexity of antibody-based therapeutics to study core afucosylation, and as an initial step in the GlyCLICK and TransGLYCIT platforms.

- 


Human IgG1-4, Fc-fusion proteins, IgG from mouse, rabbit, rat, monkey, sheep, goat, cow and horse
- 

Hydrolyzes the  $\beta$ 1,4-linkage between the two innermost GlcNAc residues
- 

30 min reaction
- 

Requires native IgG fold

### Antibody Deglycosylation Workflow

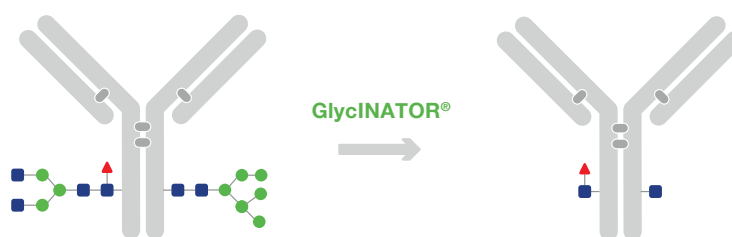

The GlycINATOR enzyme requires a native fold of the Fc in order to trim the N-glycan. It hydrolyzes the  $\beta$ 1,4-linkage between the two innermost GlcNAc residues of

the Fc glycan, leaving the core GlcNAc intact on the Fc. The Fc glycosylation site is conserved among many species, and GlycINATOR deglycosylates antibodies

from several species and subclasses. It effectively hydrolyzes all IgG glycoforms, including complex, high-mannose, hybrid and bisected glycans.

### Product Formats



**GlycINATOR®**  
Lyophilized enzyme for deglycosylation of IgG



**Immobilized GlycINATOR®**  
Immobilized enzyme for deglycosylation of IgG in spin columns

### GlycINATOR®




PRODUCT	DESCRIPTION	ID
GlycINATOR, 2000 units	Deglycosylates 2 mg IgG	A0-GL1-020
GlycINATOR LE (low endotoxin), 2000 units	Deglycosylates 2 mg IgG	A0-GL8-020

## Hydrolysis of Fc glycans using GlycINATOR®



To characterize the activity of GlycINATOR on a therapeutic antibody, trastuzumab was incubated with GlycINATOR for 30 min at 37°C. The antibody was then digested using the FabRICATOR (IdeS) enzyme to generate F(ab')<sub>2</sub> and Fc/2 fragments suitable for LC-MS analysis. The MS results demonstrate hydrolysis of the glycan structure after the core GlcNAc with or without the core fucose when the antibody is incubated with GlycINATOR.


**LC-MS characterization of Fc glycan hydrolysis by GlycINATOR.** Native and GlycINATOR-treated trastuzumab was digested with FabRICATOR to generate Fc/2 fragments before LC-MS analysis.

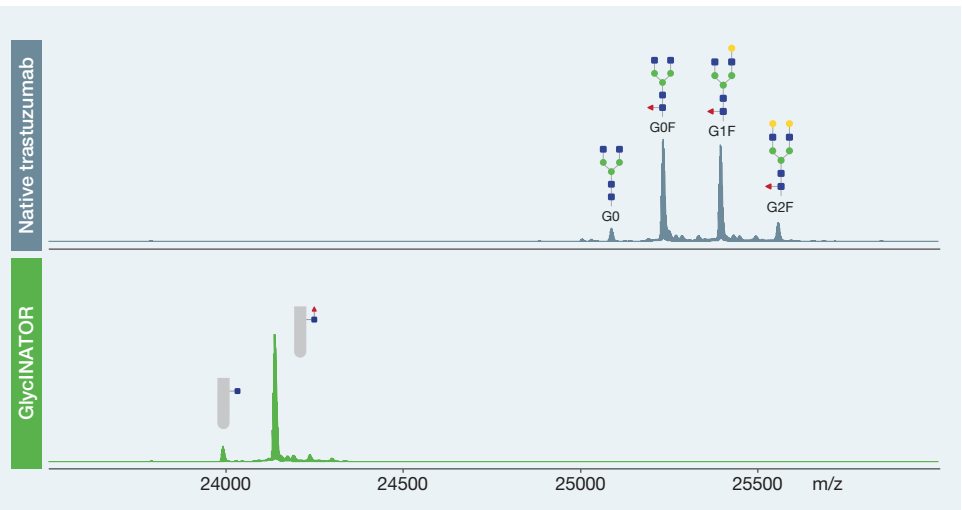

## Antibodies Deglycosylated with GlycINATOR® do Not Elicit ADCC



Afucosylation of the Fc N-glycans found on IgG can impact the binding to Fc receptors and thus the antibody-dependent cellular cytotoxicity (ADCC). Here, the ability of native and deglycosylated trastuzumab to elicit ADCC was analyzed using the cell based iLite® ADCC assay. The native antibody showed an expected increase in signal when more antibody was added to the assay. However, the antibody deglycosylated by GlycINATOR showed no signal even at larger doses, indicating an abolished ADCC activity after GlycINATOR treatment.

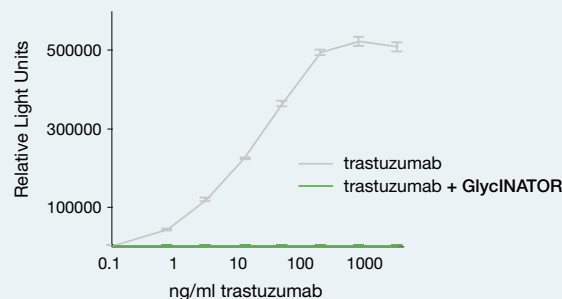

**Study of ADCC by native and deglycosylated trastuzumab.** The effect of trastuzumab was analyzed using iLite ADCC Effector (V) Assay Ready Cells (BM4001) and iLite ADCC Target HER2 (+) Assay Ready Cells (BM4011). Effector and Target cells were diluted and mixed with the antibody solutions. Incubation was performed at 37°C and 5% CO<sub>2</sub>. Luminescence reading was performed after 6 h.




## Immobilized GlycINATOR®



PRODUCT	DESCRIPTION	ID
Immobilized GlycINATOR, Microspin 2 × 0.5 mg	Deglycosylates 2 × 0.5 mg IgG	A0-GL6-010
Immobilized GlycINATOR, Microspin 5 × 0.5 mg	Deglycosylates 5 × 0.5 mg IgG	A0-GL6-025
Immobilized GlycINATOR, Microspin 10 × 0.5 mg	Deglycosylates 10 × 0.5 mg IgG	A0-GL6-050
Immobilized GlycINATOR, Midispin 1-10 mg	Deglycosylates 1-10 mg IgG	A0-GL6-100
Immobilized GlycINATOR, Maxispin 10-100 mg	Deglycosylates 10-100 mg IgG	A0-GL6-1000

# IgZERO®

## Hydrolysis of Complex Fc glycans

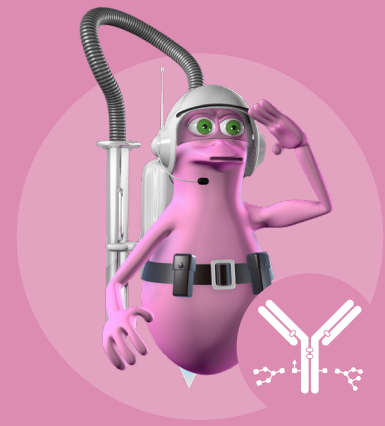

**IgZERO (EndoS)** is an IgG-specific endoglycosidase hydrolyzing complex N-glycans at the Fc glycosylation sites. The enzyme deglycosylates IgG after the core GlcNAc and display limited activity on high-mannose and hybrid-type glycans. The enzyme is used to rapidly reduce sample complexity, to inactivate Fc-mediated effector functions in immunoassays and as a tool to improve imaging by reducing Fc interactions.

- 


Human IgG1-4, Fc-fusion proteins, IgG from mouse, rabbit, rat, monkey, sheep, goat, cow and horse
- 

Hydrolyzes the  $\beta$ 1,4-linkage between the two innermost GlcNAc residues
- 

30 min reaction
- 

Requires native IgG fold

## Antibody Deglycosylation Workflow

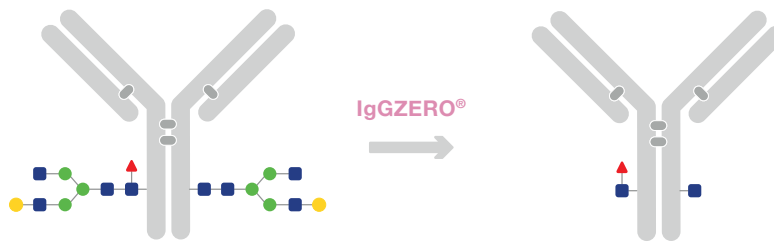

The IgZERO enzyme hydrolyzes N-glycans specifically at the Fc glycosylation site on IgG. In contrast to GlycINATOR,

IgZERO has a limited activity on high-mannose and hybrid-type glycoforms. The hydrolysis of complex glycans is fast

and carried out under native and physiological conditions, since the enzyme requires native IgG fold for activity.

## Product Formats



**IgZERO®**  
Lyophilized enzyme for deglycosylation of IgG



**deGlycIT™**  
Immobilized enzyme for deglycosylation of IgG in spin columns

## IgZERO®



PRODUCT	DESCRIPTION	ID
IgZERO, 1000 units	Deglycosylates 1 mg IgG	A0-IZ1-010
IgZERO, 5000 units	Deglycosylates 5 mg IgG	A0-IZ1-050
IgZERO LE (low endotoxin), 2000 units	Deglycosylates 2 mg IgG	A0-IZ8-020





## Comparison of Glycosidases for N-glycans

GlycINATOR, IgGZERO and EndoH are endoglycosidases acting between the core GlcNAc of the N-glycan, whereas PNGase F is a general amidase acting

on the bond between the core GlcNAc and the peptide chain. GlycINATOR and IgGZERO are specific for the Fc glycans of IgG, and EndoH is a general endoglyco-

sidase for high-mannose type glycoforms. PNGase F acts broadly on all mammalian N-glycoforms but cannot hydrolyze glycans with  $\alpha$ 1,3 core fucosylation.

	GLYCINATOR	IGGZERO	PNGASE F	ENDO H
Digestion site				
Works on native IgG	Yes	Yes	Yes	Yes/No
pH optimum	7.4	7.4	7.5	5-6
Reaction time	30 min	30 min	15 min +	1-16 h
IgG specific	Yes*	Yes	No	No
Glycoform specificity	All Fc glycoforms	Complex Fc glycoforms, limited activity on high-mannose and hybrid-type glycans	All glycoforms except core $\alpha$ 1,3 fucosylation	High-mannose and some hybrid-type glycans

\* GlycINATOR has one more known substrate,  $\alpha$ -1-acid glycoprotein.


## Degree of Fucosylation in Clone Selection using IgGZERO®



Absence of a core fucose on the Fc N-glycan of a therapeutic antibody is associated with an increased antibody-dependent cellular cytotoxicity (ADCC), a desirable property for many therapeutic antibodies. Since IgGZERO rapidly reduces sample complexity by removing the heterogeneous Fc glycans and leaves only the potentially fucosylated, innermost GlcNAc residue, interpretation of MS data is greatly simplified. A technique for studying the degree of fucosylation using IgGZERO has been developed by researchers at Institute Pasteur and LFB Biotechnologies. Briefly, the cell supernatant was concen-

trated, treated with FabRICATOR and IgGZERO, reduced, and then analyzed using LC-MS. This allowed for the easy

determination of the fucosylation levels of the antibody and enhanced clone selection (Henninot et. al., 2015).

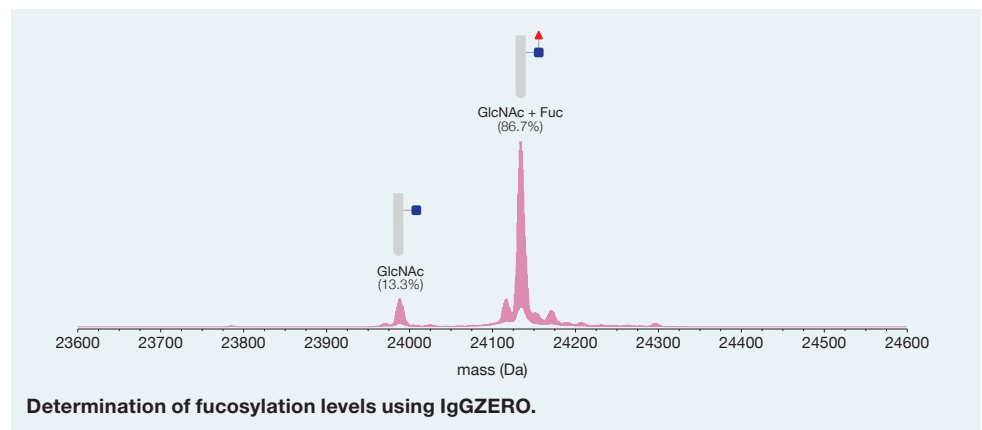

### References

Henninot, A. et al. *Characterization of monoclonal antibodies by a fast and easy LC- MS ToF analysis on culture supernatant.* Anal Biochem 1-8 (2015).

## deGlyt™




PRODUCT	DESCRIPTION	ID
deGlycIT, Microspin 2 × 0.5 mg	Deglycosylates 2 × 0.5 mg IgG	A0-IZ6-010
deGlycIT, Microspin 5 × 0.5 mg	Deglycosylates 5 × 0.5 mg IgG	A0-IZ6-025
deGlycIT, Microspin 10 × 0.5 mg	Deglycosylates 10 × 0.5 mg IgG	A0-IZ6-050
deGlycIT, Midispin 1-10 mg	Deglycosylates 1-10 mg IgG	A0-IZ6-100
deGlycIT, Maxispin 10-100 mg	Deglycosylates 10-100 mg IgG	A0-IZ6-1000


# GlyCLICK®


## Site-Specific Conjugation of Native IgG




**GlyCLICK is a three-step conjugation technology for site-specific and quantitative conjugation of IgG from several species and subclasses. Specific Fc glycan hydrolysis allows for site-specific conjugation on the core GlcNAcs using robust click-chemistry and results in a degree of label (DOL) or antibody-drug ratio (DAR) of 2.0. The reliable performance ensures quantitative conjugation and intact immunoreactivity for sensitive applications, and the technology is available in a range of kit formats.**

- 


Human IgG1-4, IgG from mouse, rabbit, rat, monkey, sheep, goat, cow and horse
- 

Conjugation occurs at the Fc N-glycan sites
- 

2-3 day protocol
- 

Alexa Fluor®, biotin, DFO, MMAE, PNU, Azide activation

### Conjugation Workflow

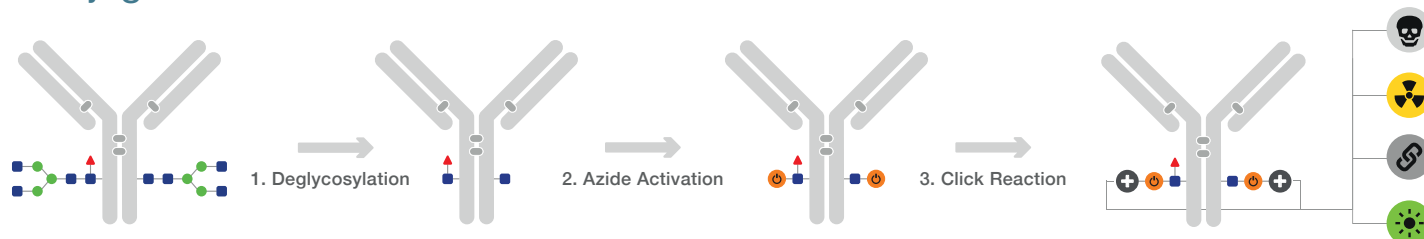

**Deglycosylation:** The Fc-specific endoglycosidase GlycINATOR (EndoS2) hydrolyzes the Fc glycans to the innermost GlcNAc moiety on several subclasses and species of IgG. The enzyme removes all glycoforms, including high-mannose, hybrid, complex, and bisecting type glycans.

**Azide Activation:** GalNAz is enzymatically attached to the exposed GlcNAc using  $\beta$ -1,4-Galactosyltransferase (Y289L)<sup>1</sup> and UDP-GalNAz, to generate an azide-activated antibody that is reactive with a cyclooctyne-functionalized label of choice such as sDIBO or DBCO-reactive labels.

**CLICK Reaction:** The azide-activated antibody is conjugated with a cyclooctyne-functionalized label of choice in a bioorthogonal reaction through strain-promoted copper-free click chemistry (SPAAC) to form a covalent conjugation via a stable triazole.

### Product Formats



#### GlyCLICK® Fluorophore

Site-specific conjugation of IgG with Alexa Fluor® 488, 555 or 647



#### GlyCLICK® Affinity

Site-specific conjugation of IgG with biotin



#### GlyCLICK® Chelator

Site-specific conjugation of IgG with DFO (desferrioxamine)



#### GlyCLICK® ADC

Site-specific conjugation of IgG with MMAE or PNU<sup>2</sup>




#### GlyCLICK® Azide Activation

Site-specific conjugation of IgG with any alkyne reactive click-chemistry



## A New Generation of Antibody Conjugates using GlyCLICK®

Conjugation is conventionally performed using random strategies by labeling at the available amines on lysine residues, or at the solvent accessible inter-chain cysteines. Although versatile, such methods offer limited control of the degree of labeling and site-occupancy and are thereby prone to influence stability, immunoreactivity and reproducibility of the final product. Site-specific strategies have become an attractive option to overcome these limitations, but often require genetic engineering to control and stabilize the process. Conjugation at the Fc glycan sites using GlyCLICK is not

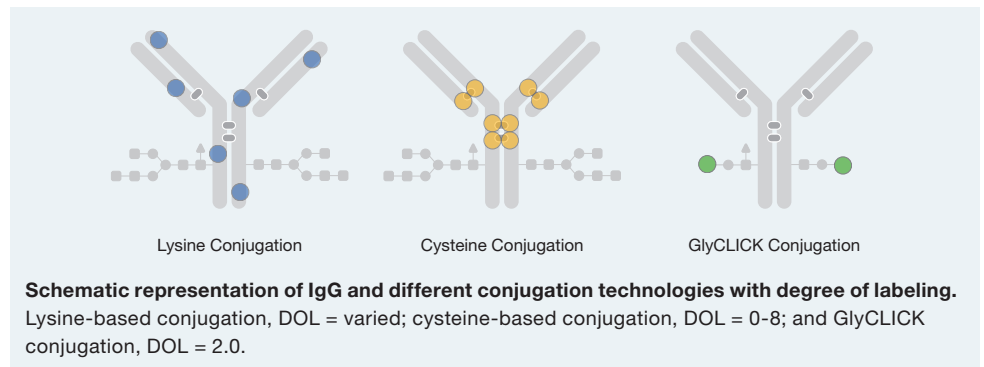

only site-specific, but also decreases the off-target effects of the IgG by


disrupting the binding to Fc receptors (FcγRs) on immune cells.

## Customized Labels for Versatile Applications using GlyCLICK®



The combination of enzymatic remodeling and copper-free click chemistry (SPAAC) enables GlyCLICK to site-specifically conjugate antibodies using any cyclooctyne-label of choice such as sDIBO or DBCO labels. To demonstrate the homogenous material obtained using GlyCLICK, panitumumab was conjugated with different labels and analyzed at the F(ab')<sub>2</sub> and Fc/2 subunit level by RP-HPLC. The resulting separation shows shifts in retention detected only for the labeled Fc/2 subunits when compared to the native counterpart, demonstrating that GlyCLICK specifically conjugates at the Fc glycan sites.

**Analysis of Fc glycan conjugation using GlyCLICK.** RP-HPLC analysis of unmodified panitumumab and panitumumab conjugated with desferrioxamine (DFO), Alexa Flour® 488, biotin or monomethyl auristatin E (MMAE) after digestion with FabRICATOR.

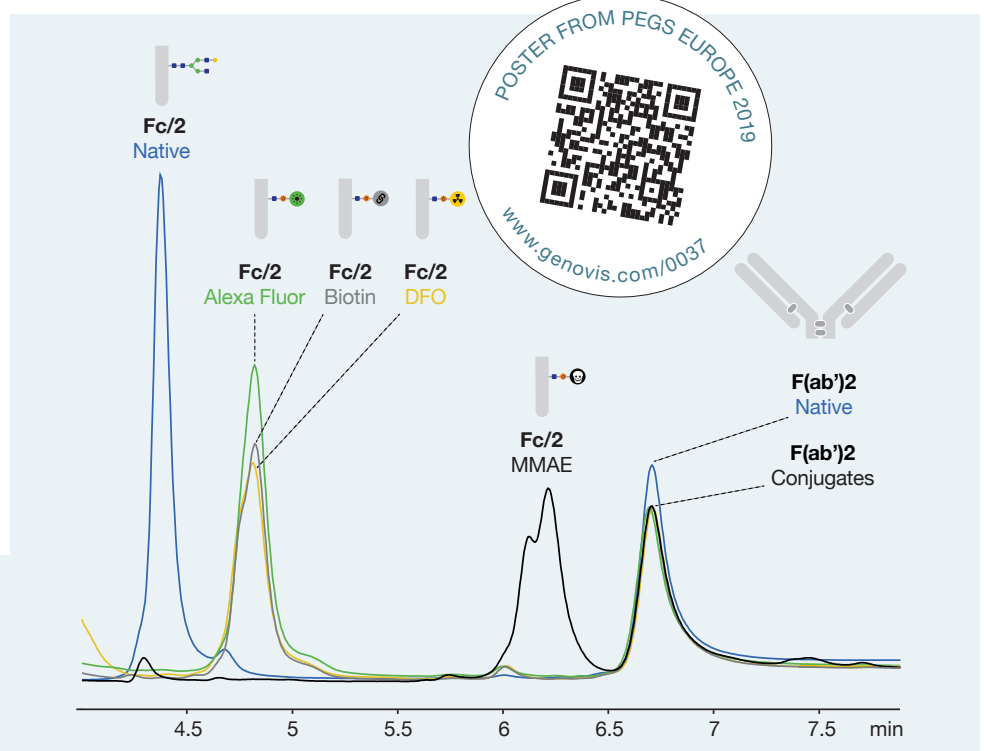

1. SiteClick™ is provided under an intellectual property license from Life Technologies Corporation.

The trademark SiteClick™ is the property of Life Technologies Corporation.

2. Linker-toxin Payloads included in GlyCLICK® are provided under an intellectual property license from Glykos Finland Ltd.


## Reliable Antigen Detection with GlyCLICK®



Immunostaining using randomly labeled antibodies can influence the image quality through inadequate signals and unspecific binding due to excessive labeling or impaired immunoreactivity from labels residing in the antigen-binding Fab region.

To evaluate the immunoreactivity after labeling, trastuzumab was site-specifically conjugated to Alexa Fluor® 488 using GlyCLICK or by random labeling at the lysines using NHS-activated DyLight® 488. Surface plasmon resonance (SPR)

response curves show that only 20% of the randomly labeled material retained binding capacity while the GlyCLICK conjugated antibodies displayed fully preserved immunoreactivity for sensitive and reliable detection.

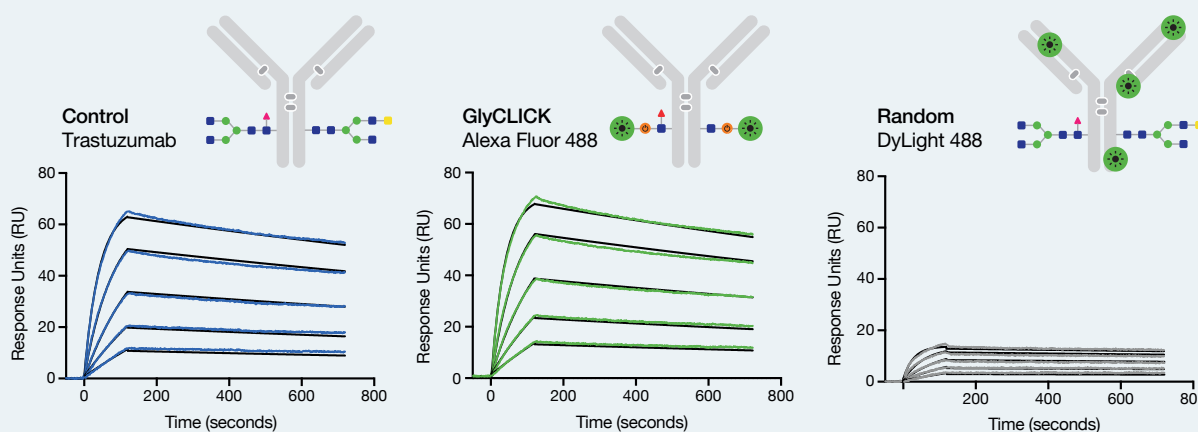

**Affinity analysis of native and conjugated trastuzumab.** Anti-human IgG (Fc) was used as the capturing molecule for trastuzumab. Native (control), Alexa Fluor® 488 (GlyCLICK, DOL=2), DyLight® 488 (random, DOL=10). HER2 was injected in a range to ensure sufficient curvature. All data were fitted against a 1:1 mathematical model.

## In Vivo Performance of GlyCLICK® Conjugated Antibodies




Antibodies radiolabeled by using the site-specific GlyCLICK technology or produced by random conjugation at lysines were used in PET/CT imaging to evaluate their performance *in vivo*. Once chelated, the tracers were labeled with the standard radioisotope Zirconium-89 (<sup>89</sup>Zr) before

injected into tumor-bearing mice. Analysis of the injected mice was carried out using PET/CT imaging to measure the biodistribution and tumor uptake of the tracers.

The results demonstrate a superior tumor uptake and significantly longer circulation time of the site-specific

GlyCLICK tracer, as compared to the randomly labeled tracer. The consistent number of labels attached per antibody using the GlyCLICK kit (DOL = 2) does not only improve performance *in vivo*, but also allows for better quantitation possibilities in PET-imaging experiments.

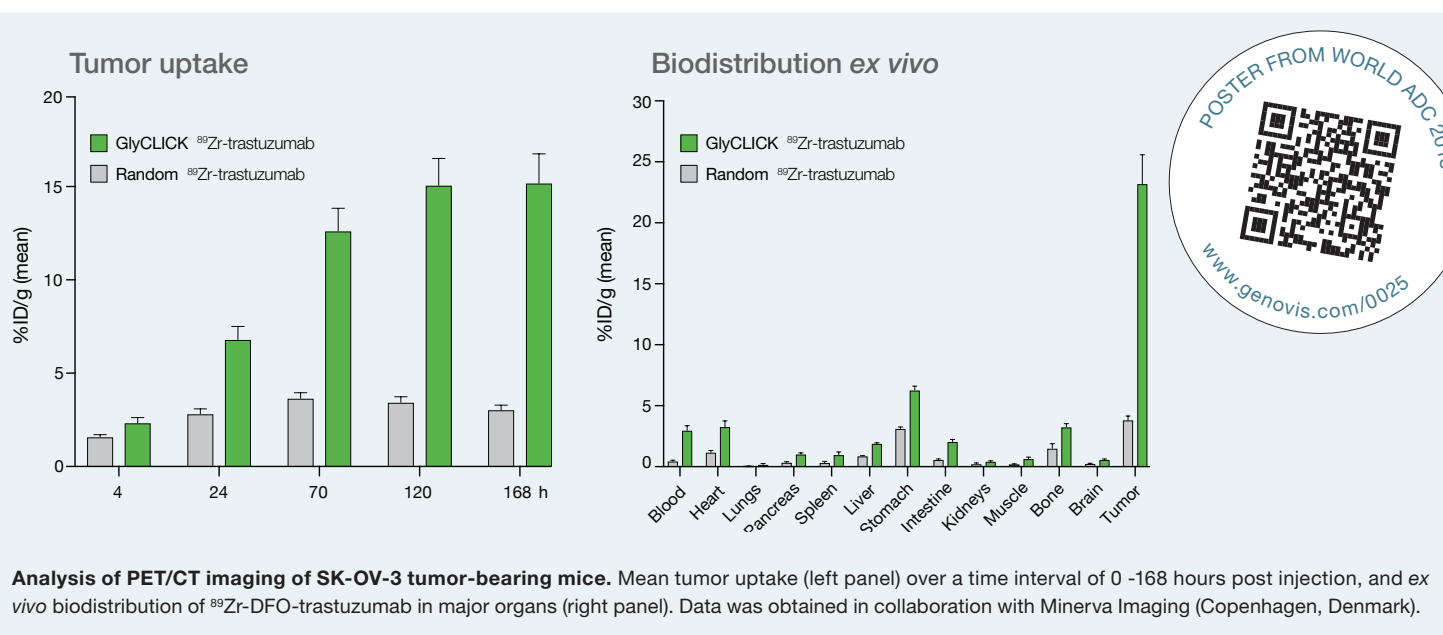

**Analysis of PET/CT imaging of SK-OV-3 tumor-bearing mice.** Mean tumor uptake (left panel) over a time interval of 0-168 hours post injection, and *ex vivo* biodistribution of <sup>89</sup>Zr-DFO-trastuzumab in major organs (right panel). Data was obtained in collaboration with Minerva Imaging (Copenhagen, Denmark).






# Generation of Highly Potent ADCs using GlyCLICK® ADC Kit



The GlyCLICK ADC kits are designed to combine any native antibody with two-step cleavable linker-payloads carrying either the MMAE or PNU toxin. The site-specific GlyCLICK technology

results in complete conjugation that incorporates linker-payloads to an antibody giving a DAR of 2.0. The kits generate homogenous ADCs with two-step cleavable linker-payloads for

intracellular drug release, and produce conjugates that show a cytotoxic effect in the pM range, indicating functional and highly potent ADCs.

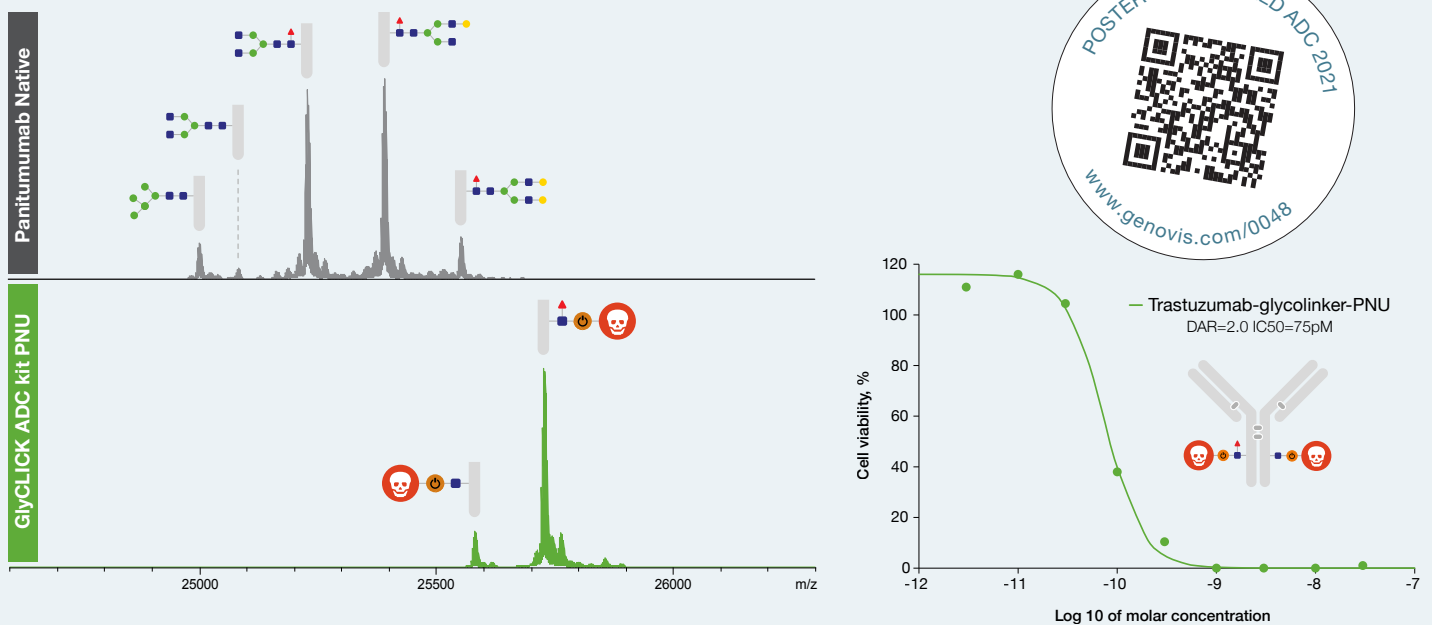

## GlyCLICK® Fluorophore



PRODUCT	DESCRIPTION	ID
GlyCLICK Alexa Fluor® 488	Conjugates 250 µg IgG	L1-F01-025
GlyCLICK Alexa Fluor® 488	Conjugates 1 × 2 mg IgG	L1-F01-200
GlyCLICK Alexa Fluor® 555	Conjugates 250 µg IgG	L1-F02-025
GlyCLICK Alexa Fluor® 555	Conjugates 1 × 2 mg IgG	L1-F02-200
GlyCLICK Alexa Fluor® 647	Conjugates 250 µg IgG	L1-F03-025
GlyCLICK Alexa Fluor® 647	Conjugates 1 × 2 mg IgG	L1-F03-200

## GlyCLICK® Affinity



PRODUCT	DESCRIPTION	ID
GlyCLICK Biotin	Conjugates 250 µg IgG	L1-A01-025
GlyCLICK Biotin	Conjugates 1 × 2 mg IgG	L1-A01-200

## GlyCLICK® Chelator



PRODUCT	DESCRIPTION	ID
GlyCLICK DFO	Conjugates 250 µg IgG	L1-C01-025
GlyCLICK DFO	Conjugates 1 × 2 mg IgG	L1-C01-200

## GlyCLICK® ADC



PRODUCT	DESCRIPTION	ID
GlyCLICK ADC kit, MMAE	Conjugates 1 × 2 mg IgG	L1-T02-200
GlyCLICK ADC kit, PNU	Conjugates 1 × 2 mg IgG	L1-T01-200


## GlyCLICK® Azide Activation



PRODUCT	DESCRIPTION	ID
GlyCLICK Azide Activation	Activates 250 µg IgG	L1-AZ1-025
GlyCLICK Azide Activation	Activates 1 × 2 mg IgG	L1-AZ1-200
GlyCLICK Azide Activation	Activates 1 × 10 mg IgG	L1-AZ1-100

# TransGLYCIT™


## Transglycosylation of IgG

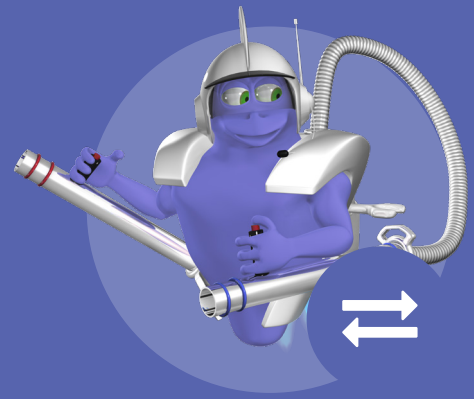

TransGLYCIT is a platform that enables efficient and site-specific human IgG glycan remodeling and conjugation. With TransGLYCIT, antibodies with defined and homogenous glycoforms or site-specific azide activation are generated using fast and robust enzymatic workflows. With the optional FucosEXO™ 16 enzyme in the glycan remodeling workflow, afucosylated antibodies carrying the Fc glycoform of choice (G0, G1, G2 or G2S2) can be obtained for comparison of antibodies with or without core fucose. With the azide activation workflow, antibodies are prepared for conjugation with any label of choice using click-chemistry.

-  Human IgG1, IgG2 and IgG4
-  Transglycosylation occurs at the Fc N-glycan sites
-  3-6 hour workflow
-  G0, G1, G2, G2S2 or azide-activated glycoform is included

### The Transglycosylation Workflow

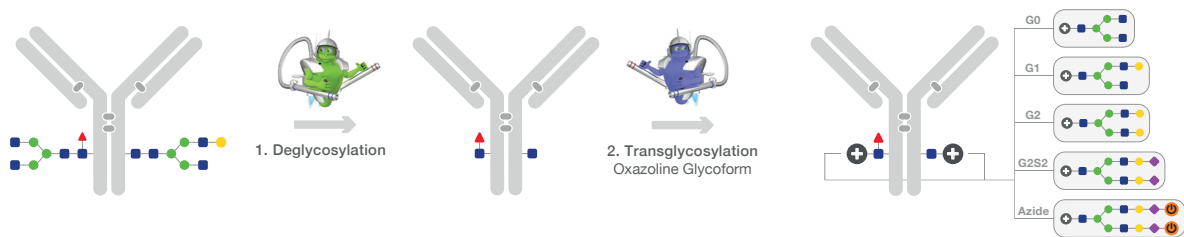

**Deglycosylation:** The Fc N-glycans are trimmed to the core GlcNAc using the IgG-specific Immobilized GLYCINATOR® (EndoS2) enzyme that hydrolyses all Fc glycoforms, including high-mannose, hybrid, complex and bisecting glycans. The degree of fucosylation is the same as on the original IgG.

Note that for the glycan remodeling products, it is possible to obtain afucosylated G0, G1, G2 or G2S2 glycoforms using TransGLYCIT Afucosylated, where the Immobilized FucosEXO™ 16 enzyme hydrolyzes the  $\alpha$ 1,6-linked core fucose prior to the transglycosylation step.

**Transglycosylation:** The engineered glycosynthase TransINATOR™ catalyzes the transglycosylation reaction between the oxazoline reactive G0, G1, G2, G2S2 or azide-activated glycoform and the core GlcNAc.

### Product Formats



#### TransGLYCIT™ Azide Activation

Azide activation of human IgG1/IgG4 or IgG2



#### TransGLYCIT™

Glycan remodeling of human IgG with the G0, G1, G2 or G2S2 glycoform



#### TransGLYCIT™ Afucosylated

Glycan remodeling of human IgG with the G0, G1, G2 or G2S2 glycoform and core afucosylation




## Generating Agalactosylated, Galactosylated, or Sialylated Fc Glycoforms using TransGLYCIT™

Glycoengineering of the IgG Fc N-glycan profile is important for the development of next-generation therapeutic antibodies with enhanced or silenced effector functions. TransGLYCIT is used to generate antibodies carrying agalactosylated (G0), galactosylated (G1, G2) or sialylated (G2S2) glycan structures. With the option to generate glycan structures lacking the core fucose using the TransGLYCIT

Afucosylated kit, antibodies that show increased binding to activating FcγIIIa receptors and thus an elevated ADCC response can be obtained for direct comparison between fucosylated and afucosylated antibodies.

Here, the glycan profile of the therapeutic antibody trastuzumab was remodeled to carry either fucosylated or afucosylated G0, G1, G2 or G2S2

glycoforms. The resulting mass spectra show the heterogeneous glycan profile of native trastuzumab and the mass shifts after complete transglycosylation of the antibody to generate G0, G1, G2 or G2S2 using TransGLYCIT. Homogenous afucosylated glycan profiles were generated using the TransGLYCIT Afucosylated kits.

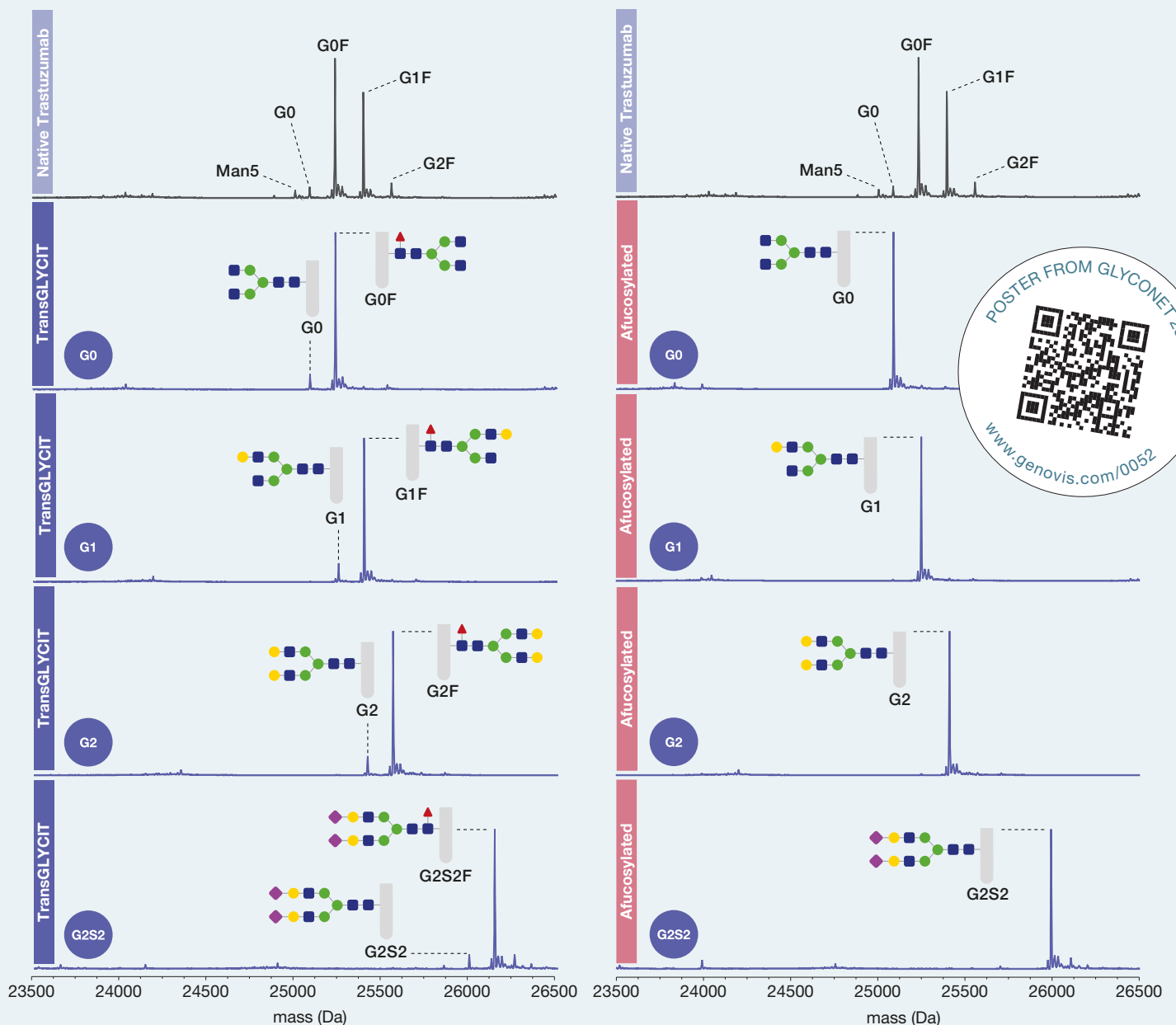

**Transglycosylation using TransGLYCIT.** Deconvoluted mass spectra of the Fc/2 fragment of native trastuzumab (top) and after transglycosylation with TransGLYCIT (left panel) or TransGLYCIT Afucosylated (right panel). The mAb was digested with FabRICATOR and the subunits were analyzed by reversed-phase LC-MS on a Waters™ BioAccord™ system equipped with a Waters™ BioResolve™ RP mAb column (2.1 × 50 mm).


## Site-specific Azide Activation of Human IgG using TransGLYCIT™ Azide Activation



Site-specific Fc N-glycan conjugation is a powerful tool with applications ranging from biomedical research to the development of diagnostic methods and antibody-drug conjugates (ADCs). A reliable conjugation technology is crucial to obtain a homogenous and

reproducible conjugate while preserving the stability and immunoreactivity of the final product. Using the TransGLYCIT Azide Activation technology, human IgG is azide-activated by transglycosylation within a few hours. Here, TransGLYCIT Azide Activation was used for site-specific

azide activation of trastuzumab using the modified N-glycan oxazoline carrying two azide functionalities. The reaction generated an IgG with four azide-activated sites for subsequent conjugation, two sites on each Fc/2.

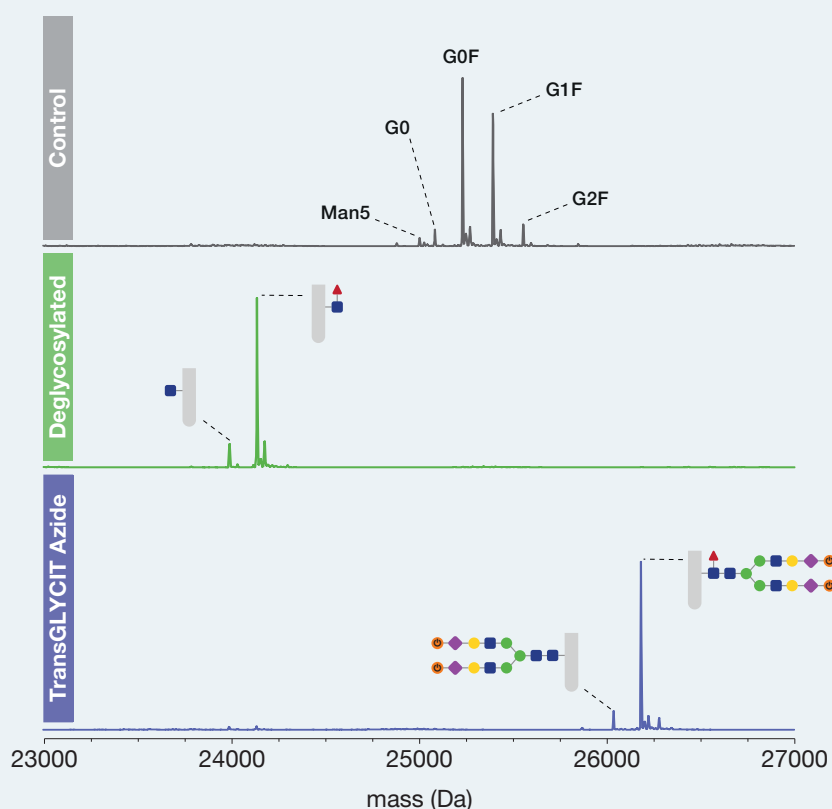

**Azide activation of trastuzumab using the TransGLYCIT Azide Activation workflow.** Deconvoluted mass spectra of the Fc/2 fragment of native trastuzumab (top), after deglycosylation by GlycinATOR (middle) and after transglycosylation by TransINATOR with the oxazoline glycoform azide as the substrate (bottom). The mAb was digested with FabRICATOR and the subunits were analyzed by reversed-phase LC-MS on a Waters™ BioAccord™ system equipped with a Waters™ BioResolve™ RP mAb column (2.1×50mm).

### TransGLYCIT™

PRODUCT	DESCRIPTION	ID
TransGLYCIT G0, 1 mg	Generates 1 mg human IgG with the G0 glycoform	T1-G0F-010
TransGLYCIT G1, 1 mg	Generates 1 mg human IgG with the G1 glycoform	T1-G1F-010
TransGLYCIT G2, 1 mg	Generates 1 mg human IgG with the G2 glycoform	T1-G2F-010
TransGLYCIT G2S2, 1 mg	Generates 1 mg human IgG with the G2S2 glycoform	T1-S2F-010

### TransGLYCIT™ Azide Activation

PRODUCT	DESCRIPTION	ID
TransGLYCIT Azide Activation, 100 µg	Activates 100 µg human IgG1 or IgG4	T1-AZ1-001
TransGLYCIT Azide Activation hIgG2, 100 µg	Activates 100 µg human IgG2	T1-AZ2-001






## Comparison of Site-specific and Random Conjugation Technologies

The azide-activated Fc N-glycans on antibodies modified using the TransGLYCIT Azide Activation workflow can be used for site-specific conjugation with a label of choice via click chemistry, resulting in an antibody homogeneously labeled with four labels per molecule. Using site-specific conjugation instead of random conju-

gation technologies is favorable for preserving affinity and generating conjugates where quantitative properties are desired. To illustrate the impact of different labeling methods on the degree of labeling (DOL), we here analyzed trastuzumab conjugated with a fluorophore using either TransGLYCIT, GlyCLICK or

random labeling by NHS-chemistry. The TransGLYCIT and GlyCLICK workflows resulted in homogeneously labeled samples with defined DOL 4 or DOL 2, respectively, whereas the random technique resulted in a heterogeneously labeled sample with a DOL varying from 2 to 10.

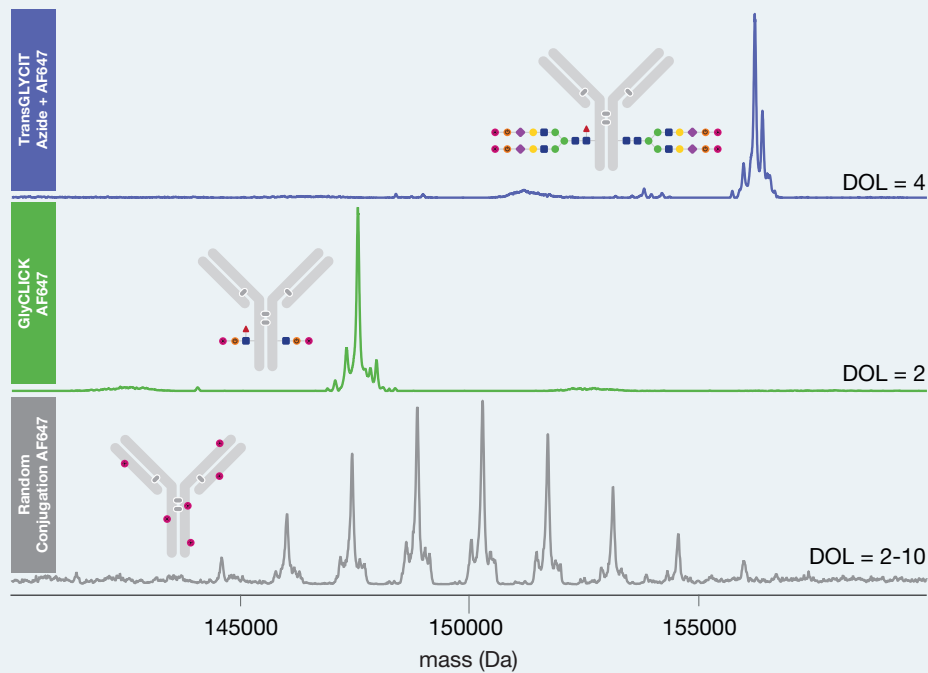

**Conjugation of trastuzumab using site-specific and random technologies.** Deconvoluted mass spectra of intact trastuzumab modified using TransGLYCIT Azide Activation followed by site-specific conjugation by strain-promoted azide alkyne click chemistry using DBCO-AF647 (DOL=4; top), conjugated using GlyCLICK AlexaFluor®647 (DOL=2; middle), and randomly labeled using NHS-activated AlexaFluor 647 resulting in heterogenous labeling (DOL=2-10; bottom). The randomly labeled sample was analyzed after deglycosylation by GlycINATOR. The intact masses of the trastuzumab conjugates were analyzed using a Bruker Impact II ESI-QTOF MS.

### TransGLYCIT™ Afucosylated

PRODUCT	DESCRIPTION	ID
TransGLYCIT G0 Afucosylated, 1 mg	Generates 1 mg afucosylated human IgG with the G0 glycoform	T1-G0A-010
TransGLYCIT G1 Afucosylated, 1 mg	Generates 1 mg afucosylated human IgG with the G1 glycoform	T1-G1A-010
TransGLYCIT G2 Afucosylated, 1 mg	Generates 1 mg afucosylated human IgG with the G2 glycoform	T1-G2A-010
TransGLYCIT G2S2 Afucosylated, 1 mg	Generates 1 mg afucosylated human IgG with the G2S2 glycoform	T1-S2A-010

# GingisREX®

## Arginine-Specific Protein Digestion



**GingisREX (RgpB)** is an arginine-specific protease that digests proteins C-terminally of arginine residues, including arginine – proline linkages that are difficult to digest with other enzymes. The enzyme can be used in peptide mapping, *de novo* peptide sequencing and analysis of post-translational modifications.

- Digests any peptide or protein containing Arg. Specific for Arg-X motifs
- Digests C-terminally of Arg residues
- 60 min reaction
- Active in 6M urea and 0.1% SDS

### Protein Digestion Workflow



GingisREX is a cysteine protease that specifically digests peptide bonds C-terminally of arginine residues, including sites next to proline. Longer peptides are generated, as compared to trypsin digestion, which can be resolved and

identified by high-resolution MS. This offers more options in sample preparation for bottom-up approaches to achieve alternative digestion profiles and increased sequence coverage. GingisREX does not have activity at

lysines, as commonly observed using Arg-C. The enzyme is active at a broad pH range of 5.0-9.0, and requires cysteine but is inhibited by guanidine hydrochloride.

### Product Formats



#### GingisREX®

Lyophilized enzyme for digestion C-terminally of arginine residues

### GingisREX®

PRODUCT	DESCRIPTION	ID
GingisREX, 5 µg enzyme	Enzyme:Protein ratio: 1:20 - 1:200	B0-GRX-005






## Arginine-Specific Digestion using GingisREX®

Interpretation of bottom-up mass spectrometry data is often based on specific assumptions about the specificity of the protease used to digest the sample. Unspecific activities create uncertainties that make data interpretation difficult and increase both analysis time and the risk for misinterpretations. We compared the activities of GingisREX and Arg-C

using oxidized insulin  $\beta$ -chain as a very simple model substrate. Oxidized insulin  $\beta$ -chain contains only one arginine as well as one lysine residue and should therefore only be digested at one single site by Arg-specific proteases. When analyzing the peptides generated by GingisREX and Arg-C by RP-HPLC and MS, GingisREX showed specific digestion at the arginine

residues with no activity on the lysine residue and no other unspecific activities even after prolonged incubation over-night with a high enzyme to substrate ratio (1:5). In contrast, Arg-C displayed a primary activity at the arginine residue, but digestion could also be observed at the lysine residue as well as a number of other minor sites.

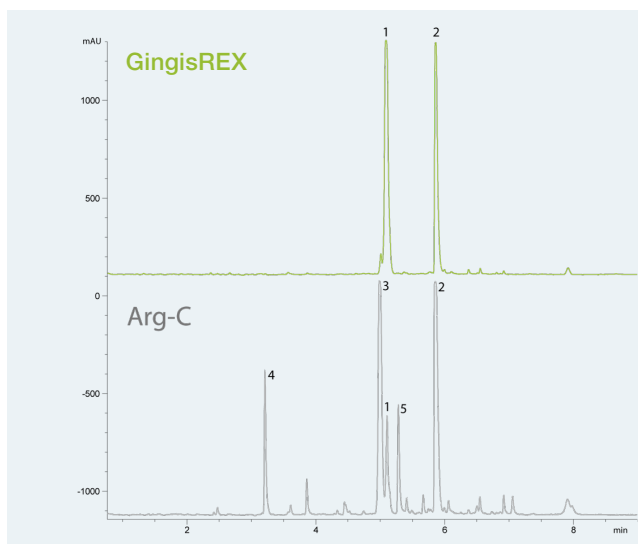

### Digestion of oxidized $\beta$ -chain of insulin with GingisREX and Arg-C.


The digestion was performed O/N at 37°C, with an enzyme to substrate ratio of 1:20 (w/w), 20 mM cysteine in buffers at pH 7.4 (GingisREX) and pH 7.6 (Arg-C). The peptides were separated on RP-HPLC.

PEPTIDE NO.	AMINO ACID SEQUENCE
Intact protein	FVNQHLGSHLVEALYLVCGERGFFYTPKA
1	GFFYTPKA
2	FVNQHLGSHLVEALYLVCGER
3	GFFYTPK
4	FVNQHLGSH
5	LVEALYLVCGER + Na

Sequences of oxidized insulin  $\beta$ -chain digested by GingisREX or Arg-C. Green indicates arginine residues and red indicates lysine residues.

## Peptide Mapping of Therapeutic Antibodies using GingisREX®


On a large and complex sample, such as a therapeutic antibody, the digestion at arginine residues gives larger peptides and results in fewer peaks in the peptide map. Here, the GingisREX and Arg-C digestion profiles of trastuzumab are presented. The digestion by GingisREX at Arg sites resulted in fewer peaks and the peptide map was less complicated, a feature that is beneficial for data interpretation in mass fingerprint analysis.

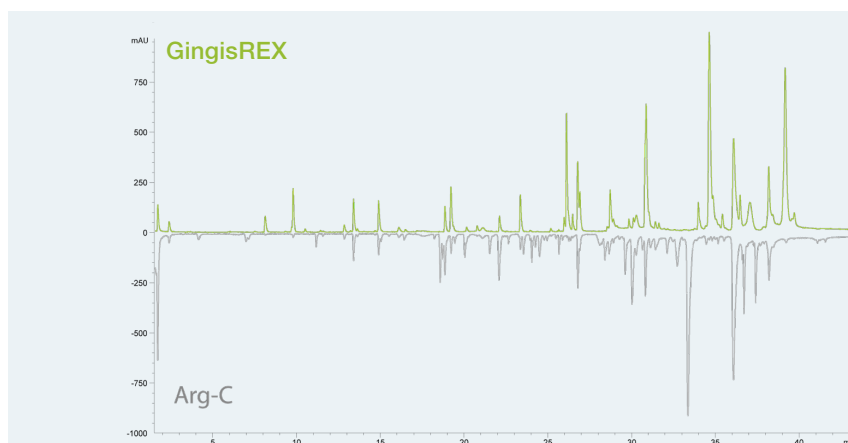

Peptide maps of trastuzumab after digestion using GingisREX or Arg-C.

# PNGase F

## Hydrolysis of N-glycans

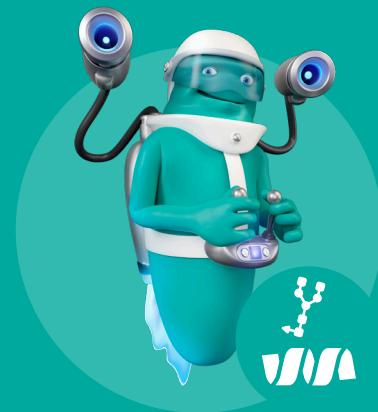

**PNGase F (Peptide N-glycosidase F)** is a glycoamidase hydrolyzing the amide bond between the polypeptide asparagine and the innermost GlcNAc of all mammalian asparagine-linked complex, hybrid, or high-mannose oligosaccharides. The enzyme is widely used for sample preparation prior to MS analysis – to reduce the protein heterogeneity and enable released glycan analysis – and to study the functional role of the N-glycan.

-  N-linked glycans on glycoproteins
-  Hydrolyzes the glycosidic bond between N-glycans and Asn
-  15 + min reaction time
-  *RapiGest™* included in the Immobilized PNGase F Denaturing kit

### Deglycosylation Workflow



PNGase F hydrolyzes the glycosidic bond between the N-glycans and asparagine. During the reaction, the asparagine residue

from which the glycan is removed is deamidated to aspartic acid. The released glycan is left intact and can be used for

further analysis. PNGase F is active under both native and denaturing reaction conditions.

### Product Formats



**PNGase F**  
Lyophilized enzyme for  
deglycosylation of glycoproteins



**Immobilized PNGase F**  
Hydrolysis of N-glycans from glycoproteins in  
spin columns under native conditions



**Immobilized PNGase F Denaturing**  
Hydrolysis of N-glycans from glycoproteins in  
spin columns under denaturing conditions


## Removal of N-glycans under Native Conditions using PNGase F



Removing the Fc N-glycans with PNGase F under native conditions enables characterization of the free N-glycans as well as the function and structure of the deglyco-

sylated antibody. Trastuzumab was used to demonstrate the efficient removal of N-glycans by PNGase F and Immobilized PNGase F under native conditions.

The mass shift demonstrates successful removal of the Fc N-glycans with no enzyme interfering in the analysis of the sample processed with Immobilized PNGase F.

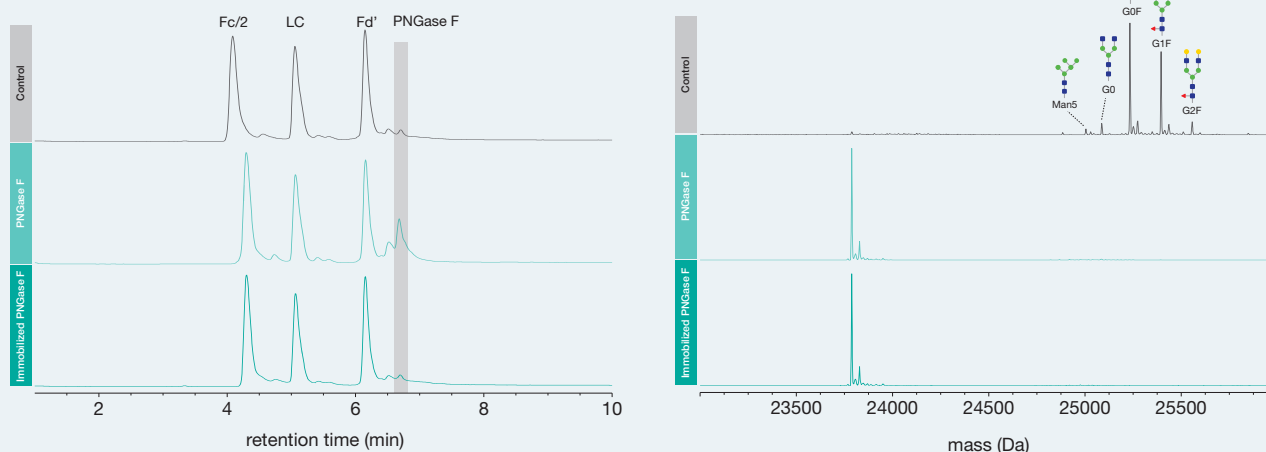

**Removal of N-glycans under native conditions.** TIC chromatogram (left) and deconvoluted mass spectra of the Fc/2 fragment (right) of intact trastuzumab (top panel) and trastuzumab treated with Immobilized PNGase F (bottom panel) or PNGase F (middle panel) for 1 h at 37°C. The antibody samples were analyzed after digestion with FabRICATOR.


## Released Glycan Analysis using PNGase F



The N-glycosylation profile of a therapeutic protein is a critical quality attribute. It affects both safety and efficacy of the biopharmaceutical and therefore needs to be characterized and monitored during development and production. Here, the N-glycans of the therapeutic antibody trastuzumab, the Fc-fusion protein etanercept and the glycoprotein RNase B were analyzed using a released glycan approach. Trastuzumab and etanercept were deglycosylated with PNGase F under native conditions, while RNase B required denaturation and reduction to be fully deglycosylated. The resulting released glycans were labelled with *RapiFluor-MS*<sup>™</sup> (Waters Corporation) and analyzed by HILIC UPLC-FLD-MS.

### Released glycan analysis of glycoproteins.

*RapiFluor-MS*<sup>™</sup> labelled N-glycans released from trastuzumab (top), etanercept (middle) and RNase B (bottom) were analyzed by HILIC-FLD-MS using a Waters<sup>™</sup> BioAccord<sup>™</sup> system equipped with a Waters<sup>™</sup> ACQUITY Premier Glycan BEH Amide column (130Å, 2.1 × 150mm). The fluorescence chromatograms are shown, and the glycan structures were annotated using a combination of glucose unit (GU) library search and MS confirmation.

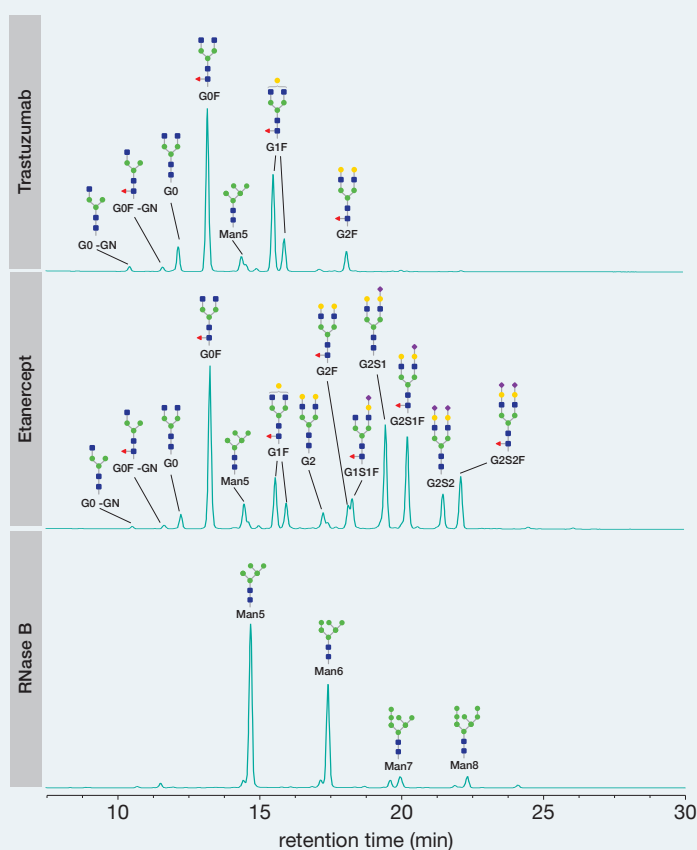



## N-glycan Removal using PNGase F



	TRASTUZUMAB		ETANERCEPT	
Protein type	IgG1		Fusion protein	
N-glycans	2		6	
O-glycans	-		26	
Product format	PNGase F	Immobilized PNGase F	PNGase F	Immobilized PNGase F
Reaction conditions	Native 37°C, 1 h	Native 37°C, 1 h	Native 37°C, 1 h	Native 37°C, 1 h
Preparation of samples for gel	Reduced	FabRICATOR digested	Reduced	Reduced


### Key characteristics and SDS-PAGE assay of a selection of glycoproteins deglycosylated using PNGase F and Immobilized PNGase F.

The table specifies important characteristics of the examined glycoproteins, the reaction conditions used for the deglycosylation and the sample preparations performed before SDS-PAGE analysis. The mass shifts on the gels show the successful removal of N-glycans after PNGase F or Immobilized PNGase F treatment from the various glycoprotein substrates.

## Rapid N-glycan Removal using Immobilized PNGase F Denaturing



Some glycosylation sites are poorly accessible for PNGase F and the deglycosylation reaction is slow or inhibited by steric hindrance. To demonstrate the increase in reaction speed under denaturing conditions, the fusion protein abatacept was processed with Immobilized PNGase F using either the native or the denaturing workflow. Native conditions allowed for removal of N-glycans without any additives or increased temperature but required overnight (18 h) incubation for complete deglycosylation. Removal of the N-glycans using denaturing conditions was complete within 30 minutes. As can be seen from the MS spectra, all N-glycans were removed while the O-glycans were left intact. The intact protein is complex to analyze, but after Immobilized PNGase F treatment, the complexity of the sample was greatly reduced.

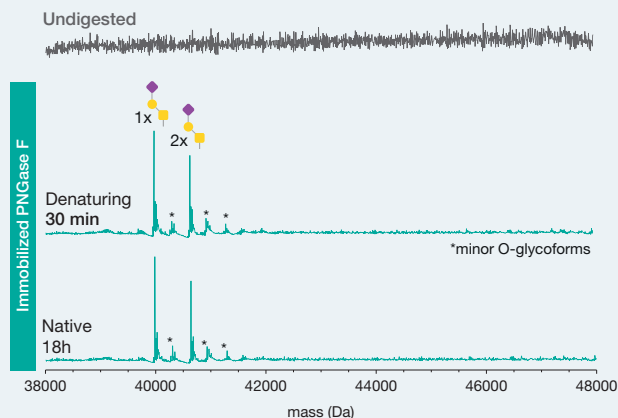

**Deconvoluted mass spectra of abatacept (top panel) treated with Immobilized PNGase F under denaturing (middle panel) or native conditions (bottom panel).** The protein was incubated for 18 h with Immobilized PNGase F using native reaction conditions or 30 min using denaturing conditions. The deglycosylated protein samples were reduced and analyzed by reversed-phase LC-MS on a Waters™ BioAccord™ system equipped with a Waters™ BioResolve™ RP mAb column (2.1 × 50 mm). The mass shifts show the complete removal of N-glycans, leaving only the O-glycoforms.

CETUXIMAB		ABATACEPT	AFLIBERCEPT	RNASE B (BOVINE PANCREAS)							
IgG1		Fusion protein	Fusion protein	Glycoprotein							
4		6	10	1							
-		8	-	-							
PNGase F		Immobilized PNGase F	Immobilized PNGase F	Immobilized PNGase F	PNGase F	Immobilized PNGase F					
Denaturing, 50°C, 10 min		Native (N) 37°C, 1 h Denaturing (D) 50°C, 15 min		Denaturing 50°C, 30 min		Denaturing 50°C, 30 min		Denaturing and reducing, 50°C, 10 min		Denaturing and reducing, 50°C, 15 min	
Reduced		FabRICATOR digested		Reduced		Reduced		None		None	
↓		↓		↓		↓		↓		↓	

## PNGase F



PRODUCT	DESCRIPTION	ID
PNGase F, 1000 units	Deglycosylates 1 mg glycoproteins	G1-PF1-010
PNGase F, 5 × 1000 units	Deglycosylates 5 × 1 mg glycoproteins	G1-PF1-050

## Immobilized PNGase F



PRODUCT	DESCRIPTION	ID
Immobilized PNGase F	5 × 0.2 mg microspin columns	G1-PF6-010
Immobilized PNGase F	10 × 0.2 mg microspin columns	G1-PF6-020


## Immobilized PNGase F Denaturing



PRODUCT	DESCRIPTION	ID
Immobilized PNGase F Denaturing	5 × 0.2 mg microspin columns + 5 × 1 mg <i>RapiGest</i> <sup>™</sup> SF	G2-PDK-010
Immobilized PNGase F Denaturing	10 × 0.2 mg microspin columns + 10 × 1 mg <i>RapiGest</i> <sup>™</sup> SF	G2-PDK-020


# OpeRATOR®


## O-glycan-specific Protein Digestion

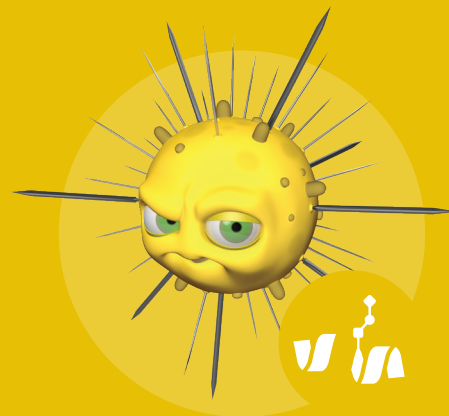

OpeRATOR is an O-glycan-specific protease that digests proteins carrying mucin-type O-glycans N-terminally of the glycosylated Ser and Thr residues. Glycopeptides carrying O-glycans are generated, and the enzyme enables O-glycan profiling, O-glycopeptide mapping and site occupancy determination, as well as middle-level approaches using MS analysis.

- 


Native proteins with mucin-type O-glycosylation
- 

Digests N-terminally of O-glycosylated Ser and Thr residues
- 

2h to O/N (16-18h) reaction
- 

Desialylation using SialEXO (included) increases performance

### O-glycoprotein Digestion Workflow

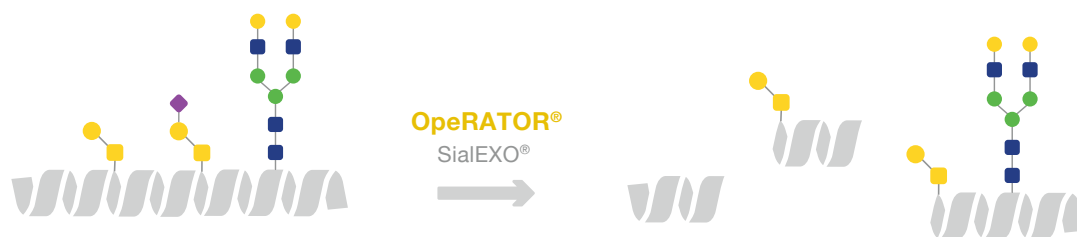

With high specificity, OpeRATOR digests mucin-type O-glycoproteins and peptides N-terminally of the O-glycans, at serine or threonine residues. It will not

digest glycoproteins at N-glycan sites. The enzyme is most active towards sites with asialylated core 1 O-glycans. It also digests sialylated core 1 and

core 3, but to a much lower extent. OpeRATOR is compatible with SialEXO (a sialidase mix), which is included in the purchase.

### Product Formats



#### OpeRATOR®

Lyophilized enzyme for digestion of mucin-type O-glycoproteins and peptides

### OpeRATOR®




PRODUCT	DESCRIPTION	ID
OpeRATOR, 2000 units	Digests 2 mg O-glycoprotein	G2-OP1-020

## O-glycan Site-specific Digestion using OpeRATOR®

OpeRATOR can be used to map O-glycosylation sites at the middle-level. Indeed, for proteins carrying only one or a few O-glycans, the O-glycosylation site can be determined directly from identifying the site of OpeRATOR digestion. Erythropoietin (EPO) is a ~30 kDa glycoprotein with a single O-glycosylation site. After removal of the N-glycans by PNGaseF and desialylation using SialEXO, the native protein was digested N-terminally of the O-glycan site using OpeRATOR. The resulting protein fragments were analyzed by reversed-phase LC-MS.


**Specific digestion N-terminally of the O-glycosylation site.** The reduced fragments were separated on a reversed-phase C4 column followed by ESI-QTOF MS detection. The EPO protein carrying one core 1 O-glycan was hydrolyzed at a single specific site N-terminally of the O-glycosylated serine using OpeRATOR.

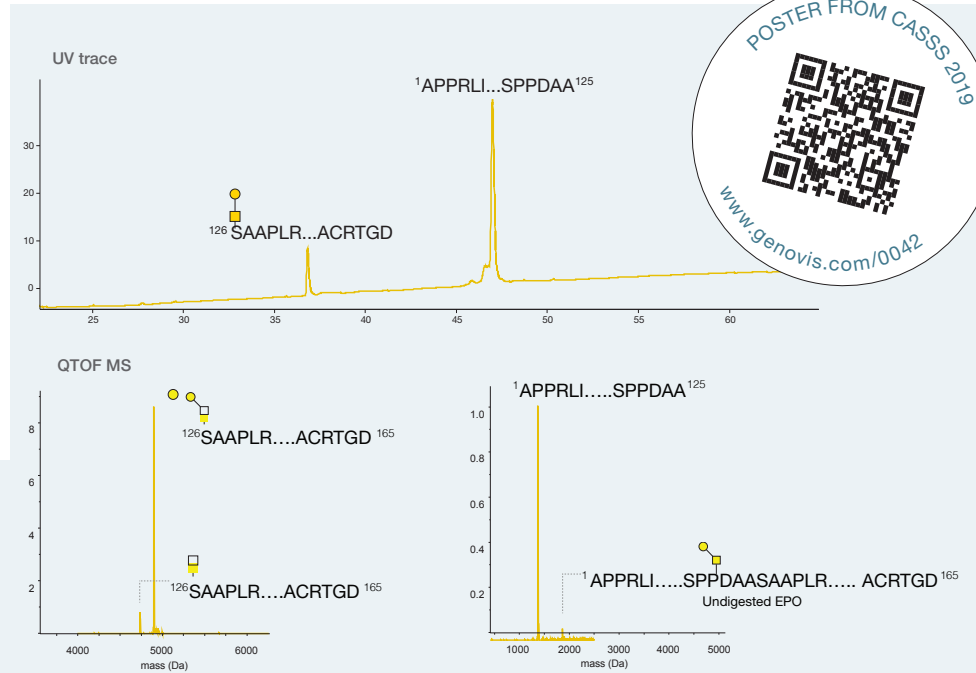

## A Complete Map of Etanercept O-glycan Sites using OpeRATOR®

Etanercept is an Fc-fusion protein with a highly O-glycosylated hinge region that is difficult to analyze by conventional assays. The protein was here incubated with OpeRATOR, and the resulting glycopeptides were analyzed using LC-MS/MS (top panel). Due to the heterogeneity in the O-glycan pattern of the protein and the OpeRATOR specificity for O-glycan structures, overlapping peptides were formed that made it possible to acquire a complete map of the O-glycan sites (middle panel). For comparison, a standard tryptic digest of the same material (bottom panel), shows only few and very large glycopeptides that were difficult to analyze.

**In-depth analysis of the O-glycosylation sites of etanercept originator.** HILIC-MS analysis (top panel) of etanercept peptides generated by digestion with OpeRATOR. The middle panel shows the peptide map of the O-glycosylation sites of etanercept based on the results from the OpeRATOR-based workflow. The identified O-glycosylation sites are marked in yellow. Two sites (marked with an asterisk) could only be inferred from the OpeRATOR digestion pattern without direct identification of a peptide containing the glycosylated amino acid. The bottom panel shows a RP-MS analysis of etanercept digested with trypsin. Data was obtained in collaboration with Thermo Fisher Scientific (UK).

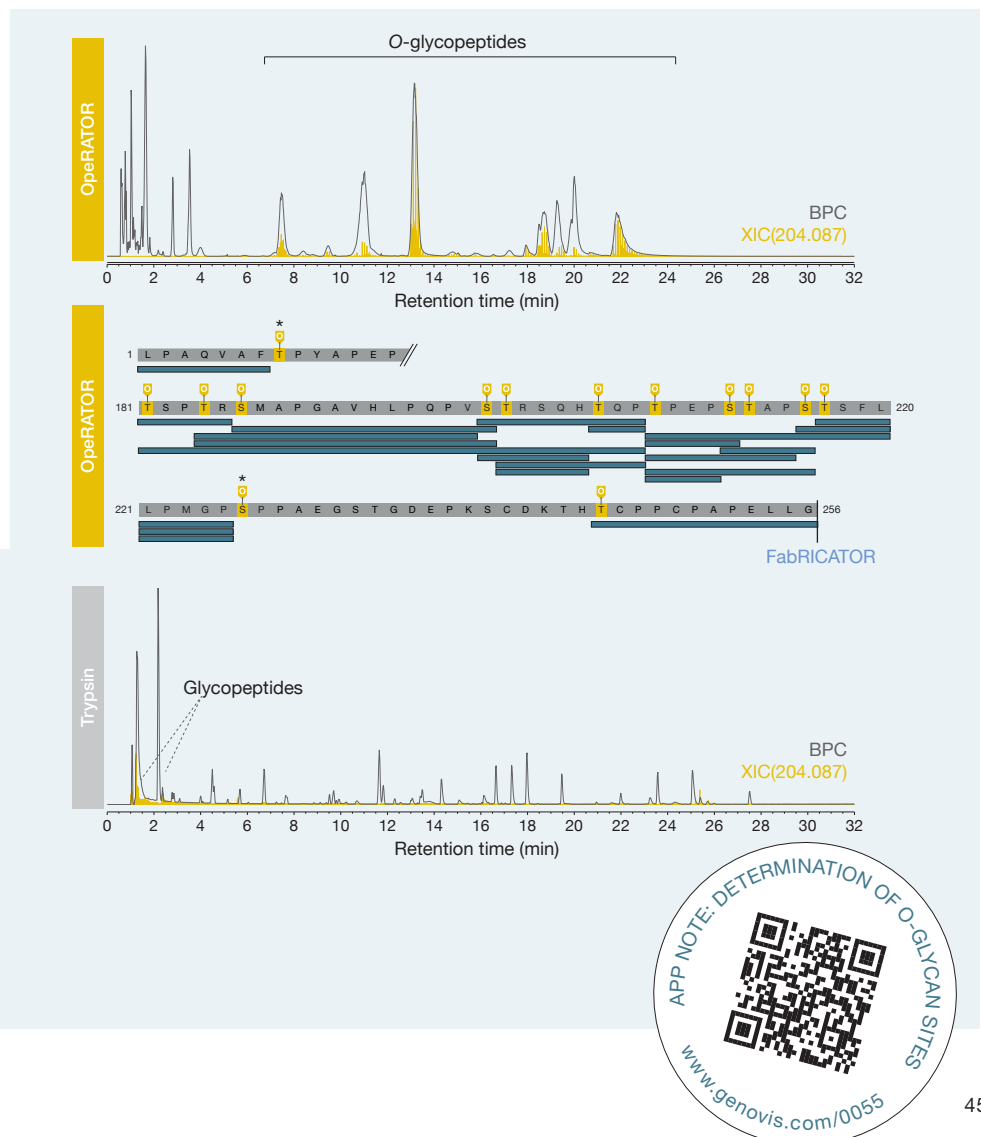

# OglyZOR®

## Hydrolysis of Core-1 O-glycans



**OglyZOR is an endoglycosidase (O-Glycosidase) that specifically hydrolyzes core 1 O-glycan disaccharides on native glycoproteins. OglyZOR is used for removal of O-glycans for glycan analysis, confirmation of O-glycan presence and reduction of sample heterogeneity.**

-  O-glycans on glycoproteins
-  Hydrolyzes Core 1 disaccharides
-  2-4 h reaction
-  Requires removal of terminal sialic acids (SialEXO included)

### O-glycoprotein Deglycosylation Workflow



OglyZOR is an endoglycosidase (endo- $\alpha$ -N-acetylgalactosaminidase) that hydrolyzes core 1 (Gal- $\beta$ 1-3-GalNAc), and to a limited extent core 3 (GlcNAc- $\beta$ 1-3-

GalNAc) O-glycan disaccharides on native glycoproteins. The enzyme does usually not require denaturation of the substrate. Since the terminal sialic acids

of the O-glycans need to be removed for OglyZOR activity, the sialidase mix SialEXO is included. Optimal activity is obtained at pH 6.5 to 7.5 and at 37°C.

### Product Formats



#### OglyZOR®

Lyophilized enzyme for hydrolysis of core 1 O-glycans on glycoproteins

#### OglyZOR®



PRODUCT	DESCRIPTION	ID
OglyZOR, 2000 units	Deglycosylates 2 mg glycoprotein	G2-OG1-020






## Deglycosylation of Biopharmaceuticals using OglyZOR®

Etanercept is a highly glycosylated recombinant Fc-fusion protein and its analysis at the intact level is very challenging. Digestion with FabRICATOR (IdeS) allows for the two subunits (TNFR and Fc/2) to

be analyzed separately in a middle-up LC-MS approach. While the untreated TNFR domain was too complex to be analyzed, the O-glycan composition could be resolved and compared after

treatment with PNGase F. Complete glycan removal provided information of the primary structure of the TNFR part.

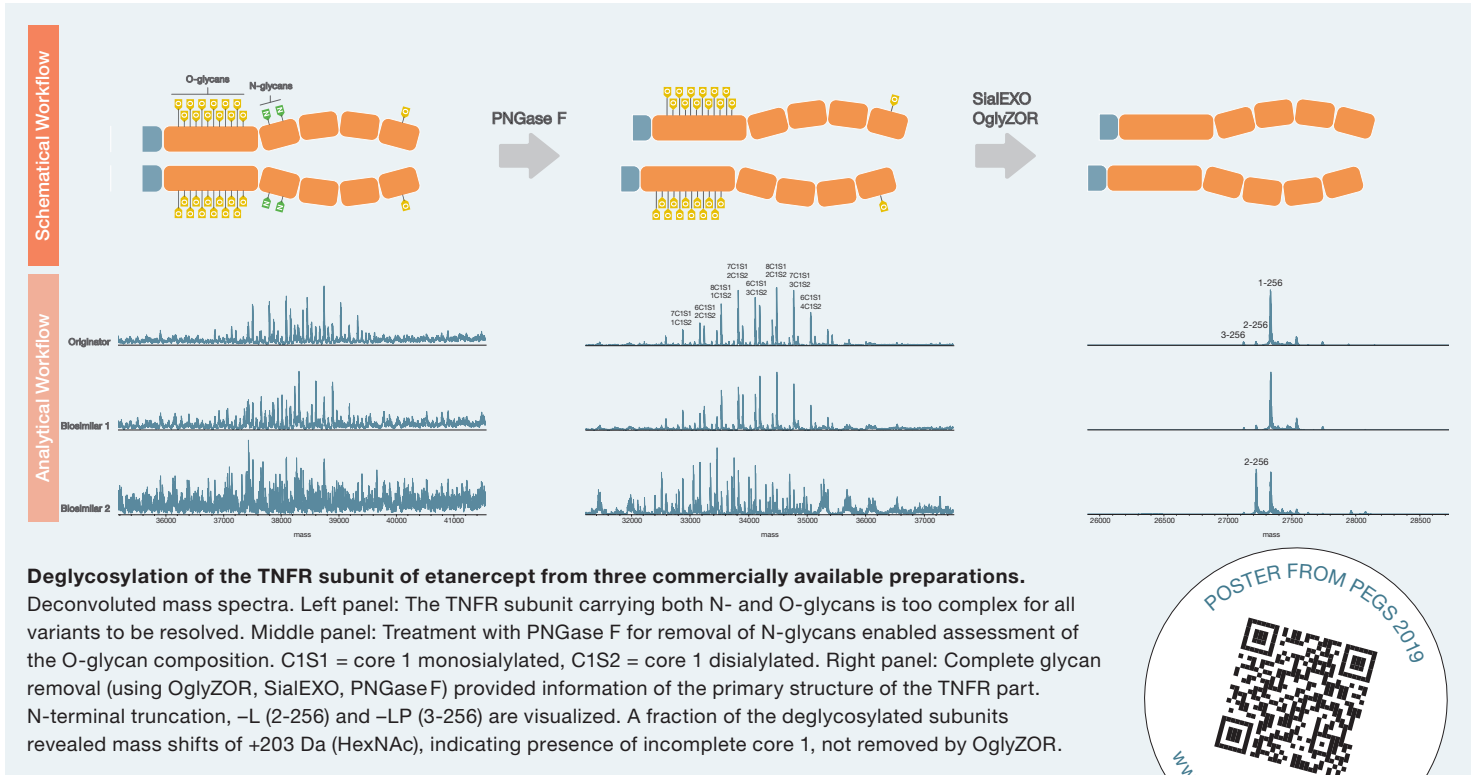


## Benchmarking of OglyZOR® Activity

The enzymatic activities of OglyZOR and SialEXO are here compared to other commercially available endoglycosidases and sialidases. When the sialic acids were removed, the OglyZOR enzyme efficiently

hydrolyzed the O-glycans from TNFR and the intact etanercept molecule. OglyZOR efficiently hydrolyzes core 1 disaccharides on native glycoproteins, making the enzyme suitable for prepara-

tion of samples prior to LC-MS analysis for identification of other PTMs or confirmation of amino acid sequences.

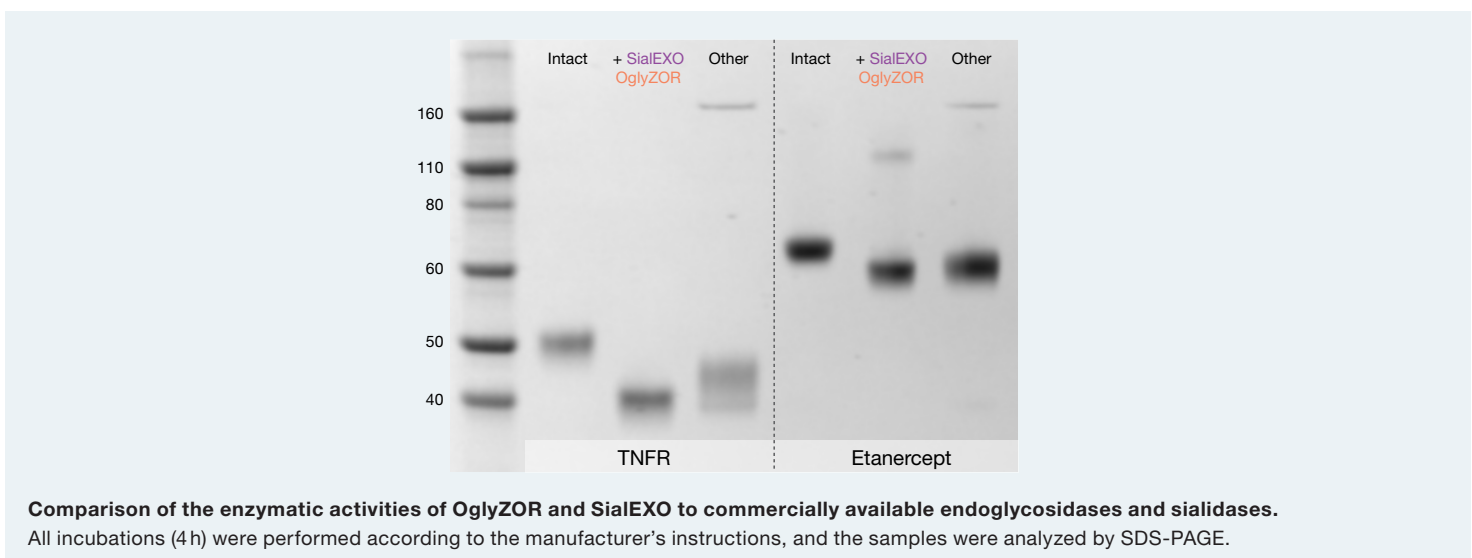

# SialEXO®

## Hydrolysis of Sialic Acids




The SialEXO enzymes are sialidases for the removal and analysis of sialic acids. The enzymes are active on both N- and O-linked glycans present on native glycoproteins or released glycan structures. SialEXO is used to pretreat O-glycosylated proteins prior to digestion with OperATOR or OglyZOR. Other applications include reduction of sample complexity, charge variant analyses and exoglycosidase arrays.

- 

Sialic acids on N- and O-linked glycans
- 

$\alpha$ 2-3,  $\alpha$ 2-6 and  $\alpha$ 2-8-linked sialic acids
- 

30 min – 2 h reaction
- 

Requires no co-factors

### Desialylation Workflow



The inherent negative charge of sialic acids might complicate analytical workflows and mask other important

modifications. The removal of sialic acids therefore facilitates the study of underlying variants in the protein.

SialEXO hydrolyzes glycoproteins under native conditions and displays a high activity in a broad pH range, 6.5 to 9.

### Product Formats



#### SialEXO®

Lyophilized enzyme mix for desialylation of glycoproteins



#### SialEXO® 23

Lyophilized enzyme for hydrolysis of  $\alpha$ 2-3-linked sialic acids



#### Immobilized SialEXO®

Immobilized enzymes for desialylation of glycoproteins in spin columns

#### SialEXO®



PRODUCT	DESCRIPTION	ID
SialEXO, 2000 units	Desialylates 2 mg glycoprotein	G1-SM1-020

#### SialEXO® 23




PRODUCT	DESCRIPTION	ID
SialEXO 23, 500 units	Desialylates 0.5 mg glycoprotein	G1-SD2-005

## Analysis of Released Sialic Acids using SialEXO® and SialEXO® 23



Sialic acids present on biopharmaceuticals are primarily attached in  $\alpha$ 2-3 or  $\alpha$ 2-6-linkages. The broad acting sialidase mix SialEXO hydrolyzes all sialic acid linkages, including  $\alpha$ 2-3,  $\alpha$ 2-6, and  $\alpha$ 2-8. SialEXO 23 is an enzyme that display high specificity for  $\alpha$ 2-3-linked sialic acids. Here, hydrophilic interaction liquid chromatography (HILIC) was used to analyze two glycan libraries. One library contained released glycans with  $\alpha$ 2-3-linked sialic acids and the other contained  $\alpha$ 2-6-linked sialic acids. The shift in retention times clearly showed the release of  $\alpha$ 2-3 and  $\alpha$ 2-6-linked sialic acids from the glycans by SialEXO 23 and SialEXO respectively.

**Analysis of released sialic acids.** HILIC analysis of released glycans with  $\alpha$ 2-3-linked (left) or  $\alpha$ 2-6-linked (right) sialic acids after incubation with SialEXO 23 or SialEXO respectively. Sialylated glycan structures are shaded in purple.

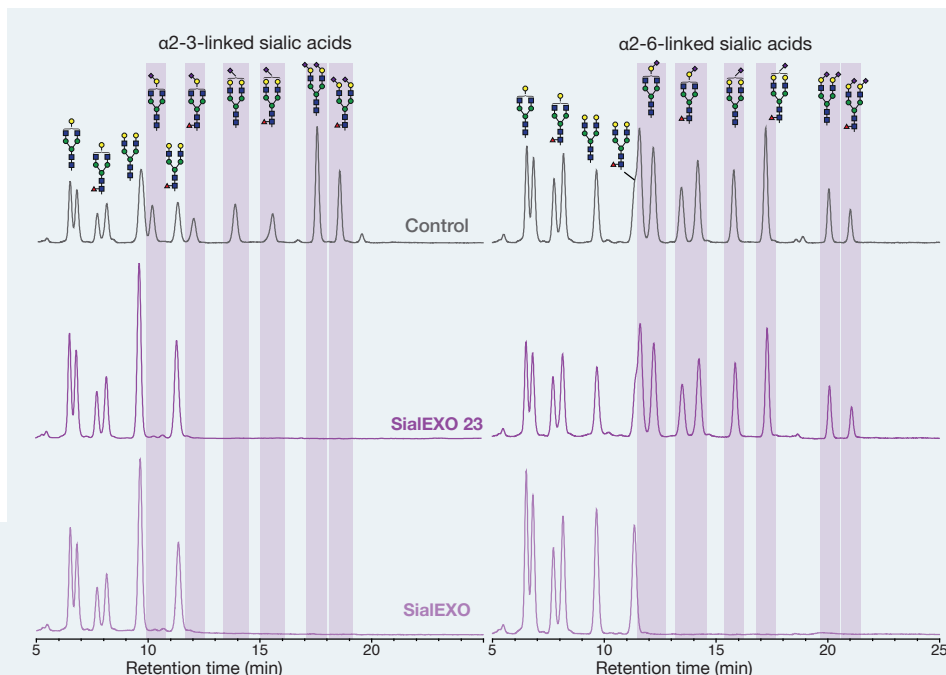


## Charge Variant Analysis using Immobilized SialEXO®



Capillary isoelectric focusing is commonly used to determine charge variants during characterization and quality control of biologics. Etanercept is a highly sialylated Fc-fusion protein which makes the analysis of its charge variants difficult.

Immobilized SialEXO microspin columns contain SialEXO covalently coupled to agarose beads for desialylation of native glycoproteins in 30 minutes without any residual enzyme in the final preparation. To study the underlying charge variants on etanercept, the protein was here desialylated using Immobilized SialEXO and analyzed by imaged capillary isoelectric focusing using Maurice from ProteinSimple.

The desialylation removes the charge heterogeneity originating from the sialic acids and allows for the analysis of the underlying charge variants of the protein.

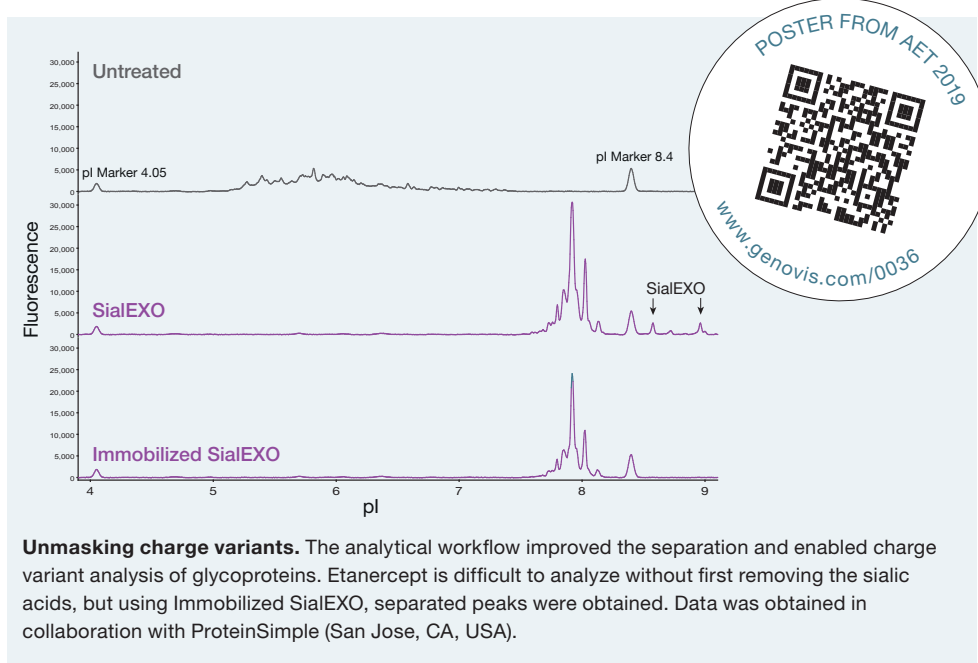

**Unmasking charge variants.** The analytical workflow improved the separation and enabled charge variant analysis of glycoproteins. Etanercept is difficult to analyze without first removing the sialic acids, but using Immobilized SialEXO, separated peaks were obtained. Data was obtained in collaboration with ProteinSimple (San Jose, CA, USA).

## Immobilized SialEXO®




PRODUCT	DESCRIPTION	ID
Immobilized SialEXO, Microspin 2 × 0.5 mg	Desialylates 2 × 0.5 mg glycoprotein	G1-SM6-010
Immobilized SialEXO, Microspin 5 × 0.5 mg	Desialylates 5 × 0.5 mg glycoprotein	G1-SM6-025
Immobilized SialEXO, Microspin 10 × 0.5 mg	Desialylates 10 × 0.5 mg glycoprotein	G1-SM6-050

# FucosEXO™

## Hydrolysis of $\alpha$ 1-2,3,4 Fucose



**FucosEXO is a mix of  $\alpha$ -fucosidases for efficient hydrolysis of  $\alpha$ 1-2,  $\alpha$ 1-3 and  $\alpha$ 1-4-linked fucose residues on native N- and O-glycosylated proteins or free oligosaccharides. It is a valuable tool for glycan structure analysis on N- and O-glycoproteins and in exoglycosidase arrays on oligosaccharides.**

-  Fucose on N- and O-glycosylated proteins
-  Hydrolyzes  $\alpha$ 1-2,3,4-linked fucose
-  1-2h reaction
-  No need for co-factors

### Defucosylation Workflow



FucosEXO hydrolyzes fucose on glycoproteins under native conditions and displays high activity in a pH ranging from 6 to 8. Since the enzymatic activity

is not dependent on any co-factors or additives, the enzyme is compatible with most sample preparations, including for LC-MS analysis. Depending on the nature

of the substrate, FucosEXO hydrolyzes fucose on glycoproteins within 1 to 2 hours, while longer incubation may be required for very complex samples.

### Product Formats



**FucosEXO™**  
Lyophilized enzyme for defucosylation of glycoproteins



**Immobilized FucosEXO™**  
Immobilized enzyme for defucosylation of glycoproteins in spin columns

### FucosEXO™




PRODUCT	DESCRIPTION	ID
FucosEXO, 2000 units	Defucosylates 2 mg glycoprotein	G1-FM1-020



## Substrate Specificity of FucosEXO™

Fucose residues are attached to both N- and O-glycans with different linkages.  $\alpha$ 1-2,  $\alpha$ 1-3, and  $\alpha$ 1-4-linked fucose most commonly occurs in O-glycans and as antenna fucosylation of N-glycans, whereas the core  $\alpha$ 1-6-linked fucose is found as a modification of the N-glycan core. A panel of oligosaccharide substrates representing different linkages was here used to measure the release of fucose after FucosEXO incubation. FucosEXO was found to efficiently release  $\alpha$ 1-2,  $\alpha$ 1-3, and  $\alpha$ 1-4-linked fucose, whereas no activity was observed on  $\alpha$ 1-6-linked fucose residues.

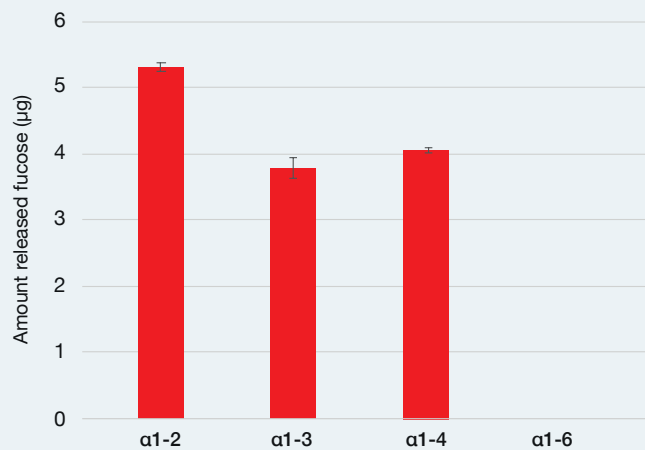

**FucosEXO substrate specificity.** The substrate specificity of FucosEXO was analyzed on equal molar amounts of synthetic oligosaccharides;  $\alpha$ 1-2 (2'-fucosyllactose),  $\alpha$ 1-3 (3-fucosyllactose),  $\alpha$ 1-4 (Lewis a) and  $\alpha$ 1-6 ( $\alpha$ 1-6 fucosylated chitobiose). Substrates were incubated with FucosEXO for 30 min at 37°C and the amount of released fucose was measured spectrophotometrically using an L-Fucose Assay Kit (Megazyme).

## Defucosylation of Native Glycoproteins using Immobilized FucosEXO™




Fucosylation of O-glycans is involved in the synthesis of functionally important glycan epitopes, including blood group antigens and the Lewis structures. The analysis of glycoproteins modified with such complex glycans can be challenging and requires specific and efficient enzymatic tools. Here, a glycoengineered TNFR protein carrying up to 11 O-glycans decorated with  $\alpha$ 1-2 and  $\alpha$ 1-3-linked fucose was treated with Immobilized FucosEXO, and the activity was compared to a commercially available  $\alpha$ -fucosidase. Within 1 hour, a complete removal of fucose was achieved using FucosEXO, while treatment with the  $\alpha$ -fucosidase only had a minor impact on the fucoses.

### Defucosylation using Immobilized FucosEXO.

Glycoengineered TNFR with core 1 and core 2 O-glycans decorated with both  $\alpha$ 1-2 and  $\alpha$ 1-3 fucose was incubated with bovine kidney fucosidase (1 h, 37°C) or Immobilized FucosEXO (1 h, RT). The processed protein was analyzed by reversed-phase LC-MS on a Waters™ BioAccord™ system equipped with a Waters™ BioResolve™ RP mAb column (2.1 × 50 mm).

### Glycoengineered TNFR (PNGase F treated)

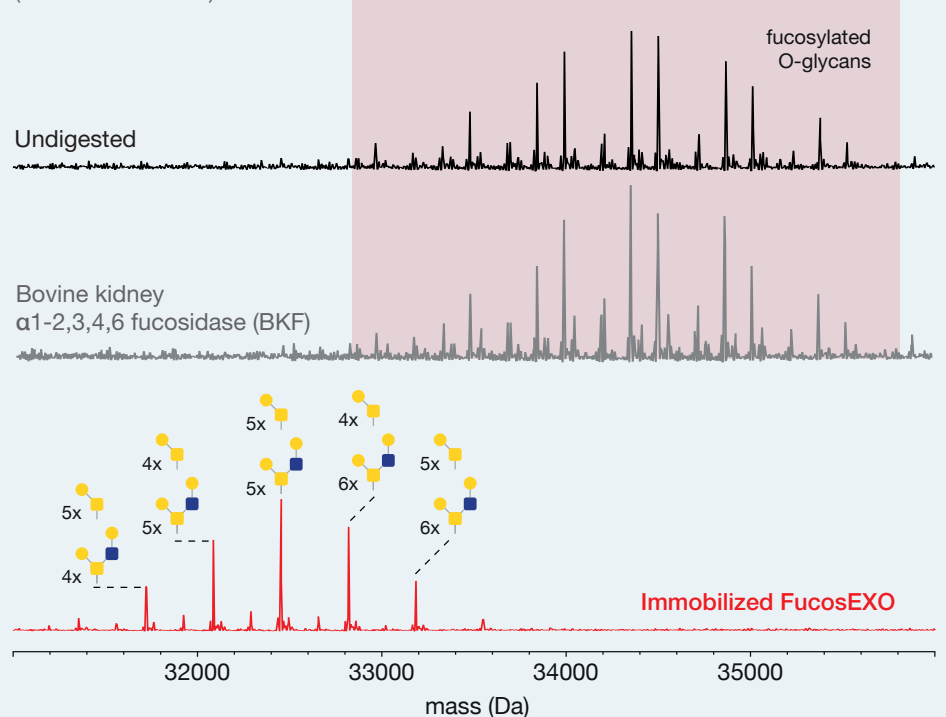

## Immobilized FucosEXO™



PRODUCT	DESCRIPTION	ID
Immobilized FucosEXO, Microspin 5 × 0.5 mg	Defucosylates 5 × 0.5 mg glycoprotein	G1-FM6-025
Immobilized FucosEXO, Microspin 10 × 0.5 mg	Defucosylates 10 × 0.5 mg glycoprotein	G1-FM6-050




# GalactEXO™

## Hydrolysis of $\beta$ 1-3,4-linked Galactose



GalactEXO is a  $\beta$ -galactosidase mix that hydrolyses  $\beta$ 1-3,4-linked terminal galactose in N- and O-glycosylated proteins. The enzyme is applicable for trimming of galactose on released glycan structures for exoglycosidase arrays or for the generation of antibodies with homogenous G0/G0F glycan profiles.

-  Galactose on N- and O-glycosylated proteins
-  Hydrolyzes  $\beta$ 1-3 and  $\beta$ 1-4-linked galactose
-  2h reaction
-  No co-factors needed

### Degalactosylation Workflow



GalactEXO is a mix of  $\beta$ -galactosidases that facilitates efficient hydrolysis of  $\beta$ 1-3,4-linked galactose on glycoproteins. The GalactEXO enzyme acts on both

N- and O-glycosylated structures for complete hydrolysis of terminal galactose within two hours. The enzyme can be used for intact glycoproteins as well as

on released glycans for applications in glycan sequencing or trimming of antibody Fc glycans.

### Product Formats



#### GalactEXO™

Lyophilized enzyme for hydrolysis of  $\beta$ 1-3,4-linked galactose



#### Immobilized GalactEXO™

Immobilized enzyme for hydrolysis of  $\beta$ 1-3,4-linked galactose in spin columns

### GalactEXO™



PRODUCT	DESCRIPTION	ID
GalactEXO, 2000 units	Degalactosylates 2 mg glycoprotein	G1-GM1-020





## Released Glycan Analysis with GalactEXO™

Released glycan analysis using exoglycosidases requires specific and complete hydrolysis to minimize errors in data interpretation. A labeled N-glycan library

was incubated with the SialEXO sialidase and GalactEXO for 1 h and analyzed with UHPLC. The HILIC-FLD separations show the complete removal of sialylated

and galactosylated structures, allowing the remaining peaks to be easily identified as G0, G0F, or G0F with bisecting GlcNAc.

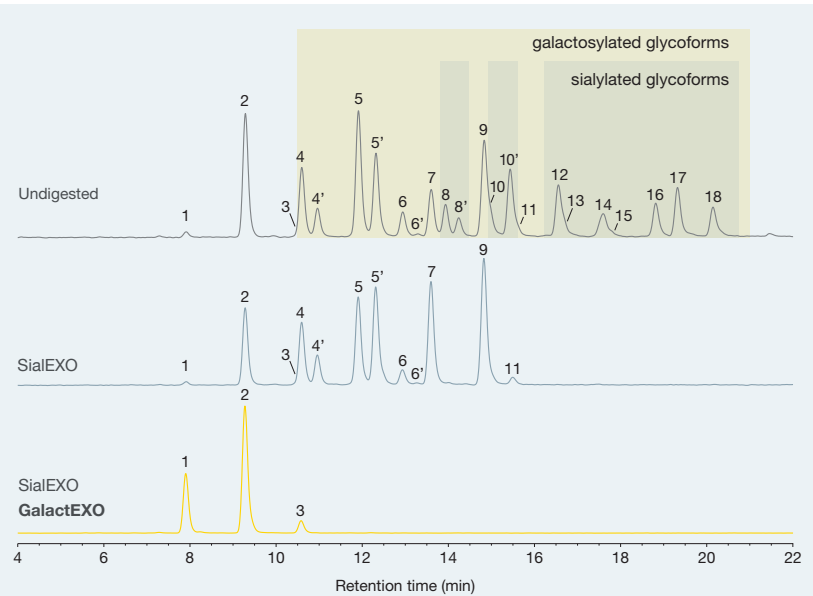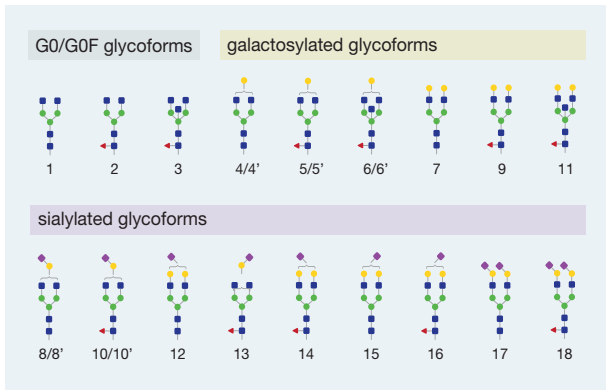


**Released glycan analysis.** HILIC-FLD UHPLC chromatograms of a 2-AB labeled glycan library analyzed undigested (top), after treatment with SialEXO (middle) or SialEXO and GalactEXO (bottom). Analyzed on Thermo Scientific Vanquish Duo UHPLC system equipped with Thermo Scientific Accurcore 150 Amide HILIC column (2.1 × 150 mm).

## G0 Antibodies with Immobilized GalactEXO™



The presence of galactose impacts antibody effector functions, and to study this mode of action, a clean removal of galactose residues may be necessary. The  $\beta$ 1-4 galactosidase activity of Immobilized GalactEXO was demonstrated on trastuzumab with a 30 min incubation prior to subunit digestion using FABRICATOR and analysis by LC-MS. The UV chromatograms clearly show the absence of the GalactEXO enzymes (yellow shading) in the Immobilized GalactEXO sample. The shift to G0F glycoforms can be seen in the mass spectra.

**Degalactosylation of trastuzumab.** Deconvoluted mass spectra of the Fc/2 fragment of trastuzumab treated with immobilized GalactEXO (bottom panel). The mAb was digested with FABRICATOR, separated by RP-HPLC (Waters™ BioResolve™ RP mAb (2.1 × 50 mm) and analyzed by ESI-Q-TOF MS (Bruker Impact II).

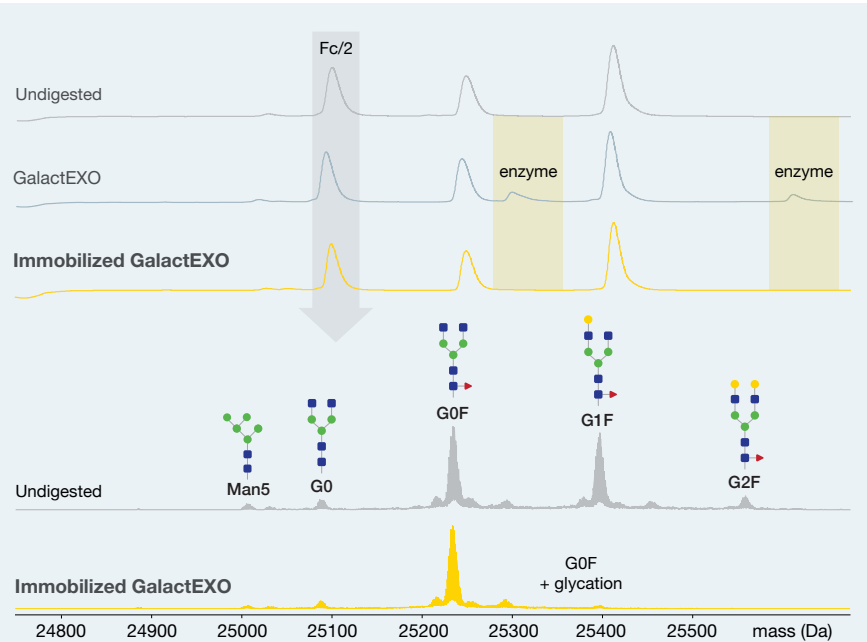

## Immobilized GalactEXO™



PRODUCT	DESCRIPTION	ID
Immobilized GalactEXO, Microspin 5 × 0.5 mg	Degalactosylates 5 × 0.5 mg glycoprotein	G1-GM6-025
Immobilized GalactEXO, Microspin 10 × 0.5 mg	Degalactosylates 10 × 0.5 mg glycoprotein	G1-GM6-050

# GalNAcEXO™

## Hydrolysis of $\alpha$ -linked GalNAcs



The GalNAcEXO enzyme is an  $\alpha$ -N-Acetylgalactosaminidase for efficient hydrolysis of terminal GalNAc residues on glycoproteins. GalNAcs linked to serine or threonine, referred to as Tn antigen, are quickly and efficiently hydrolyzed by GalNAcEXO, and the enzyme also displays activity on  $\alpha$ 1-3-linked terminal GalNAcs. The enzyme is a valuable tool that reduces sample heterogeneity for the analysis of complex O-glycoproteins that carry  $\alpha$ -linked GalNAc residues as immature truncated core 1.

- 

$\alpha$ -linked GalNAc residues on O-glycoproteins
- 

Tn antigen and  $\alpha$ 1-3-linked terminal GalNAcs
- 

4 h reaction
- 

No need for co-factors

## Removal of GalNAc Residues



GalNAcEXO hydrolyzes GalNAcs on glycoproteins under native conditions and is active in pH ranging from 6.0 to 7.6.

No co-factors or special buffers are required. Depending on the nature of the substrate, the GalNAcEXO reaction on

glycoproteins is complete in two hours while more complex samples may require overnight incubation.

## Product Formats



**GalNAcEXO™**  
Lyophilized enzyme for hydrolysis of GalNAc residues



**Immobilized GalNAcEXO™**  
Immobilized enzyme for hydrolysis of GalNAc residues in spin columns

### GalNAcEXO™




PRODUCT	DESCRIPTION	ID
GalNAcEXO, 2000 units	Hydrolyzes GalNAcs from 2 mg glycoprotein	G1-NA1-020



## GalNAcEXO™ Performance on Complex Substrates

Etanercept is a complex biopharmaceutical carrying 13 O-glycan sites, of which an average of 9 to 10 are occupied. After N-glycan removal by PNGase F and trimming of O-glycans with SialEXO and GalactEXO, GalNAc-related peaks were observed in the mass spectra (top). After incubation with the GalNAcEXO enzyme, over 95% of the remaining α-GalNAcs were removed (bottom).

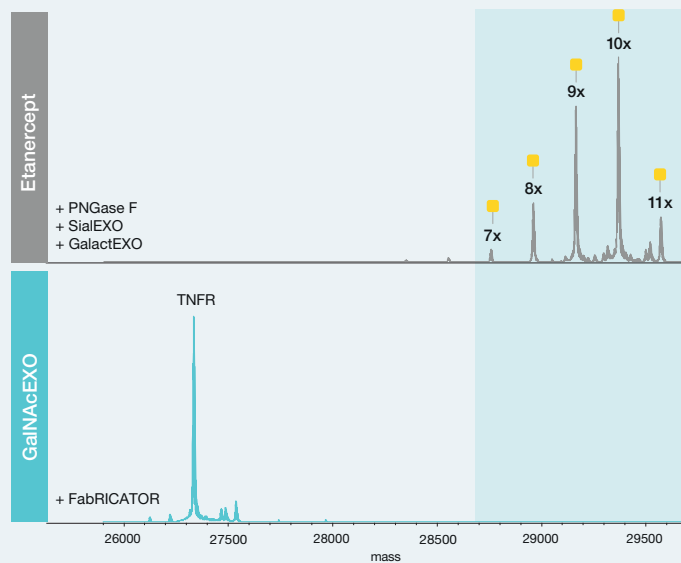

**GalNAcEXO activity on etanercept with truncated O-glycans.** Deconvoluted mass spectra of partially deglycosylated etanercept (top) and after treatment with GalNAcEXO (O/N at 37°C, bottom). Samples were FabRICATOR-digested, separated by RP-HPLC (Waters™ BioResolve™ RP mAb, 2.7 μm, 2.1 × 100 mm) and analyzed by ESI-Q-TOF MS (Bruker Impact II).


## Deglycosylation of the C1 Inhibitor using GalNAcEXO™



During the production of O-glycosylated biopharmaceuticals, Tn antigens may appear as a result of incomplete processing of the core GalNAc residues. This appears as a repeating mass shift with a 203 Da HexNAc unit difference. To confirm the presence of Tn antigens on the heavily (26 sites) O-glycosylated

C1 inhibitor, the following workflow was applied: first, the N-glycans were removed using PNGase F and core 1 O-glycans were removed using OglyZOR and SialEXO. A pattern of repeating peaks was observed in the mass spectra. The remaining GalNAc residues were completely removed using GalNAcEXO,

leaving one single peak. This workflow allowed for quantification of the number of core GalNAcs and reduced the remaining heterogeneity, facilitating the confirmation of the intact mass of the protein as well as revealing truncated variants.

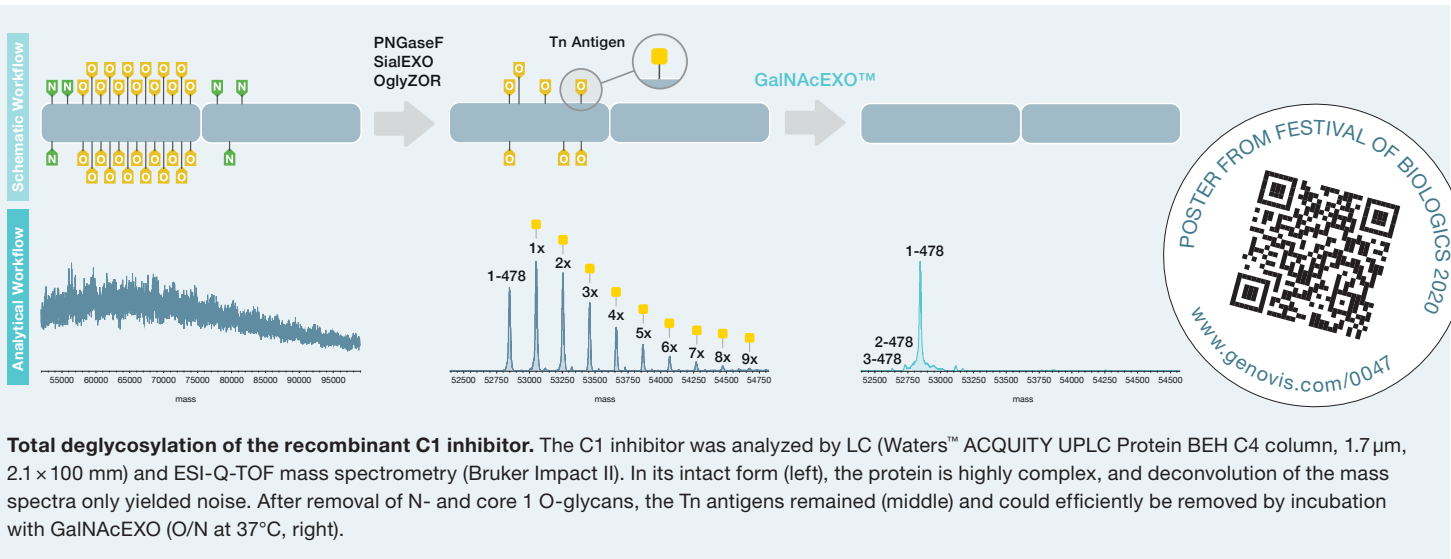

**Total deglycosylation of the recombinant C1 inhibitor.** The C1 inhibitor was analyzed by LC (Waters™ ACQUITY UPLC Protein BEH C4 column, 1.7 μm, 2.1 × 100 mm) and ESI-Q-TOF mass spectrometry (Bruker Impact II). In its intact form (left), the protein is highly complex, and deconvolution of the mass spectra only yielded noise. After removal of N- and core 1 O-glycans, the Tn antigens remained (middle) and could efficiently be removed by incubation with GalNAcEXO (O/N at 37°C, right).


## Immobilized GalNAcEXO™







PRODUCT	DESCRIPTION	ID
Immobilized GalNAcEXO, Microspin 5 × 0.5 mg	Hydrolyzes GalNAcs from 5 × 0.5 mg glycoprotein	G1-NA6-025
Immobilized GalNAcEXO, Microspin 10 × 0.5 mg	Hydrolyzes GalNAcs from 10 × 0.5 mg glycoprotein	G1-NA6-050
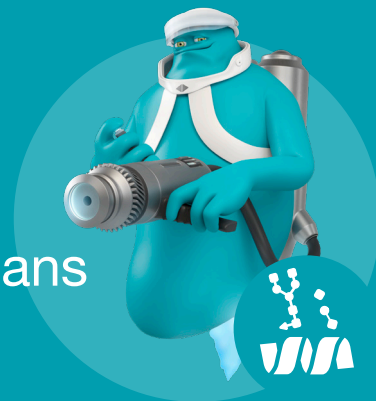
# OmniGLYZOR™

## Hydrolysis of N- and Mucin-type O-glycans



OmniGLYZOR contains a mixture of immobilized enzymes for the removal of N- and simple mucin-type O-glycans on antibodies, fusion proteins and other glycosylated proteins. Removal of glycans is widely used to reduce heterogeneity to facilitate analysis of the protein by for example mass spectrometry. Deglycosylation can also be used to study the functional role of the glycans.

-  N- and simple mucin-type O-glycans on glycoproteins
-  Hydrolysis of N- and O-glycans from glycoproteins
-  1-4 hour reaction
-  RapiGest™\* and PNGase F included

### Deglycosylation Workflow



OmniGLYZOR contains enzymes required for the removal of N-glycans and the most commonly occurring mucin-type O-glycans, namely mono- and disialyl core 1 and Tn antigen ( $\alpha$ -GalNAc). PNGase F, O-glycosidase, sialidase and  $\alpha$ -GalNAcase activities are included in OmniGLYZOR.

The glycoprotein sample is incubated with the OmniGLYZOR resin in a microspin column for 1-4 h under native reaction conditions. The deglycosylated glycoprotein is then easily collected by a centrifugation step. Some N-glycosylation sites are poorly or not accessible to

PNGase F unless the substrate protein is denatured. Remaining N-glycans can therefore be removed by an additional deglycosylation step using lyophilized PNGase F under denaturing conditions. Lyophilized PNGase F and RapiGest SF is included in the OmniGLYZOR kit.

### Product Formats



#### OmniGLYZOR™

A mix of immobilized enzymes in spin columns for deglycosylation of glycoproteins carrying N- and simple O-glycans

#### OmniGLYZOR™

PRODUCT	DESCRIPTION	ID
OmniGLYZOR, Microspin 5 x 50-100 $\mu$ g	Deglycosylates 5 x 50-100 $\mu$ g glycoprotein	G3-OM6-005
OmniGLYZOR, Microspin 10 x 50-100 $\mu$ g	Deglycosylates 10 x 50-100 $\mu$ g glycoprotein	G3-OM6-010




## Efficient Removal of N- and Mucin-type O-glycans

We demonstrate the performance of OmniGLYZOR using two therapeutic proteins as substrates, the heavily glycosylated Fc-fusion protein etanercept and erythropoietin (EPO). Due to the glycan heterogeneity, intact analyses resulted in complex mass spectra. By incubating the

samples with OmniGLYZOR for 1 hour at 37°C, the N- and O-glycans were efficiently removed as indicated by single peaks corresponding to the unmodified protein. For EPO, a minor amount of O-glycans modified with acetylated sialic acids were left on the protein since


those structures were inefficiently hydrolyzed by OmniGLYZOR. Both substrate proteins were also treated with another commercially available deglycosylation product according to the manufacturer's recommendation (O/N incubation at 37°C), and the data are shown for comparison.

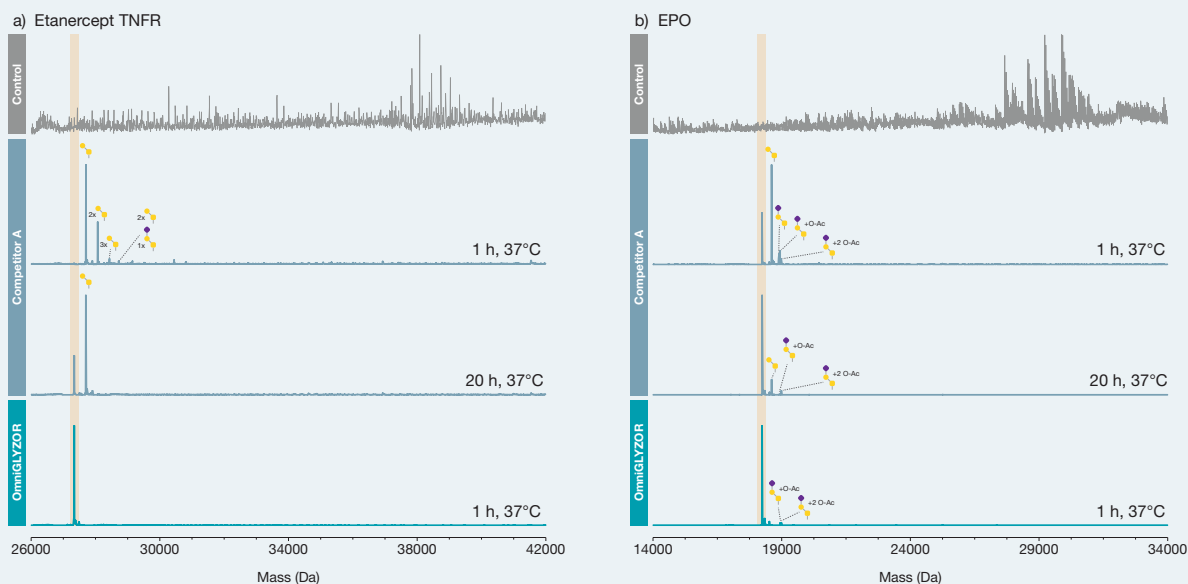

**Complete hydrolysis of N- and O-glycans by OmniGLYZOR.** Deglycosylation of a) etanercept and b) EPO. The substrate proteins were incubated on OmniGLYZOR Microspin columns for 1 h at 37°C, or with another commercially available deglycosylation product for 1 h and 20 h at 37°C. To simplify the analysis of etanercept, the deglycosylated protein was digested with FabRICATOR to separate the O-glycosylated TNFR domain from the Fc fragment. The resulting subunits were analyzed by reversed-phase LC-MS using a Waters™ BioAccord™ LC-MS system. EPO was analyzed in the same way in its intact state. The peaks corresponding to the fully deglycosylated substrate proteins are shaded in orange.

## Deglycosylation under Denaturing Reaction Conditions

The C1-inhibitor is a human plasma-derived biotherapeutic that is modified with 6 N-glycans and up to 28 O-glycans consisting of mostly sialyl core1 structures. Without pretreatment of this highly heterogeneous protein, reversed-phase LC-MS analysis yielded a complex mass spectrum impossible to interpret in detail. Using OmniGLYZOR Microspin columns under native conditions, all O-glycans were efficiently removed within 3 hours. However, between 1 and 3 N-glycans remained on the protein. These inaccessible N-glycans were removed by an additional deglycosylation step under denaturing conditions using the lyophilized PNGase F and MS-compatible RapiGest™ SF surfactant included in the OmniGLYZOR kit. A complete removal of all glycans was observed, with the exception of the minor amount of core 2 O-glycans present on the molecule.

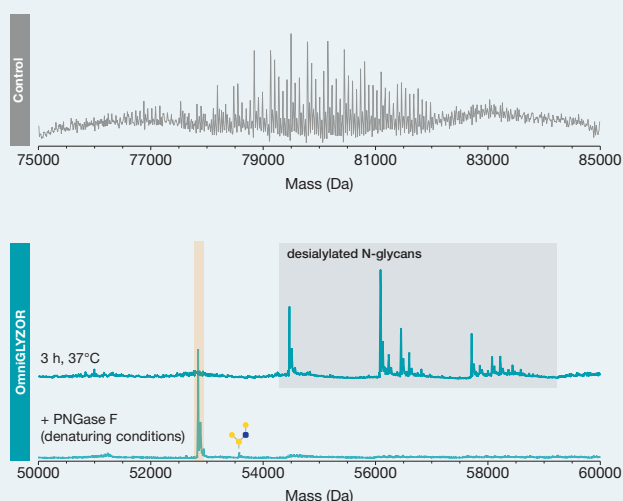

**Complete deglycosylation of plasma-derived human C1 inhibitor.** The protein was analyzed by reversed-phase LC-MS using a Waters™ BioAccord™ LC-MS system in its untreated state (top), after deglycosylation on an OmniGLYZOR Microspin column under native conditions (middle) and after additional deglycosylation using PNGase F under denaturing conditions (bottom). The peak corresponding to the deglycosylated substrate protein is shaded in orange.





\* RapiGest™ SF Surfactant from Waters Corporation is included in OmniGLYZOR™. RapiGest™ is a trademark of Waters Corporation.

# GlycOCATCH®

## Enrichment of O-glycoproteins



**GlycOCATCH is an affinity resin for purification of mucin-type O-glycosylated proteins and peptides. The affinity resin is based on inactive OpeRATOR enzyme that has been engineered to bind O-glycosylated proteins and peptides with high affinity. The applications of GlycOCATCH include specific enrichment or removal of O-glycoproteins and peptides, glycomics, studies of complex samples and characterization of biopharmaceuticals.**

-  Mucin-type O-glycosylated proteins and peptides
-  Binds to glycoproteins and peptides carrying mucin-type O-glycans.
-  2h workflow
-  Desialylation using SialEXO (included) increases the binding

### Affinity Purification Workflow



GlycOCATCH is provided in a spin column format to enable easy-to-use enrichment of mucin-type O-glycoproteins. The GlycOCATCH protein binds to O-glycoproteins with high affinity, and

due to the strong interaction, elution is performed using 8M urea. Alternatively, the elution can be performed with the included OpeRATOR enzyme, that digests the bound O-glycoproteins N-terminally of

the O-glycan sites. In this case, peptides are recovered under native conditions. Since non-sialylated O-glycans are bound much more efficiently to the resin, the sialidase mix SialEXO is included.

### Product Formats



#### GlycOCATCH®


Enrichment resin for affinity purification of mucin-type O-glycoproteins

### GlycOCATCH®

PRODUCT	DESCRIPTION	ID
GlycOCATCH	Enrichment of 0.2 mg O-glycoprotein	G3-OC6-002

## Selective Purification of O-glycosylated Peptides using GlycOCATCH®

Here, the O-glycan specificity of GlycOCATCH was studied at the peptide level. An O-glycosylated peptide was spiked into a mixture of non-glycosylated peptides and loaded onto GlycOCATCH, the resin was washed, and the bound peptides were eluted with 8 M urea. After purification, the samples were separated and analyzed by RP-LC-MS. The data show selective purification of the glycodrosocin peptide, carrying a core 1 O-glycan.

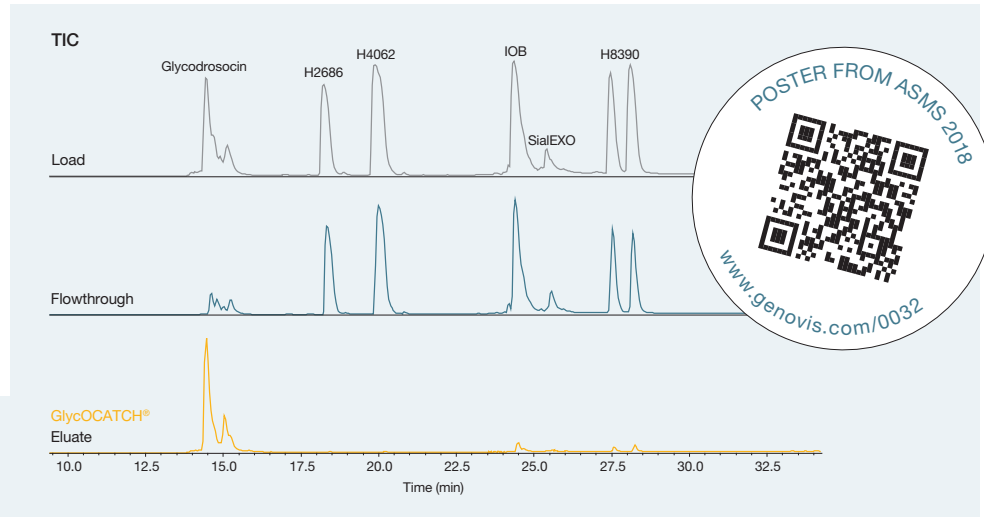

### Affinity purification using GlycOCATCH.


Selective purification of an O-glycosylated peptide (glycodrosocin) carrying a core 1 O-glycan.

## Affinity Purification of O-glycosylated proteins using GlycOCATCH®

GlycOCATCH can be used for specific enrichment of O-glycosylated peptides and proteins prior to MS analysis. The

strong binding ability of GlycOCATCH was effectively used to enrich O-glyco-


sylated proteins from a complex matrix such as human serum.

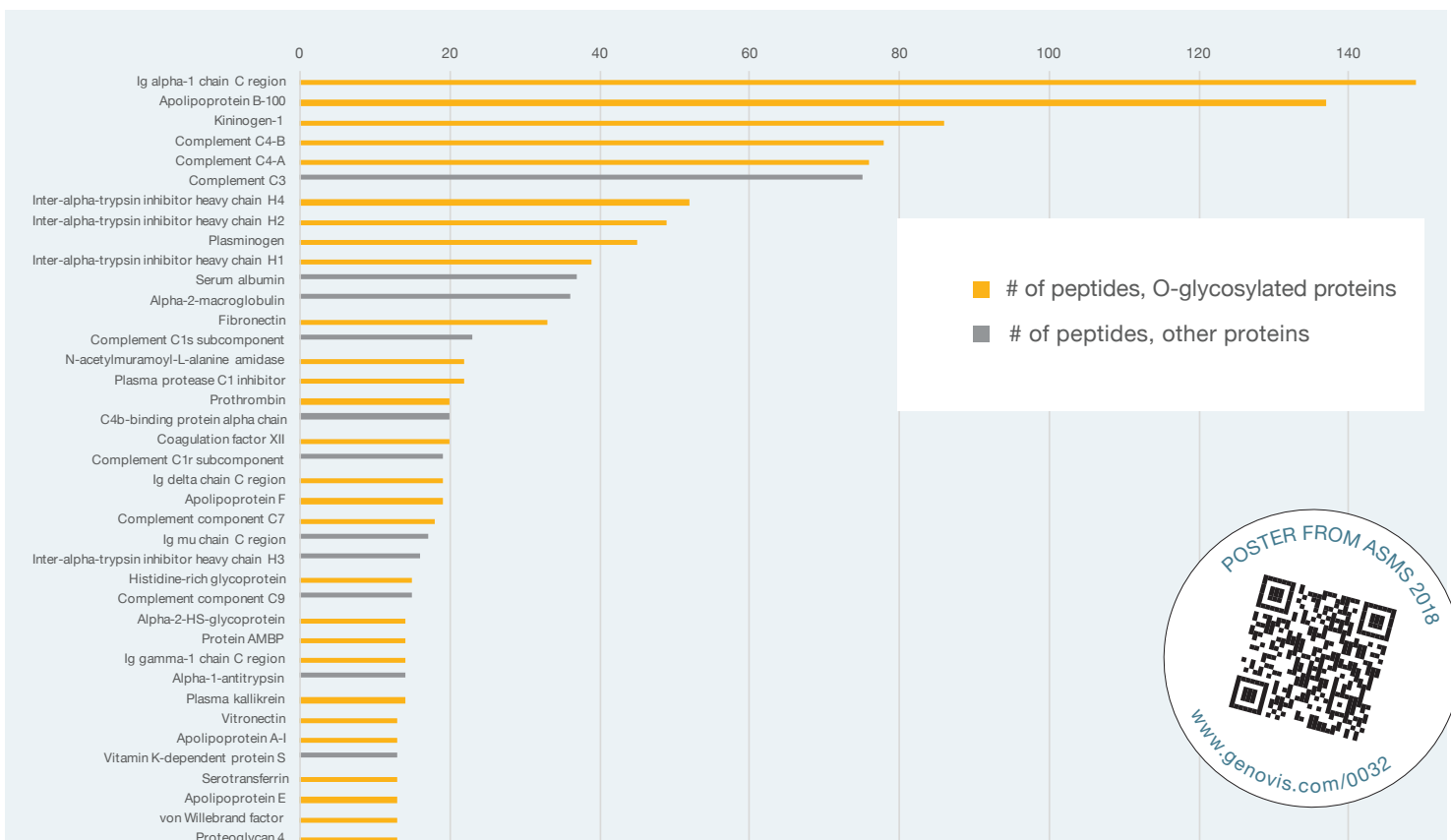

**Enrichment of O-glycoproteins from human serum.** Human serum was treated with SialEXO and incubated with the GlycOCATCH resin. The eluted proteins were denatured, reduced, alkylated and trypsin digested. Peptides were analyzed using RP-LC-MS/MS, and searched against the Swiss Prot database. Identified O-glycosylated proteins and other proteins with > 12 matching peptides and a MASCOT score > 200 are listed.

# Anti-FabRICATOR®

## Detection of the FabRICATOR® Enzyme

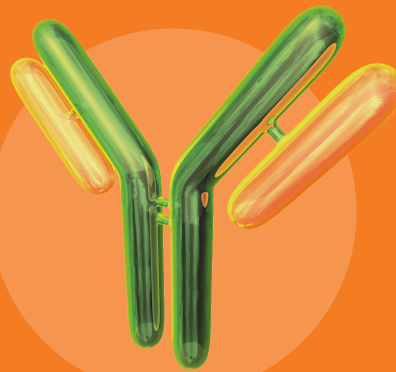

**Anti-FabRICATOR (Anti-IdeS) is a goat polyclonal antibody detecting the FabRICATOR (IdeS) enzyme. It is available in three formats; affinity purified, biotin-conjugated affinity purified and purified on protein G. Anti-FabRICATOR can be used to detect the IdeS enzyme in Western blot and ELISA assays, as for example in preclinical gene therapy applications.**

-  Detects the FabRICATOR (IdeS) enzyme
-  Source: goat
-  Polyclonal
-  Concentrations ranging from 0.5 mg/ml to 4 mg/ml

### Product Formats



**Anti-FabRICATOR®  
Affinity Purified**  
Affinity purified goat polyclonal antibody



**Anti-FabRICATOR®  
Affinity Purified Biotin Conjugated**  
Affinity purified, biotin-conjugated goat polyclonal antibody




**Anti-FabRICATOR®  
Protein G Purified**  
Protein G-purified goat polyclonal antibody

### Anti-FabRICATOR® Affinity Purified Performance

The Anti-FabRICATOR products are sensitive reagents for detecting the FabRICATOR (IdeS) enzyme. The Anti-FabRICATOR Affinity Purified displays a higher affinity to the FabRICATOR analyte as compared to the Anti-FabRICATOR Protein G Purified. This was demonstrated by a sandwich ELISA, where the plate

was coated with Anti-FabRICATOR Protein G Purified or Anti-FabRICATOR Affinity Purified respectively, prior to the binding of the FabRICATOR enzyme. The biotinylated Anti-FabRICATOR Affinity Purified was used as the secondary antibody, and the readout was detected by HRP-labeled streptavidin.

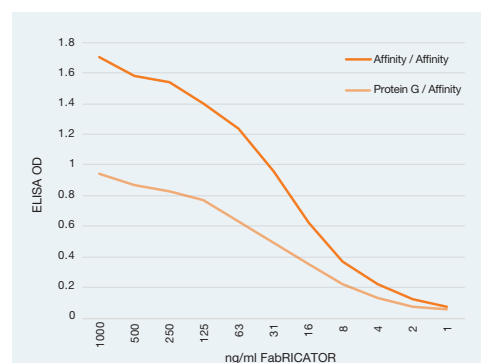

**Performance of Anti-FabRICATOR Affinity Purified.** To display the difference in performance between the Anti-FabRICATOR Protein G Purified and the Anti-FabRICATOR Affinity Purified, wells were coated with either Anti-FabRICATOR Protein G Purified or Anti-FabRICATOR Affinity Purified at 2 µg/ml. The FabRICATOR enzyme was added in dilutions from 1 µg/ml to 1 ng/ml. Anti-FabRICATOR Affinity Purified Biotin Conjugated was then added at 1 µg/ml. HRP-Streptavidin was added in a 1:2000 dilution, followed by visualization using ABTS.

### Anti-FabRICATOR®

PRODUCT	DESCRIPTION	ID
Anti-FabRICATOR Affinity Purified	Frozen solution of 0.1 ml, concentration 1 mg/ml	A3-AF2-010
Anti-FabRICATOR Affinity Purified, Biotin Conjugated	Frozen solution of 0.1 ml, concentration 0.5 mg/ml	A3-AF3-005
Anti-FabRICATOR Protein G Purified	Frozen solution of 0.1 ml, concentration 4 mg/ml	A3-AF1-010

## Legal and Disclaimers

All rights reserved. Genovis products may be covered by one or more patents, trademarks and copyrights owned or controlled by Genovis AB. For more information about commercial rights, please contact the Genovis team at [licensing@genovis.com](mailto:licensing@genovis.com).

Genovis products are intended for research use only. They are not intended to be used for therapeutic or diagnostic purposes in humans or animals. All goods and services are sold subject to Genovis' General Terms and Conditions of Sale.

### **FabRICATOR®**

This product is provided under an exclusive world-wide intellectual property license from Hansa Biopharma AB derived from international publication WO03051914, including granted US Patent No US 7,666,582 and granted European Patent No 1458861. The license encompasses IdeS from *Streptococcus pyogenes* for biotechnical industrial applications which are neither therapeutic nor diagnostic, other than the following exception which is included within the license: digesting IgG *in vitro* in clinical samples for diagnostic purposes.

### **POROS® included in FabRICATOR® HPLC**

This product is provided under an intellectual property license from Life Technologies Corporation. The transfer of this product is conditioned on the buyer using the purchased product solely in research or for analytical testing during manufacturing, all conducted by the buyer, excluding contract research or any fee for service research, and the buyer must not (1) use this product or its components

for (a) diagnostic, therapeutic or prophylactic purposes; (b) testing, analysis or screening services, or information in return for compensation on a per-test basis; and/or (c) sell or transfer this product or its components for resale, whether or not resold for use in research. For information on purchasing a license to this product for purposes other than as described above, contact Thermo Fisher Scientific, 5826 Newton Drive, Carlsbad, CA 92008 USA or [outlicensing@thermofisher.com](mailto:outlicensing@thermofisher.com). The trademark POROS® trademark is the property of Life Technologies Corporation.

### **FabRICATOR® Magic**

The trademarks Thermo Scientific™, KingFisher™ and BindIt™ are the property of Thermo Fisher Scientific Inc. and its subsidiaries.

### **SiteClick™ included in GlyCLICK®**

This product is provided under an intellectual property license from Life Technologies Corporation. The transfer of this product is conditioned on the buyer using the purchased product solely in research conducted by the buyer, excluding contract research or any fee for service research, and the buyer must not (1) use this product or its components for (a) diagnostic, therapeutic or prophylactic purposes; (b) testing, analysis or screening services, or information in return for compensation on a per-test basis; or (c) manufacturing or quality assurance or quality control, and/or (2) sell or transfer this product or its components for resale, whether or not resold for use in research. For information on purchasing a license to this product for purposes other than

as described above, contact Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or [outlicensing@thermofisher.com](mailto:outlicensing@thermofisher.com). The trademark SiteClick™ is the property of Life Technologies Corporation.

### **Linker-toxin Payloads included in GlyCLICK®**

This product is provided under an intellectual property license from Glykos Finland Ltd. The transfer of this product is conditioned on the buyer using the purchased product solely in research conducted by the buyer. The buyer must not (1) use this product or its components for (a) diagnostic, therapeutic or prophylactic purposes; or (b) manufacturing or quality assurance or quality control, and/or (2) sell or transfer this product or its components for resale, whether or not resold for use in research.

### **CaptureSelect™**

Thermo Scientific™ CaptureSelect™ resin from Thermo Fisher Scientific Inc. and its subsidiaries. Thermo Scientific and CaptureSelect are trademarks of Thermo Fisher Scientific Inc. and its subsidiaries.

### **TransGLYCIT™**

TransGLYCIT contains a technology licensed from GlycoT Therapeutics, LLC.

### **Immobilized PNGase F Denaturing**

RapiGest™ SF Surfactant from Waters Corporation is included in Immobilized PNGase F Denaturing. RapiGest™ is a trademark of Waters Corporation.

© Genovis AB



The SmartEnzymes portfolio consists of enzymes and technologies designed to improve the efficacy and throughput in analytical or preparative workflows for complex biopharmaceuticals such as antibodies, Fc-fusion proteins, complex glycoproteins and antibody conjugates.



**FabRICATOR®**  
Below hinge  
digestion of IgG



**FabALACTICA®**  
Above hinge  
digestion of human IgG1



**FabDELLO™**  
Above hinge  
digestion of human IgG1



**FabRICATOR® Z**  
Below hinge  
digestion of mouse IgG



**FabULOUS™**  
Above hinge  
digestion of IgG



**GingisKHAN®**  
Above hinge  
digestion of human IgG1



**GlycINATOR®**  
Hydrolysis of  
IgG Fc glycans



**IgGZERO®**  
Hydrolysis of  
complex Fc glycans



**GlyCLICK®**  
Site-specific  
conjugation of native IgG



**TransGLYCIT™**  
Transglycosylation  
of IgG



**GingisREX®**  
Arginine-specific  
protein digestion



**PNGase F**  
Hydrolysis of  
N-glycans



**OpeRATOR®**  
O-glycan-specific  
protein digestion



**OglyZOR®**  
Hydrolysis of  
core-1 O-glycans



**SialEXO®**  
Hydrolysis of  
sialic acids



**FucosEXO™**  
Hydrolysis of  
 $\alpha$ 1-2,3,4 fucose



**GalactEXO™**  
Hydrolysis of  $\beta$ 1-3,4-  
linked galactose



**GalNAcEXO™**  
Hydrolysis of  
 $\alpha$ -linked GalNAcs



**OmniGLYZOR™**  
Hydrolysis of N- and  
Mucin-type O-glycans



**GlycOCATCH®**  
Enrichment of  
O-glycoproteins



**Anti-FabRICATOR®**  
Detection of the  
FabRICATOR enzyme

**EMEA & Asia – Genovis AB**

Box 790, SE-220 07 Lund, Sweden  
Customer service: +46 46 10 12 30  
Order phone: +46 46 10 12 30  
Order fax: +46 46 12 80 20  
order@genovis.com

**US & Canada – Genovis Inc.**

245 First Street, Suite 1800, Cambridge, MA 02142, USA  
Customer service: 617-444-8421  
Order phone (toll free): 855-782-0084  
Order fax: 858-524-3006  
orders.us@genovis.com