

Cell culture media For better lives

CELLiST™ Glycyl-L-Tyrosine Dihydrate

Enhance Your Biologics Manufacturing with Superior Solubility and Cost-Efficiency

Shortcut to Official Channel






# **CELLIST™ Glycyl-L-Tyrosine Dihydrate**

Enhance Your Biologics Manufacturing with Superior Solubility and Cost-Efficiency

## **Overview**

CELLISTTM offers a comprehensive solution for all your biologics manufacturing requirements. CELLISTTM Glycyl-L-Tyrosine Dihydrate (Gly-L-Tyr) provides exceptional solubility, enhancing the overall available tyrosine source, which is a critical component for cell growth and productivity. Gly-L-Tyr solubility surpasses that of traditional L-Tyrosine. Incorporating CELLISTTM Gly-L-Tyr into the manufacturing processes removes the need for two separate feeding lines, resulting in a simple, robust and cost-effective workflow. This product was developed by Ajinomoto Group, leveraging years of experience in amino acids science and bioprocess design, as well as high-quality raw materials production capabilities.

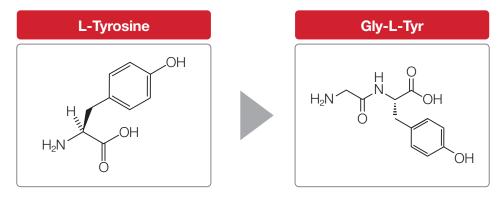

# **Advantages**


#### Gly-L-Tyr is an effective tyrosine source with high solubility

CELLIST<sub>TM</sub> Glycyl-L-Tyrosine Dihydrate is a highly-soluble and effective tyrosine source, which is an essential ingredient for biologics manufacturing. Unlike other tyrosine sources, such as L-Tyrosine disodium salt dihydrate (L-Tyr+2Na), which requires an alkaline pH and presents challenges during dissolution, Gly-L-Tyr provides optimal solubility without compromising pH conditions.

	L-Tyrosine	Gly-L-Tyr
Solubility at 25°C	0.479 g/L	36.9 g/L


### **Streamlined Process and Cost-Cutting**

Choosing CELLISTTM Gly-L-Tyr simplifies your manufacturing process. Unlike L-Tyrosine or L-Tyr+2Na, which requires an additional, separate feed source at higher pH due to reduced solubility, Gly-L-Tyr allows for a one-line feeding process, resulting in substantial decrease in running costs. Gly-L-Tyr can be easily supplemented to the main feed, dissolved at neutral pH and provides consistent tyrosine source to your cell culture, without introducing any pH fluctuations to the culture.

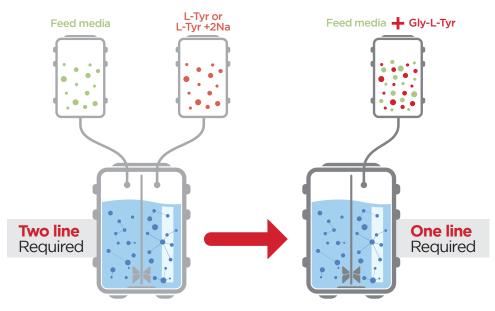

### Optimized tyrosine source for improved productivity

The addition of Gly-L-Tyr to the process optimizes tyrosine source supply to the cell culture and results in significant improvement to the performance and productivity of the feb-batch process. The graph below shows the increase in productivity (titer) for a fed-batch CHO cell culture.

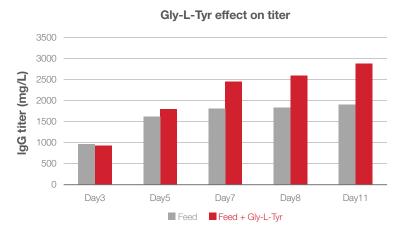

### **Animal Origin Free (AOF) Guarantee**

CELLIST<sub>TM</sub> prioritizes the quality and safety of our products. Our Gly-L-Tyr is guaranteed to be Animal Origin Free (AOF). To ensure the highest standards, we employ a fermentation process for L-Tyrosine, avoiding the use of animal-origin materials found in extraction processes, such as those utilized by other manufacturers.

### **Shorter Lead-time and Stable Supply**

We value our customers' time constraints and production needs. With a lead time of just 2 weeks for in-stock items, and 6 months for out-of-stock products, we are able to provide you with extremely rapid delivery and support.

# **Specifications**

Model No.	Glycyl-L-Tyrosine Dihydrate
Format	Powder
Storage Condition	Room Temperature
Shelf Life	24 months
Item Description	Powder, chemically defined and animal derived compound free





# **Global CELLIST Solution Centers**



CELL ST | Solution Center

Ahmedabad, India Mail: parshadg.hirapara.ut3@asv.ajinomoto.com CELL ST | Solution Center

Shanghai, China

Mail: mingqi\_fan@ajinomoto.com.cn

#### **North America**

JSR Life Sciences, LLC. 1280 North Matilda Ave. Sunnyvale, CA 94089 Tel: +1-408-543-8800 E-Mail: bioprocess.us@jsrlifesciences.com

#### Europe

JSR Micro NV. Technologielaan 8 B-3001 Leuven, Belgium
Tel: +32 (0)16 832 832 E-Mail: bioprocess.eu@jsrlifesciences.com

Eat Well, Live Well.

